

# **Table of Contents**

Everything you need for a well stocked laboratory or educational environment



- 4 Laboratory Equipment
- 16 Sampling & Testing Instruments
- 44 Workstation Essentials
- 55 Glassware
- 62 Sample Collection & Transport
- 72 Liquid Sample Handling
- 86 Sample Handling
- 97 Tubes & Vials

- 100 Tube Racks & Storage
- 119 Cold Products & Storage
- 135 Microscopes, Slide Storage and Accessories
- 146 Safety Products & Dispensers
- 155 Sealing & Labeling
- 163 Charts & References

# **Laboratory Equipment**

High-quality equipment and instruments for a vast range of protocols



Centrifuges

**Stirrers** 

Mixers

**Shakers** 

# **Clinical Centrifuge**

# Powerful clinical centrifugation in compact bench-top size

Powerful bench-top centrifugation for clinical use in stat labs, medical clinics and offices, veterinarian offices, hospitals and education.

### **Superb Function**

- Powerful, Adjustable Speed (500 to 6,000 RPM) can be set in RPM or RCF
- Programmable with memory for up to 8 programs
- Adjustable time (30 seconds to 60 minutes or continuous)
- Rapid acceleration and gentle electronic braking (deceleration time of <60 seconds)
- Digitally controlled brushless motor for accurate speed control and no motor maintenance required

### **Optimal Design**

- Easy to read Digital back-lit display
- Tool-less rotor removal for convenient loading and cleaning
- Clear lid/Viewing porthole to allow non-contact speed verification by non-contact tachometer
- Bench-top size, compact footprint for workspace use allows user to see rotor movement for added visual safety
- Includes automatic shut-off feature for sample imbalance and tilt

### **Built-in Safety**

- Safety features help ensure it will not operate when the lid is open, unit is jammed or unit is out of balance
- Manual lid unlock allows users to open lid in case of a power failure and prevents loss of samples

### **IVDR** Designation

- Class A, Rule 5a (Article 47, Annex VIII) per IVDF 2017/746

### **Applications include:**

- Clinical chemistry

120775

Grey/Purple ■ ■

- Blood collection prep -

Immunology

Clinical microbiology - Blood b

Blood banking


- Pathology



Rotor fits a multitude of tube sizes with the included tube adapters and spacers

- 3 10 mL Blood Collection Tubes
- 5 mL MicroCentrifuge Tubes
- 15 mL Conical Centrifuge Tubes



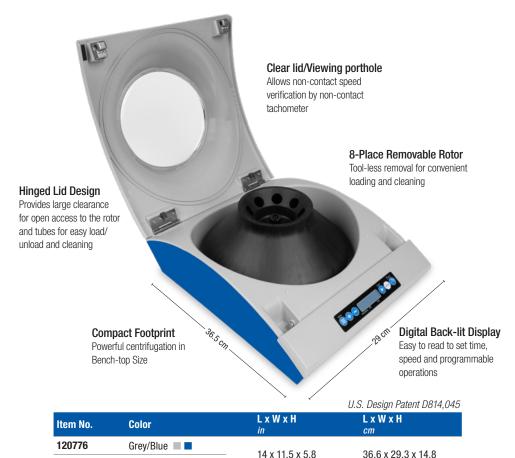

Specifications			
Speed Range	500 – 6,000 ı	pm	
RCF	4,668 x g (max RCF calculated) with typical 15 mL tube		
Capacity	8-place 15 mL tube rotor: 120 mL total maximum liquid volume		
Run Time Ranges	30 seconds – 60 minutes or continuous		
Weight	4.5 kg (9.93 l	b)	
Operating Temp.	Range 2°C to	40°C	
Max. Relative Humidity	Non-condensing, 80% for temp up to 31°C decreasing linearly to 50% R.H. at 40°C		
Electrical	Meets global power and efficiency requirements with grounded cord sets for NA, Japan, EU, UK and Australia		
Power Adapter	Input	100-240 VAC, 2.5A, 50/60 Hz	
Rating	Output	12 VDC, 12.5 A, 150W	
	Input	12 VDC	
Centrifuge Rating	Current (Under load)	Approx. 7.7 A	
	Power	Approx. 93W	
Certifications	C€	(SGS) (1 √RoHS UK X	

	_	Compilant	
Package Includes			
Clinical Centrifuge			
8-Place Rotor			
8 Tube Adapters			
24 Tube Spacers			

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs









# **MagFuge®**

A high-speed mini centrifuge and magnetic stirrer in one unit





Two Units

in One!



**High-Speed Mini Centrifuge** 

High-speed Centrifuge is easy to use and designed to deliver maximum speed (12,500rpm/9,783x g) and flexibility with rapid acceleration and fast controlled breaking

Ideal for most protocols requiring fast spins including PCR, micro- filter cell separation and HPLC. Versatility to handle 0.2, 0.5, 1.5, 2.0 and 5.0 mL tube sizes

### Intuitive One-Dial Programming

- Set speed and time with ease:
- Adjustable speed from 500 to 12,500 rpm
- Speed setting for RPM or RCF
- Adjustable time from 30 seconds to 60 minutes

### Safe to Use

- Safety features help ensure it will not operate when the lid is open, a rotor is jammed, or the unit is out of balance

Item No.	Color	L x W x H in	L x W x H cm
120581	Grey/Purple ■■	0.0 7.0 4.7	23.5 x 18.5 x 12
120582	Grey/Blue ■	9.2 x 7.3 x 4.7	

Replacement parts			
120622	MagFuge® Rotor, 12 x 1.5/2.0 mL tube		
120623	MagFuge® Rotor, 6 x 5.0 mL tube		
120624	MagFuge® Stirrer Rotor		
120682	MagFuge® Silicone Mat		
120527	Power Adapter		
120015	Tube Adapters, 0.2 mL, 0.5 mL, Pack of 12 (6 ea.)		
HS100502	Universal Power Cord Set		

### Patent Information

Community Design No. 003515337-0001-0002 Euro. Pat. App. 3 246 088

US Design Patent No. D814,045 US Patent No. 10471439



High-performance Magnetic Stirrer is designed using high quality rare earth magnets creating a robust magnetic coupling that significantly reduce the chance of spin-outs

**Magnetic Stirrer** 

Perfect for work with sensitive samples and protocols up to 3L requiring aggressive mixing for full dilution

### Easy to Use

- Set stir mode, speed and time with easy one-dial programming:
- Select Mode: clockwise, c. clockwise and oscillate
- Adjustable speed from 50 to 2,500 rpm
- Continuous and adjustable time settings from 30 seconds to 60 minutes

### **Durable Accessories Included**

 Re-usable, easy to clean silicone mats included to keep vessel stable and prevent spills and two NdFeB stir bars (40 x 12 mm) with chemical resistant PTFF coating included

Г	THE COating included		
Package Includes			
Centrifuge/Stirrer	Stirrer Rotor 12 Tube Adapters 0.2 mL & 0.5 mL		
5.0 mL Tube Rotor	2 Stir Bars, 35 mm x 12 mm, cylindrical shape		
1.5/2.0 mL Tube Ro	otor 2 Silicone Mats for Stirrer Lid		
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4			
Interchangeable Plugs			
Specifications: Cer	ntrifuge		
Speed Range	500-999 rpm (increments of 1); $1,000-5,000$ rpm (increments of 5); $5,000-12,500$ rpm (increments of 25)		
RCF	16 – 9,783 x g (max RCF calculated)		
D Ti D	00		

Kun Time Kanges	30 seconds – 60 minutes (5 second increments)		
Specifications: Ma	anotic Stirror		
•	·		
Speed Range	50 – 999 rpm (increments of 1); 1,000 – 2,500 rpm (increments of 5)		
Run Time Range	Stirrer Mode	30 seconds to 60 minutes or continuous (5 second increments)	
	Oscillating Stirrer Mode	1 second – 60 minutes or continuous	
	Oscillating Period	30 seconds – 15 minutes	
Canacity	3 Linius Diameter of vessel not to exceed 130 mm		



requirements, press and hold to

operate, release to end spin cycle

LxWxH

# **Gusto® High-Speed Mini Centrifuge**

**Maximum speed with minimal space** 

This compact centrifuge is ideal for protocols calling for fast spins with a higher RCF requirement

Made of durable ABS with UV-Resistant coating

- Intuitive Push button programming and bright LCD screen make it easy to use

### **Quiet Operation**

- Brush-less motor enables quiet operation and low vibration

### **Safety Features**

Item No.

- Lid safety feature prevents opening while rotor is turning and clear polycarbonate allows user to see rotor movement for added visual safety
- Includes automatic shut-off feature for sample imbalance and tilt

Color

HS10050	Blue	8 x 6.7 x 4.5	20.3 x 17.1 x 11.4	
Replaceme	nt parts			
HS100502	Universal Power Cord Set			
HS100503	Gusto® Rotor 12 x 1.5/2.0 mL Tube (with Cover and Knob)			
HS100504	Gusto® Rotor 4 x 8 Strip 0.2 mL Tube (with Knob)			
120015	Tube Adapters, 0.2, 0.5 mL, Pack of	12 (6 ea.)		
HS100506	Gusto® Rotor Cover (for 12 x 1.5/2.0	mL Tube Rotor Only)		
120527	Power Adapter			

LxWxH

Specifications	
Speed Range	500-12,500 rpm
Run Time Ranges	1 to 30 minutes
RCF	15 - 9,783 x g (max RCF calculated)









HS CA	>
-------	---

Package Includes
Gusto® High-Speed Mini Centrifuge
12-Place 1.5/2.0 mL Tube Rotor
32-Place 0.2 mL Tube Rotor (singles, or 4 x 8 Strip)
12 Tube Adapters 0.2 mL
12 Tube Adapters 0.5 mL
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs



# Sprout® Plus Mini Centrifuge

Mini but Mighty - ideal for small runs

### Easy to Use

Tool-free SnapSpin rotors and quick-spin options make the centrifuges fast and easy to use

# **Ergonomic Design**

- The lid was functionally designed to fit in the palm of the hand for one handed lock and unlock allowing for quick tube extraction
- Large lid size accommodates longer and skirted tubes and
- The housing is made from ABS plastic and lid is polycarbonate plastic enabling the unit to be easily cleaned

### Safety Measures

Sprout® Plus

- Fitted with a braking function, the unit quickly stops spinning when powered off or when lid is opened
- A weighted, rubberized base eliminates potential movement and dampens vibration and noise

sprout	sprout
Specifications	

Specifications	
Speed Range	Approx. 6,000 rpm (fully loaded)
RCF	2,000 x g



Item No.	Color	LxWxH in	L x W x H cm
120610	Green	- 6 x 5 x 4.3	15.3 x 12.8 x 10.9
120611	Blue	- 0 X 3 X 4.3	13.3 x 12.0 x 10.9
Replacement part	ts		
HS100011	Sprout® Plus Rotor SnapSpin 6 x 1.5/2.0 mL Tube		
HS100012	Sprout® Plus Rotor SnapSpin 2 x 8 Strip, 0.2 mL		
120148	Sprout® Plus Lid O-Ring, <i>Pack of 10</i>		
120491	Sprout® Plus Power Adapter		
120015	Tube Adapters, 0.2 mL, 0.5 mL, Pack of 12 (6 ea.)		

# Package Includes

Sprout® Plus Mini Centrifuge

6-Place 1.5/2.0 mL Tube Rotor

16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)

Accessory Kit Included

6 Tube Adapters 0.2 mL

6 Tube Adapters 0.5 mL

4 Lid 0-Rings

Storage Case for Rotors and Adapters

Tube Rack

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs







# **Need tubes?**

See our portfolio of 0.2 mL, 0.5 mL, 1.5 mL Microcentrifuge Tubes





# **Portable Centrifuge Kit**

Spin samples on-the-go

The Portable Centrifuge Kit offers a convenient way to spin down and prepare samples outside of a conventional laboratory environment

### **Immediate Field Processing**

- Process samples in remote areas before deterioration can take place
- Allows for immediate on-site testing for faster results
- Made of durable ABS plastic for easy clean in field environments

### Easy to Use

- Centrifuge is equipped with tool-less SnapSpin rotors for easy change out
- Simple operation with lock/unlock activation

### **Portable Convenience**

- All items in kit are housed in a handled DuraPorter® Transport Box for portability
- The handled box has an air-tight seal for optimal sample protection

See page 64 to learn more about the DuraPorter® Transport Box

Item No.	
120568	Portable Centrifugation Kit, 100-240 VAC, 50/60Hz Universal Plug + 12V Car Adapter
120569	Portable Centrifuge Transport Case with Inserts - Only
120570	12V Sprout® Plus Centrifuge Car Adapter*
	*Adaptor can be used with Sprout® Dlue Mini Contribuge

\*Adapter can be used with Sprout® Plus Mini Centrifuge

# Package Includes

Sprout® Plus Mini-Centrifuge

6-Place 1.5/2.0 mL Tube Rotor

16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)

6 Tube Adapters 0.2 mL

6 Tube Adapters 0.5 mL

4 Lid 0-Rings

Storage Case for Rotors and Adapters

**Tube Rack** 

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs









# DuraPorter® Transport Box

**Tube Rack Insert** 

(12-Place 1.5/2.0 mL, 14-Place 0.2 mL)

Transport Box Dividers

Sprout® Plus Cradle

Portable 12V Sprout® Plus Centrifuge Car Adapter

# **Magnetic Induction Stirrer**

Strong, motor-less and maintenance free

### **Superior Functional Design**

- Advanced magnetic coil technology drives the spinning of the stir bar maintaining continuous, constant speed for reproducible results
- Magnetic coupling is optimized by using alternating current (AC) as its driving force strong enough to stir viscous solutions
- Motor-less operation generates minimal heat ideal for stirring temperature sensitive samples and environments for long periods
- With no moving parts, produces less vibration and quieter operation over traditional magnetic stirrers

### **Intuitive Operation**


- Multiple stirring modes are easy to program
- Bright digital display precisely controls time and speed
- Quickly meets stirring speeds with rapid acceleration and braking (less than 10 seconds)




Item No.	Color	L x W x H	L x W x H
120584	Grey/Purple	9.7 x 7.7 x 3	24.5 x 19.5 x 7.5
120585	Grev/Blue	9.1 x 1.1 x 3	24.0 X 19.0 X 7.0

Replacement parts		
120738	MI Stirrer Silicone Mat	
120739	MI Stirrer Universal Power Adapter	



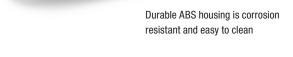



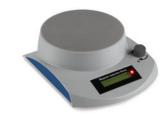

Specifications		
Speed Range	50-2,000 rpm variable, 1 rpm increments	
Run Time Range	30 seconds to 60 minutes, variable, 1 second increments or continuous	
Stirrer Modes	Single clockwise	
	Single counter clockwise	
	Auto reverse rotation 1 to 99 cycles 50% duty	
Capacity	3 L plus	
Mixing Surface	165 mm diameter, silicone mat 160 mm diameter	









Package Includes	
Rare Earth, NdFeB Magnetic Induction Stirrer	
2 Stir Bars 35 mm x 12 mm Cylindrical Shape	
1 Silicone Anti-Skid Plate Cover	

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs



The Mini Magnetic Stirrer is an affordable, dependable solution ideal for general research, academic, life sciences, industrial and clinical lab settings

### **Durable Design**

- Constructed of chemical resistant ABS and PC blend
- Compact footprint saves valuable space on benchtop
- Lightweight for easy transport

### **Ultimate Function**

- Powerful magnetic coupling ensures accurate alignment, reducing spin-outs during acceleration or deceleration
- Performs routine stirring applications with low vibration and quiet operation

U.S.and E.U. Utility Patent Pending U.S. Design Patent No. D814,045

Item No.	Color	L x W x H in	L x W x H cm
120155	Grey/Purple ■■	7.2 x 5.2 x 1.85	18.3 x 13.2 x 4.7
120595	Grev/Blue	7.2 X 3.2 X 1.03	10.3 X 13.2 X 4.7

Replacement	Part		
120815	Mini Stirrer Universal Power Adapter		



# **Specifications**

**Speed Range** 200-2,200 rpm (Infinitely adjustable)

Capacity Up to 3 L



C E VROHS COMPILIANT CA
-------------------------

<sup>ac</sup>	kage	e Inc	lud	es

Mini Magnetic Stirrer

2 Stir Bars 35 mm x 9 mm, Octagonal Shape

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs









# **Stir Bars & Stir Bar Retriever**

- Magnetic bar is made from ALNICO, one of the strongest permanent magnets available
- The sintered manufacturing process reduces the chance of decoupling (spin outs), encapsulation in PTFE protects both magnet and what it is stirring
- Chemically inert and easy to clean
- Can withstand temperatures up to 250°C, making them suitable for use with both unheated and heated stirrers
- Available in two surfaces:

Stir Bars

Item No.

120547

120548

- Low friction octagon shape with a center pivot ring helps create a strong mixing vortex, even with more viscous solutions
- Smooth shape for better coupling when using concave bottom glassware

LxW (in)

0.4 x 0.2

 $0.8 \times 0.3$ 

Pack of 10

L x W (mm)

10 x 5

20 x 8

25 x 8

28 x 8

30 x 8

35 x 9

10 x 5

25 x 8

30 x 8

35 x 8

 The retrieval bar can recover stir bars even in hot solutions when stirring and mixing is no longer required

Octagonal Shape	Smooth Shape
35 x 9 mm	
35 x 8 mm	35 x 8 mm
30 x 8 mm	30 x 8 mm
28 x 8 mm	
25 x 8 mm	25 x 8 mm
f 20 x 8 mm	
(10 x 5 mm)	(10 x 5 mm)

Stir Bar Retriever

120549 Octagonal 1 x 0.3 120550 with pivot ring 1.1 x 0.3 120551 1.2 x 0.3 120552 1.4 x 0.4 120894  $0.4 \times 0.2$ 120895 1 x 0.3 Smooth 1.2 x 0.3 120896 120897 1.4 x 0.3

**Shape** 

Stir Bar Retriever Pack of 5

Item No.	L (in)	L (cm)
120553	11.8	30

# **Lab Mats**

Two Lab Mat options to keep your workstation safe, quiet, clean and organized





# **Vortexer**

# Multi-purpose, maximum convenience

### **Innovative Design**

- One mixing head does it all with this sleek, modern and compact designed vortex mixer
- Hold slot configurations for multiple sizes of tubes, vessels and assay plates
- Center pad for quick hand-held mixing

### Easy to Use

- Control switch with three operation selections (intermittent, off, continuous)
- Variable Speed Rotary Dial provides pre-set speed control selection during shaking mode

### **Durable Construction**

- ABS material resist most commonly used chemicals and is easy to clean
- Stable elastomer feet dampen vibration and prevent movement during operation



Center pad area for hand held mixing



Assay plate x 1 (to be used with retention cords)



25 mL, 50 mL conical tubes x 2 (horizontally mounted)



15 mL conical tubes x 2 (horizontally mounted)



5~mL tubes x 2 (horizontally mounted in the same retention space as 15 mL)



0.5 mL microfuge tubes x 6 (vertically mounted)



1.5/2.0 mL microfuge tubes x 4 (vertically mounted)



0.2 mL microfuge tubes (strip tubes) x 24 (vertically mounted)

US Design Patent No. D787,085 US Utility Patent No. 9895670

			-	
Item No.	Speed Control	Power Input	L x W x H in	L x W x H cm
120209	Infinity Dial	110-120 VAC, 50/60 Hz, 55 W, US/NA I.I		
120212	Infinity Dial	220-240 VAC, 50 Hz, 44 W, CE	6.7 x 7.5 x 7.8	17.1 x 19 x 19.7
120214	Infinity Dial	220-240 VAC, 50 Hz, 44 W, UK	(with mixing head	(with mixing head
120210	Preset Dial	220-240 VAC, 50 Hz, 44 W, AUS	in place)	in place)
120318	Preset Dial	100-110 VAC 50/60 Hz 66 W(PSF) Japan III	_	

Replacement parts		
120327	Vortexer Multi-Head	
120328	Vortexer Central Pad, 3 pk	
120422	Vortexer Retention Cords, 3 pk	
120329	100-120 VAC US Cord 1.1	

Replacement parts		
120333	220-240 VAC AUS Cord	$\odot$
120332	220-240 VAC CE Cord	••
120331	220-240 VAC UK Cord	1
120330	100-110 VAC Japan Cord	11

Specifications	:
Speed Range	Up to 3,000 rpm
Orbit	3.7 mm
Control	3 Position Control — (top) intermittent, (middle) off, (bottom) continuous











# **Mini Vortexer**

**Small footprint with maximum power** 

Simple, touch start operation offers exceptional mixing performance with one hand **Robust Design** 

- Low profile design and small footprint saves valuable space
- Made of ABS and silicone materials that resist most commonly used chemicals
- Stable platform and elastomer feet ensure unit stays stable during operation

### **Powerful Function**


- Rapid vortex action with low activation force
- Quiet operation and low vibration

U.S. and E.U. Utility Patent Pending U.S. Design Patent No. D814,045

			, , , ,
Item No.	Color	L x W x H in	L x W x H cm
120567	Grey/Purple		
120598	Grey/Blue	5.1 x 4.3 x 3	13 x 11 x 7.5
120827	Black/Black		

Replacement parts		
120617	Mini Vortexer Silicone Cover, set of 2	
120814	Mini Vortexer Universal Power Adapter	



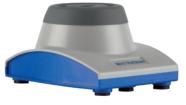



Specifications	
Speed	4,500 rpm
Orbit	6 mm
Capacity	Up to 50 mL tube









Package Includes	Packa	qe l	Incl	ud	es
------------------	-------	------	------	----	----

Mini Vortexer

2 Silicone Caps

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs









# **Digital Orbital Shaker**

# Shaker for a vast range of protocols

### Versatile Design

- Create your own matrix for vessels of all shapes and sizes with Hi-Lo Platform and Bungee Cords
- Instantly adjusts to weight load changes and unbalanced loads with nostop electronic motor feedback allowing for truly reproducible results
- Durable ABS construction suitable for use in incubators (excluding CO2) and cold rooms within the stated temperature range

# **High-Featured Operation**

- Suitable for wide range of protocols including very low start speed for sensitive sample agitation to high speeds for vigorous applications
- Soft start and shut down feature reduces spilling from full vessels
- Designed to handle heavy loads up to 10 lbs (4.5 kg)
- Overload protection shut-off prevents overloading

Item No.	LxWxH	LxWxH	Platform Int	erior W x D
	in	cm	in	ст
120460	12.4 x 7.8 x 12.2	31.4 x 19.7 x 31	11 x 11	27.9 x 27.9

Replacement parts		
120521	Digital Orbital Shaker Silicone Mat & Bumpers	
120524	Digital Orbital Shaker Bungie Cords, Pack of 8	
120491	Digital Orbital Shaker Universal Power Adapter	



Specifications	
Speed Range	20-300 rpm, 10 rpm increments
Orbit Diameter	19 mm (0.8")
Run Time Range	0 to 74 hours and 59 minutes; 1 minute increments or continuous
Maximum Load	4.5 kg (10 lbs)





Package Includes	
Digital Orbital Shaker	
Skid-Free Rubber Platform Mat	
8 Bungee Cord Set	
Universal Hi-Lo Platform with Soft Rubber Edge Bumpers	

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs

# Sampling & Testing Instruments

Reliable and accurate instrument range for diverse applications



Thermometers
Refractometers
Hydrometers & Jars
Pipette Controllers
Colony Counters & Accessories

# **THERMOMETERS**

We have a wide range of thermometers to measure and track temperatures for a wide range of applications. Our portfolio includes Dial (Bimetal), Digital and Glass Thermometer options to get the job done.

# **Dial Thermometers**

Dial thermometers don't require batteries, are enduring, and readings can be adjusted with the calibration screw on the back

Indicates temperature through a metal stem with attached circular scale. Application use varies with dial sizes, stem lengths and temperature ranges.

# **Dial Thermometers - 1" Diameter**

Ideal for many applications including food and beverage quality control, general lab use, and asphalt/concrete testing

- The small stem and face fits nicely into a shirt pocket or holder
- Includes a plastic sheath with pocket clip for storage

# Dial Thermometers - 1" Diameter

Item No.	Temperature Range	Subdivisions	Accuracy	Face Diameter (in/mm)	Stem Length (in/mm)	Stem Diameter (in/mm)
81125	25 to 125°F	1°F	±1°F			
81070	-40 to 70°C	1°C	±1°C	4 / 05	F / 407	0.1 / 0.5
81110	-10 to 110°C	1°C	±1°C	<u> </u>	5 / 127	0.1 / 3.5
81550	50 to 550°F	5°F	±5°F			

# Dial Thermometers - 1.75" / 2" Diameter

Ideal for many applications including food and beverage quality control, general lab use, and asphalt/concrete testing

- Larger dial size than common pocket thermometers for easier reading
- Features an adjustment screw for field calibration
- Includes a clip for container mounting

### Dial Thermometers - 1.75" / 2" Diameter

Item No.	Temperature Range	Subdivisions	Accuracy	Face Diameter (in/mm)	Stem Length (in/mm)	Stem Diameter (in/mm)
82125	25 to 125°F	1°F	±1°F	1.8 / 44.5	8 / 203	0.1 / 3.5
82110	-10 to 110°C	1°C	±1°C	1.8 / 44.5	8 / 203	0.1 / 3.5
82220	0 to 220°F	2°F	±2°F	1.8 / 44.5	8 / 203	0.1 / 3.5
82550	50 to 550°F	5°F	±5°F	1.8 / 44.5	8 / 203	0.1 / 3.5
82550DG	50 to 550°F /	2.5°F	±2.5°F	0 / 51	0./000	01/05
10 to 290°C	5°C	±5°C	2 / 51	8 / 203	0.1 / 3.5	



82220

# **Dial Thermometers - 3" Dial**

Ideal for many applications including soil testing, food and beverage quality control, and asphalt/concrete testing

- Constructed of stainless-steel case, bezel, and stem with a glass window
- Case is hermetically sealed to prevent lens fogging
- Includes a 1/2"NPT male thread on the back for threading into a pipeline
- Also features an adjustment screw for field calibration

### Dial Thermometers - 3" Diameter

Item No.	Temperature Range	Subdivisions	Accuracy	Face Diameter (in/mm)	Stem Length (in/mm)	Stem Diameter (in/mm)
82160-6	-40 to 160°F	2°F	±2°F	3 / 76	6 / 152	0.3 / 6.4
82160-12	-40 to 160°F	2°F	±2°F	3 / 76	12 / 305	0.3 / 6.4



Ideal for many applications including compost testing, food and beverage quality control, and asphalt/concrete testing

- Constructed of stainless-steel case, bezel, and stem with a glass window
- Case is hermetically sealed to prevent lens fogging
- Includes a 1/2"NPT male thread on the back for threading into a pipeline
- Also features an adjustment screw for field calibration

### **Long-Stem Dial Thermometers**

Long Otom	iong otom biar mormonotoro					
Item No.	Temperature Range	Subdivisions	Accuracy	Face Diameter (in/mm)	Stem Length (in/mm)	Stem Diameter (in/mm)
82200-36	0 to 200°F	2°F	±2°F	3 / 76	36 / 914	0.3 / 7.9
82100-36	0 to 100°C	1°C	±1°C	3 / 76	36 / 914	0.3 / 7.9
82200-48	0 to 200°F	2°F	±2°F	3 / 76	48 / 1220	0.3 / 7.9

# **Swivel-Head Digital Thermometers**

Digital stem thermometers are ideal for temperature monitoring in essentially any setting and application

Laboratories - easily monitoring reagents and liquids in test tubes, beakers, and flasks
 Food Processing - Refrigeration
 Horticulture

Both models feature a dual scale (°F / °C), fast 1-second update time, and a swivel head that can be moved a full  $180^{\circ}$  for easy reading in any position. Both the ON/OFF and scale selection buttons are located on the face for easy access.

The thermistor sensor stem is stainless steel, 0.138" (3.5 mm) diameter, and includes a protective case featuring a pocket clip. Includes a 1.5V alkaline battery (LR44) and features automatic power-off after 10 minutes of non-use.

# **Swivel-Head Digital Thermometers**

Item No.	Temperature Range	Resolution	Accuracy	Stem Length (in/mm)	Stem Diameter (in/mm)
83210-12	-58 to 302°F / -50 to 150°C	0.1°	±2.0°F / 1.0°C	12 / 305	0.1 / 3.5
83110	14 to 212°F / -10 to 100°C	0.1°	±1.0°F / 0.5°C	5 / 127	0.1 / 3.5





180° Head Rotation

83210-12

# **Glass Thermometers**

VEE GEE thermometers feature a SafetyBLUE fill. This is an organic performance liquid that is bio-degradable and exhibits all the qualities of mercury without the risk of contamination and hazards associated with mercury.

# **Serialized Glass Thermometers**

Precision serialized thermometers are ideal for use in laboratories, universities, food/beverage, environmental, wastewater and many other general-purpose applications

- Feature SafetyBLUE fill against a white background for easy reading
- Accuracy is equal to one subdivision

# **Serialized Glass Thermometers**

Item No.	Temperature Range	Subdivisions	Length (in/mm)	Immersion
80701	-30 to 120°F	2°F	12 / 300	Total
80702	-4 to 230°F	2°F	12 / 300	Total
80703	0 to 300°F	2°F	12 / 300	Total
80706	-35 to 50°C	1°C	12 / 300	Total
80707	-20 to 110°C	1°C	12 / 300	Total
80708	-20 to 150°C	1°C	12 / 300	Total
80901	-30 to 120°F	2°F	12 / 300	Partial
80902	-4 to 230°F	2°F	12 / 300	Partial
80903	0 to 300°F	2°F	12 / 300	Partial
80905	20 to 500°F	2°F	15.4 / 391	Partial
80906	-35 to 50°C	1°C	12 / 300	Partial
80907	-20 to 110°C	1°C	12 / 300	Partial
80908	-20 to 150°C	1°C	12 / 300	Partial
80910	-10 to 250°F	2°F	15.4 / 391	Partial

80901

# **Serialized Encapsulated Glass Thermometers**

**Encapsulation protects you from broken glass and liquid** 

- Feature FEP encapsulation; in the event of breakage the FEP coating contains the broken glass and liquid, preventing injury or contamination
- Feature SafetyBLUE fill against a white background for easy reading
- Accuracy is equal to one subdivision

### Serialized Encapsulated Glass

Item No.	Temperature Range	Subdivisions	Length (in/mm)	Immersion
80701E	-30 to 120°F	2°F	12 / 300	Total
80702E	-4 to 230°F	2°F	12 / 300	Total
80703E	0 to 300°F	2°F	12 / 300	Total
80704E	30 to 220°F	2°F	15.4 / 391	Total
80706E	-35 to 50°C	1°C	12 / 300	Total
80707E	-20 to 110°C	1°C	12 / 300	Total
80708E	-20 to 150°C	1°C	12 / 300	Total
80901E	-30 to 120°F	2°F	12 / 300	Partial
80902E	-4 to 230°F	2°F	12 / 300	Partial
80903E	0 to 300°F	2°F	12 / 300	Partial
80904E	30 to 220°F	2°F	15.4 / 391	Partial
80906E	-35 to 50°C	1°C	12 / 300	Partial
80907E	-20 to 110°C	1°C	12 / 300	Partial
80908E	-20 to 150°C	1°C	12 / 300	Partial

80702E

# **Serialized Armored Glass Thermometers**

# Armor protects the thermometer from rough handling and shock during use

- Includes nickel-plated brass armor
- Feature SafetyBLUE fill against a white background for easy reading
- Accuracy is equal to one subdivision

### **Serialized Armored Glass Thermometers**

Item No.	Temperature Range	Subdivisions	Length (in/mm)	Immersion
80701-A	-30 to 120°F	2°F	12 / 300	Total
80702-A	-4 to 230°F	2°F	12 / 300	Total
80703-A	0 to 300°F	2°F	12 / 300	Total
80706-A	-35 to 50°C	1°C	12 / 300	Total
80707-A	-20 to 110°C	1°C	12 / 300	Total
80708-A	-20 to 150°C	1°C	12 / 300	Total
80901-A	-30 to 120°F	2°F	12 / 300	Partial
80902-A	-4 to 230°F	2°F	12 / 300	Partial
80903-A	0 to 300°F	2°F	12 / 300	Partial
80905-A	20 to 500°F	2°F	15.4 / 391	Partial
80906-A	-35 to 50°C	1°C	12 / 300	Partial
80907-A	-20 to 110°C	1°C	12 / 300	Partial
80908-A	-20 to 150°C	1°C	12 / 300	Partial
80910-A	-10 to 250°F	2°F	15.4 / 391	Partial



**Armor protects the thermometer and encapsulation protects you** 

- Includes FEP encapsulation and nickel-plated brass armor
- In the event of breakage, the FEP coating contains the broken glass and liquid fill, preventing injury or contamination. The armor protects the thermometer from rough handling and shock during use
- Feature SafetyBLUE fill against a white background for easy reading
- Accuracy is equal to one subdivision

### **Armored Encapsulated Glass Thermometers**

Item No.	Temperature Range	Subdivisions	Length (in/mm)	Immersion
80701E-A	-30 to 120°F	2°F	12 / 300	Total
80702E-A	-4 to 230°F	2°F	12 / 300	Total
80703E-A	0 to 300°F	2°F	12 / 300	Total
80704E-A	30 to 220°F	2°F	15.4 / 391	Total
80706E-A	-35 to 50°C	1°C	12 / 300	Total
80707E-A	-20 to 110°C	1°C	12 / 300	Total
80708E-A	-20 to 150°C	1°C	12 / 300	Total
80901E-A	-30 to 120°F	2°F	12 / 300	Partial
80902E-A	-4 to 230°F	2°F	12 / 300	Partial
80903E-A	0 to 300°F	2°F	12 / 300	Partial
80904E-A	30 to 220°F	2°F	15.4 / 391	Partial
80906E-A	-35 to 50°C	1°C	12 / 300	Partial
80907E-A	-20 to 110°C	1°C	12 / 300	Partial
80908E-A	-20 to 150°C	1°C	12 / 300	Partial

# **Armor Cases for Glass Thermometers**

Nickel-plated brass armor cases protect liquid-in-glass thermometers from rough handling and shock during use

Armors feature a threaded screw cap with ring top for suspension where necessary

Item No.	Thermometer Length (in/mm)	Outside Diameter (in/mm)	Inside Diameter (in/mm)
83025A	12 / 305	0.4 / 9.5	0.3 / 7.5
83025N	15.5 / 393	0.4 / 9.5	0.3 / 7.5







# **Glass Thermometers: Be in the know....**

# **Column Separation**

The largest single cause for failure of liquid-in-glass thermometers in the lab is due to column separation. This can occur in transit or in the lab.

The simplest and safest method to reunite a separated column is to force the liquid down the capillary tube by using a centrifuge, if one is available with a cup deep enough to ensure that the centrifugal force is below the liquid column.

- Carefully insert the thermometer, bulb down, into the centrifuge. Place some cotton wadding at the bottom of the cup to prevent any damage to the bulb. Turn on the centrifuge - in just a few seconds all the liquid will be forced past the separation.
- If the cup is not deep enough and all the centrifugal force is not below the column, the column will split, forcing part of the liquid down. The remainder will be forced up, filling the expansion chamber.
- If a centrifuge is not available, the column can be reunited by holding the thermometer in an upright position and gingerly tapping the stem above the separation against the palm of your hand. As you gently tap the thermometer, observe the liquid above the separation until it breaks away from the wall of the capillary and runs down to join the main column.

# **Immersion Types**

Thermometers are one of the most important instruments used in the laboratory and are a relatively inexpensive means of reliably measuring temperature, provided they are properly used and not abused.

It is essential that the instruments be used in the manner prescribed if true temperature are to be ascertained. The following are some terms that should be understood and adhered to when true temperatures are to be obtained:

Total immersion - such a thermometer is designed to indicate temperatures correctly when the bulb and the entire liquid column are exposed to the temperature being measured, except for a minimal length emergent to be visible. Therefore, the thermometer is immersed into the medium up to the point of reading.

Partial Immersion - such a thermometer has a line around it at the immersion distance from the bottom. The thermometer indicates correctly when the bulb and liquid column to that line are exposed to the temperature being measured and the emergent stem is at ambient temperature. Therefore, the thermometer is immersed only up to the immersion line.

# REFRACTOMETERS

VEE GEE's benchtop and handheld refractometers provide users with highly accurate readings for a wide range of applications in the agricultural, food & beverage, chemical and manufacturing industries.

We offer a large breadth of scales for specific applications such as Brix, NaCl, Salinity, Specific Gravity of Seawater, Coolant, Battery Acid, and more.

Our portable line of optical and digital refractometers offers the convenience to perform sample analysis where and when you need. With added temperature control, our Abbe benchtop unit takes the accuracy of your Brix/Refractive Index analysis to the next level in the laboratory.

# X-Series Handheld Refractometers


**Compact & Portable** 

The X-Series handheld refractometers are exceptionally simple to operate and provide quick, accurate measurements for the concentration of nearly any aqueous solution in a wide range of applications.

Made durable with industrial-grade design features to provide years of reliable use even in the field. Rugged metal housing construction with slip-resistant rubber grip and a Non-Roll stand to protect from rolling off surfaces.

No need to recalibrate or apply correction factors in work environments with large temperature swings with the Automatic Temperature Compensation (ATC) feature included on Models BTX-1, BTX-20, STX-3, CTX-1 & CTX-2. The effective compensation range is 50-86° F / 10-30° C.



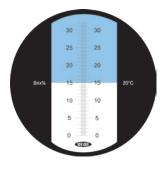

### **Accurate & Sharp Readings**

Large scale is best in class delivering clear, sharp figures, screen shading and incredible contrast for overall accuracy

# Quick Results

High quality optical glass prism and all-metal housing design allows the sample and prism to reach temperature equilibrium quickly

X-Series Specifications		
Sample Type	Transparent and Translucent	
Sample Quantity	>0.1 mL	
Prism	Optical Glass	
Prism Housing	Copper Alloy w/Nickel-Chrome Coating	
Included	Carrying Case, Plastic Transfer Pipet, Instruction Manual * ATC models include calibration screwdriver.	
Warranty	5-Year	

# **X-Series Handheld Refractometers**

For Food, Beverage, and Agriculture Applications

# **Brix Handheld Refractometers**

Item No.	Model No.	Scale	Range	Resolution	Accuracy	ATC	Approx Length (in/mm)
43001	BX-1	Brix	0 to 32%	0.2%	±0.2% Brix	No	6.5 / 165
43002	BTX-1	Brix	0 to 32%	0.2%	±0.2% Brix	Yes	6.5 / 165
43003	BX-2	Brix	28 to 62%	0.2%	±0.2% Brix	No	5.9 / 150
43004	SGX-1	Brix/Specific Gravity Beer Wort	0 to 32% 1.000 to 1.130	0.2% / 0.001	±0.2% / ±0.001	Yes	6.5 / 165
43005	BX-3	Brix	60 to 92%	0.2%	±0.2% Brix	No	5.5 / 140
43006	HXM-1	Brix,	Brix, 58.0 to 92%	0.5%	±0.5% Brix	Yes	5.8 / 148
		Honey	Water, 13 to 27%	1.0%	±1.0% Water		
		Moisture	Baume, 38 to 43°	0.5°	±0.5° Baume		
43007	BX-4	Brix	45 to 82%	0.2%	±0.2% Brix	No	5.5 / 140
43009	BX-10	Brix	0 to 10%	0.1%	±0.1% Brix	No	7.3 / 185
43011	BX-20	Brix	0 to 20%	0.1%	±0.1% Brix	No	7.3 / 185
43012	BTX-20	Brix	0 to 20%	0.1%	±0.1% Brix	Yes	7.3 / 185
43015	BX-50	Brix	0 to 50%	0.5%	±0.5% Brix	No	5.7 / 145



# **New!** Brix Field Kits

**Everything you need for in field testing at the Vineyard or Honey Bee Farm** 

Includes durable Carrying Case with Refractometer, Microfiber Cloth, Bottle for Distilled Water, Sample Cup, Sample Pipette, Calibration Screwdriver and Manual

Item No.	Model No.	Scale	Range	Resolution	Accuracy	ATC	Approx. Length (in/mm)	L x W x H of Kit (in/cm)
43002-KIT	BTX-1	Brix	0 to 32%	0.2%	±0.2% Brix	Yes	6.5 / 165	10.5 x 7 x 3 in.
43006-KIT	HXM-1	Brix,	Brix, 58.0 to 92%	0.5%	±0.5% Brix	Yes	5.8 / 148	26.7 x 17.8 x 7.6 cm
		Honey	Water, 13 to 27%	1.0%	±1.0% Water			
		Moisture	Baume, 38 to 43°	0.5°	±0.5° Baume	-		



# **Water Testing Refractometers**

For Water/Wastewater, Aquatic and Aquaculture Applications

Sodium Chloride And Salinity Handheld Refractometers

Item No.	Model No.	Scale	Range	Resolution	Accuracy	ATC	Approx.Length (in/mm)
43033	SX-1	NaCl	0 to 10%	0.1%	±0.1%	No	7.3 / 185
43035	SX-2	NaCl	0 to 28%	0.2%	±0.2%	No	6.5 / 165
43036	STX-3	Salinity	0 to 100 ppt	1.0 ppt	±1.0 ppt	Yes	7.3 / 185
		Specific Gravity of Seawater	1.000 to 1.070	0.001	±0.001	_	



43035/Model SX-2



# **New!** Water Testing Field Kit

Testing when and where you need to with this convenient kit

Includes durable Carrying Case with Refractometer, Microfiber Cloth, Bottle for Distilled Water, Sample Cup, Sample Pipette, Calibration Screwdriver and Manual

Item No.	Model No.	Scale	Range	Resolution	Accuracy	ATC	Approx. Length (in/mm)	L x W x H of Kit (in/cm)
43036-KIT	STX-3	Salinity	0 to 100 ppt	1.0 ppt	±1.0 ppt	Yes	7.3 / 185	10.5 x 7 x 3 in.
		Specific Gravity of Seawater	1.000 to 1.070	0.001	±0.001	_		26.7 x 17.8 x 7.6 cm




# **X-Series Handheld Refractometers**

For Light Industrial, Chemical, and Automotive Applications

# General Purpose and Industrial Chemical Handheld Refractometers

Item No.	Model No.	Scale	Range	Resolution	Accuracy	ATC	Approx. Length (in/mm)
43055	NDX-2	Refractive Index	1.435 to 1.520	0.001	±0.001	No	5.5 / 140
43062	CTX-1	Propylene Glycol Freezing Point	-58 to 32°F / -50 to 0°C	33.80°F 1°C	±1°C	Yes	5.9 / 150
	Ethylene Glycol Freezing Point	-50 to 0°C	1°C	±1°C			
		Battery Acid Specific Gravity	1.10 to 1.40	0.01	±0.01		
43064	CTX-2	Propylene Glycol Freezing Point	-50 to 32°F / -45.56 to 0°C -60 to 32°F / -51.11 to 0°C	2°F -16.67°C	±2°F	Yes	5.9 / 150
		Ethylene Glycol Freezing Point	-60 to 32°F	2°F	±2°F		
		Battery Acid Specific Gravity	1.10 to 1.40	0.01	±0.01		
43071	URX-60	Urea	0 to 60%	0.5%	±0.5%	No	5.7 / 145

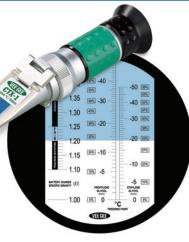

**43062/Model CTX-1**For specific scale images, visit product on our website.

# **Accessories**

Item No.	Description	Model
43000-DLXXX	Daylight Plate Kit	X-Series Handheld Refractometers
		Kit includes: Ring, Hinge, Plate, and Screw. Assembled
43RP-008	Soft Fabric Carrying Case	X-Series Handheld Refractometers
1111	Soft Vinyl Carrying Case	BX-1, BTX-1, BX-2, BX-3, BX-4, BX-50, SX-2, NDX-2, CTX-1, CTX-2
1111L	Soft Vinyl Carrying Case	BX-10, BX-20, BTX-20, SX-1, STX-3



# **Transfer Pipets**

# One-piece, unbreakable polyethylene plastic pipettes

# Transfer and deliver small volumes of liquid quickly and safely

- Bulk packed versions come in handy dispenser box
- Sterile versions available single wrapped or 20 wrapped at a time
- Lot-to-lot traceability on sterile versions



# **Transfer Pipets**

	Item No.	Vol. (mL)	Working Vol. (mL)	Approx. Length (in/mm)	Graduations (mL)	Packaging 500 per box	Sterile
	HS206371A	5	1	5.9 / 150	0.25	Bulk	
a)	HS206371B	5	1	5.9 / 150	0.25	20/bag	✓
	HS206371C	5	1	5.9 / 150	0.25	1/single wrap	✓
	HS206372A	5	2	6.1 / 155	0.5	Bulk	
b	HS206372B	5	2	6.1 / 155	0.5	20/bag	✓
	HS206372C	5	2	6.1 / 155	0.5	1/single wrap	✓
	HS206373A	7	3	6.1 / 155	0.5	Bulk	
C	HS206373B	7	3	6.1 / 155	0.5	20/bag	✓
	HS206373C	7	3	6.1 / 155	0.5	1/single wrap	✓



# MDX-Series Digital Refractometers

# **Compact, Durable and Accurate**

The MDX-Series Digital Refractometers feature IP65/67 ratings, making them ideal for use in the field and in a multitude of environments.

Delivering quick digital accuracy with high-end performance features like:

- Automatic Temperature Compensation (ATC)
- Continuous Loop Average Readings





# Easy to Clean

Prism cover and rubber grip add extra protection and are removable for easy cleaning



ATC Range	32 to 104°F / 0 to 40°C
IP Rating	IP67: detecting unit IP65: body
Prism	Optical Glass
Wavelength	589.3 nm LED
Sample Cell	Stainless Steel
Calibration Solution	Distilled Water (not included)
Min. Sample Volume	0.3 mL
Battery	1 X AAA, 1.5V
Dimensions (L x W x H in/mm)	4.8 x 2.3 x 1 / 121 x 58 x 25
Weight (lb/kg)	0.198 lb / 0.089 kg
Included	Hard shell case with protective lining, disposable pipette, (1) AAA Battery, and user manual
Warranty	1-Year Parts and Labor

# **MDX-Series Digital Refractometers**

For Food, Beverage, and Agriculture Applications

# Digital Brix/Sugar Refractometer

Item No.	Model No.	Scale	Range	Resolution	Accuracy
48101	MDX-101	Brix	0.0 to 50.0% Brix	0.1% Brix	±0.2% Brix
		Refractive Index	1.3330 to 1.4200 nD	0.0001 nD	±0.0003 nD
48102	MDX-102	Brix	0.0 to 90.0% Brix	0.1% Brix	±0.2% Brix
		Refractive Index	1.3330 to 1.5300 nD	0.0001 nD	±0.0003 nD
48103	MDX-103	Dextran	0.0 to 10.5%	0.1%	±0.2%
		Fructose	0.0 to 68.0%	0.1%	±0.2%
		Glucose	0.0 to 59.0%	0.1%	±0.2%
		Lactose	0.0 to 16.5%	0.1%	±0.2%
		Maltose	0.0 to 15.5%	0.1%	±0.2%
		Sucrose (Brix)	0.0 to 90.0%	0.1%	±0.2%

# **Digital Honey Refractometer**

Item No.	Model No.	Scale	Range	Resolution	Accuracy	
48301	MDX-301	Brix	0.0 to 90.0% Brix	0.1% Brix	±0.2% Brix	
		Moisture	5.0 to 38.0%	0.1%	±0.2%	
		Baume	33.0 to 48.0°	0.1°	±0.2°	
		Refractive Index	1.3330 to 1.5300 nD	0.0001 nD	±0.0003 nD	

# **Digital Wine Refractometer**

Item No.	Model No.	Scale	Range	Resolution	Accuracy	
48401	MDX-401	Brix	0.0 to 50.0% Brix	0.1% Brix	±0.2% Brix	
		PA (% v/v)	0.0 to 22.0%	0.1%	±0.2%	
		Oechsle (Oe)	0 to 150°	1°	±2°	
		KMW (Babo)	0.0 to 25.0°	0.1°	±0.2°	

# For Water/Wastewater, Aquatic and Aquaculture Applications

# Digital Sodium Chloride Refractometer

Item No.	Model No.	Scale	Range	Resolution	Accuracy	
48201	MDX-201	NaCl (%)	0.0 to 28.0%	0.1%	±0.2%	
		NaCl (ppt)	0 to 280 ppt	1 ppt	±2 ppt	
		Specific Gravity	1.000 to 1.217	0.001	±0.002	
		Refractive Index	1.3330 to 1.4200 nD	0.0001 nD	±0.0003 nD	

# **Digital Seawater Refractometer**

Item No.	Model No.	Scale	Range	Resolution	Accuracy	
48202	MDX-202	Salinity	0 to 100 ppt	1 ppt	±2 ppt	
		Chlorinity	0 to 57 ppt	1 ppt	±2 ppt	
		Specific Gravity	1.000 to 1.075	0.001	±0.002	
		Refractive Index	1.3330 to 1.4200 nD	0.0001 nD	±0.0003 nD	

# For Light Industrial, Chemical, and Automotive Applications

# Digital DEF (UREA)/Battery Acid Refractometer

Item No.	Model No.	Scale	Range	Resolution	Accuracy
48602	MDX-602	DEF (Urea) 0.0 to 51.0%		0.1%	±0.2%
		Refractive Index	1.3330 to 1.4200 nD	0.0001 nD	±0.0003 nD
		Battery Acid Specific Gravity	1.000 to 1.500	0.001	±0.005

# **Digital Coolant Refractometer**

Item No.	Model No.	Scale	Range	Resolution	Accuracy
48603	MDX-603	Ethylene Glycol (% v/v)	0.0 to 65.0%	0.1%	±0.5%
		Ethylene Glycol Freezing Point	-80.0 to 32.0°F / -50.0 to 0.0°C	0.1°F / 0.1°C	±1.0°F/±0.5°C
		Propylene Glycol (% v/v)	0.0 to 73.0%	0.1%	±0.5%
		Propylene Glycol Freezing Point	-100.0 to 32.0°F / -50.0 to 0.0°C	0.1°F / 0.1°C	±1.0°F/±0.5°C

# C10

# **Abbe Refractometer**

# **High Precision Benchtop Unit**

The C10 Abbe Refractometer is a high precision benchtop instrument measuring the Refractive Index (nD) of liquids, glass, plastic and film.

With an oversized, easy to read scale and Prism Temperature Control feature, the C10 delivers Brix and Refractive Index Scales and dispersion measurements with high precision.

Made of easy to clean, durable metal housing, the C10 is a long-lasting tool perfect for Education, Food, Beverage and Industrial Laboratories.

C10 Abbe Refractometer Specification	ns		
Dispersion (nF-nC)	Dial Readings w/Conversion Charts		
Ambient Temperature	41 to 95°F / 5 to 35°C		
•	32 to 158°F / 0 to 70°C		
Measuring Temperature	0 to 70°C		
Optical Wavelength	589 nm		
Prism	Optical Glass		
Prism Housing	SUS316 Stainless Steel		
Circulator Nozzles	2 for Primary Prism,		
	2 for Secondary Prism;		
	accepts 3/8" ID tubing		
<b>Dimensions</b> (L x W x H in/mm)	5.5 x 3.5 x 9.4 / 140 x 90 x 240		
Weight	9.9 lb / 4.5 kg		
Supplied With	2 x LR44, 1.5V Batteries		
	Glass Standard (1.5163 nD)		
	Bromonaphthalene Contact Liquid		
	(4 mL bottle)		
	Scale Adjustment Tool		
Warranty	1-Year Parts and Labor		



### C10 Abbe Refractometer

Item No.	Scale	Range	Resolution	Accuracy			
C10	Brix	0.0 to 95.0%	0.25% Brix	±0.15% Brix			
	Refractive Index	1.3000 to 1.7200 nD	0.0005 nD	±0.00025 nD			

### Refractometer Standards, Letter of NIST Traceability Certification

		<u> </u>		
Item No.	Туре	Brix Value (68°F / 20°C)	R.I. Value (68°F / 20°C)	Bottle Volume
45000-0S025	Oil	60.00%	1.4420 nD	7.4 mL / 0.25 oz
45000-0S200	Water	0.00%	1.33297 nD	7.4 mL / 0.25 oz
45000-0S001	Contact Solution	0.00%	1.65 nD	7.4 mL / 0.25 oz

# **Manual Polarimeter**

Rugged design and simple operation allows for heavy use at the instructional level

- Polarimeters are used to measure the purity of ingredients in pharmaceuticals, quality control of raw materials, and concentration of ingredients in sugar-based substances
- The model C25L Polarimeter is designed to measure the angle of rotation caused by passing polarized light through an
  optically active substance

### **Manual Polarimeter**

Item No.	Measuring Range	Resolution	Accuracy	Light Source	L x W x H (in)	L x W x H (mm)	Weight	Power
4601P	±180°	0.1°	±0.1°	LED with 589 nm filter	18.1 x 5.5 x 13	460 x 140 x 330	11.5 lbs / 5.2 kg	110 VAC

### Replacement parts

	represent parts				
Item No.	Description				
4600P-SC100	100 mm sample cell				
4600P-SC200	200 mm sample cell				
1200-FS4	Fuse, 1A, 250V				





# **HYDROMETERS**

VEE GEE Hydrometers are manufactured to exacting standards. Each instrument is individually tested and inspected to ensure the highest level of dependability, accuracy, and uniformity. Scales include: API, Alcohol, Battery, Baume, Brix, Calcium Chloride, Liquid Petroleum Gas, Plato, Specific Gravity, Sodium Chloride & Soil Analysis.

# **Hydrometers Explained...**


The hydrometer is an instrument which is constructed on the Archimedes principle that a solid body displaces its own weight of the liquid in which it floats. Hydrometers can be divided into two general classes; namely for liquid heavier than water and for liquids lighter than water. The base hydrometer scale is Specific Gravity, in which distilled water equals 1.000 as the initial point. Liquids lighter than water are scaled below 1.000 specific gravity and liquids heavier than water are scaled above 1.000 specific gravity.

Many other scales are commonly used, such as API, Brix, Baume, Plato, etc. All of which are convertible into specific gravity by formula.

Hydrometers are usually calibrated at 15°C/60°F. To determine the density of a liquid, the liquid should be at 15°C/60°F. If the temperature varies, the liquid will either contract or expand, depending upon the temperature. Therefore, the density fluctuates with the temperature. Where there is a variation from the standard 15°C/60°F, corrections must be applied to the hydrometer reading. To assure proper corrections, a separate accurate thermometer should be used, or a hydrometer in combination with a thermometer, which is sometimes referred to as a "thermohydrometer."

# The correct method of reading a hydrometer follows:

- A. Observe a point below the plane of the liquid surface. The surface should appear as an ellipse (Fig. 1).
- B. The line of vision is raised until the surface, seen first as an ellipse, becomes a straight line (Fig. 2).
- C. The point at which this line cuts the hydrometer scale is the reading of the instrument (Fig. 3).

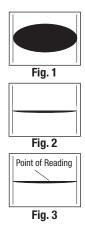



If the liquid is not sufficiently clear for readings to be made in this manner, read from above the surface and estimate as accurately as possible the point to which the liquid rises on the hydrometer stem. Since hydrometers are calibrated to give correct indications when read at the principal surface of the liquid, correct the reading just taken at the upper edge of the meniscus by an amount equal to this height above the principal surface of the liquid. The amount of correction can be determined with sufficient accuracy for most purposes by taking a few readings on the upper and lower meniscus in a clear liquid of the same character as that being tested and noting the differences.

# **Alcohol Hydrometers - Tralle and Proof Scales**

Verify the alcohol content of distilled spirits

Both models include the same ranges; model 6612-2TS also features an internal thermometer with SafetyBLUE fill for measuring solution temperature from -10 to  $+100^{\circ}$ F.

- Calibrated for use at 15°C/60°F

NOTE: Not for use with beer or wine.

Alcohol Hydrometers - Tralle and Proof Scales

Item No.	Range	Subdivisions/ Tolerance	Range	Subdivisions/ Tolerance	Approx. Length (in/mm)
		TRALLE		PR00F	
6612-1	0 to 100%	1% / ±1%	0 to 200	2 / ±2%	12 / 305
6612-2TS	0 to 100%	1% / ±1%	0 to 200	2 / ±2%	14 / 355

6612-2TS

# **Alcohol Hydrometers - Proof Scale, Internal Revenue Service Spec.**

Verify the alcohol content of your ethyl alcohol-based spirits with these hydrometers

- Tolerances match those laid out by the Internal Revenue Service specifications
- Calibrated for use at 15°C/60°F

Alcohol Hydrometers - IRS Specification

Item No.	Proof Range	IRS Size	Subdivisions	Tolerance	Approx. Length (in/mm)
6613-B	80 to 120	В	0.5	±0.5%	9 / 230
6613-F	0 to 20	F	0.2	±0.3%	12 / 305
6613-G	20 to 40	G	0.2	±0.4%	12 / 305
6613-H	40 to 60	Н	0.2	±0.4%	12 / 305
6613-I	60 to 80	I	0.2	±0.3%	12 / 305
6613-K	75 to 95	K	0.2	±0.3%	12 / 305
6613-L	90 to 110	L	0.2	±0.3%	12 / 305
6613-M	105 to 125	M	0.2	±0.3%	12 / 305
6613-N	125 to 145	N	0.2	±0.3%	12 / 305
6613-P	145 to 165	Р	0.2	±0.3%	12 / 305
6613-Q	165 to 185	Q	0.2	±0.2%	12 / 305
6613-R	185 to 206	R	0.2	±0.2%	12 / 305

Alcohol Hydrometers w/3-point NIST Certificate and PVC Case- Proof Scale, Internal Revenue Service Specification

Verify the alcohol content of your ethyl alcohol-based spirits with these hydrometers Tolerances match those laid out by the Internal Revenue Service specifications.

- Commonly used in distilleries requiring certification
- Calibrated for use at 15°C/60°F

 $\textbf{INCLUDES:} \ \textbf{Hydrometer, NIST-traceable certificate with 3 data points, and a protective PVC \ case$ 

Alcohol Hydrometers w/ NIST Cert. IRS Specs. & Case

Item No.	Proof Range	IRS Size	Subdivisions	Tolerance	Approx. Length (in/mm)
6613-B-C	80 to 120	В	0.5	±0.5%	9 / 230
6613-F-C	0 to 20	F	0.2	±0.3%	12 / 305
6613-G-C	20 to 40	G	0.2	±0.4%	12 / 305
6613-H-C	40 to 60	Н	0.2	±0.4%	12 / 305
6613-I-C	60 to 80	I	0.2	±0.3%	12 / 305
6613-K-C	75 to 95	K	0.2	±0.3%	12 / 305
6613-L-C	90 to 110	L	0.2	±0.3%	12 / 305
6613-M-C	105 to 125	M	0.2	±0.3%	12 / 305
6613-N-C	125 to 145	N	0.2	±0.3%	12 / 305
6613-P-C	145 to 165	Р	0.2	±0.3%	12 / 305
6613-Q-C	165 to 185	Q	0.2	±0.2%	12 / 305
6613-R-C	185 to 206	R	0.2	±0.2%	12 / 305





6613-F

# **Plato Hydrometers with Thermometers**

Calibrated to show the concentration of extract as a percentage by weight in brewery wort

For example, a 5°P reading indicates 5 g of extract per 100 g of wort.

- Primarily used in the brewing industry
- Calibrated for use at 68°F/20°C

# **Plato Hydrometers with Thermometers**

Item No.	Plato Range	Subdivisions	Approx. Length (in/mm)	Thermometer Range	Thermometer Subdivision	Thermometer Liquid Fill
6614TS-5	0 to 8.5°	0.1°	14 / 360	0 to 50°C	1°C	SafetyBLUE
6614TS-6	7.5 to 16°	0.1°	14 / 360	0 to 50°C	1°C	SafetyBLUE
6614TS-7	15.5 to 24°	0.1°	14 / 360	0 to 50°C	1°C	SafetyBLUE
6614TS-8	0 to 32°	0.5°	14 / 380	30 to 120°F	2°F	SafetyBLUE

6614TS-5

# **API ASTM Hydrometers with Thermometer**

**Chosen as the Standard Scale for Petroleum** 

The API Scale was selected as the standard hydrometer scale for petroleum products by the American Petroleum Institute (API), the U.S. Bureau of Mines, and the National Institute of Standards and Technology (NIST).

- The API Scale is based on the following formula:  $^{\circ}$ API @  $60^{\circ}$ F = (141.5/(Specific Gravity at  $60^{\circ}$ F) 131.5
- The thermometer scale is located in the body of the hydrometer
- Calibrated for use at 15°C/60°F

# **API ASTM Hydrometers with Thermometer**

ALIAOIM	i i y ai o i i o toi o	WILLI THOTHION	10101			ALI TAO IN TITUTO NOTO WITH THORNOLOG								
Item No.	ASTM No.	API Range	Subdivisions	Approx. Length (in/mm)	Thermometer Range	Thermometer Subdivision	Thermometer Liquid Fill							
6752HTS	52H	9 to 21°	0.1°	15 / 380										
6753HTS	53H	19 to 31°	0.1°	15 / 380										
6754HTS	54H	29 to 41°	0.1°	15 / 380	0 to 15005	005	CofobiDLUE							
6755HTS	55H	39 to 51°	0.1°	15 / 380	- 0 to 150°F	2°F	SafetyBLUE							
6756HTS	56H	49 to 61°	0.1°	15 / 380										
6757HTS	57H	59 to 71°	0.1°	15 / 380										

6756HTS

# **Battery Hydrometer - Syphon Set**

For Measuring Battery Acid

The Battery Acid Hydrometer Set is used in many industrial settings including, fleet maintenance, automotive services, solar-powered system maintenance, and general industrial applications.

- The rubber parts are made of nitrile rubber, which is chemically inert and durable
- A reinforced box provides storage and is compartmented to eliminate breakage
- Calibrated for use at 25°C/77°F

INCLUDES: Specific Gravity Hydrometer, glass syphon tube with rubber bulb, nozzle, and clamp

### **Battery Hydrometers**

, ,				
Item No.	Description	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)
6605-5	Hydrometer and Syphon Set	1.150 to 1.300	0.005	20 / 507
6605-5H	Hydrometer Only	1.150 to 1.300	0.005	8.7 / 220
6605-0SY	Glass syphon tube, rubber bulb, nozzle, and clamp only			20 / 507

6605-5

# **Baume (Heavy) Hydrometers**

These Baume Hydrometers are for liquids heavier than water

- The Baume scale is based on the following formula: Specific Gravity = 145/(145-degrees Baume)
- Calibrated for use at 15°C/60°F

### Baume (Heavy) Hydrometers

Item No.	Baume Range	Subdivisions	Length (in/mm)
6609-1	0 to 12°	0.1°	
6609-2	9 to 21°	0.1°	
6609-3	19 to 31°	0.1°	
6609-4	29 to 41°	0.1°	12 / 305
6609-5	39 to 51°	0.1°	
6609-8	0 to 15°	0.1°	
6609-9	0 to 25°	0.2°	

# Baume (Heavy) Hydrometers

Daume (meavy) myurumeters						
Item No.	Baume Range	Subdivisions	Length (in/mm)			
6609-10	0 to 35°	0.5°				
6609-11	0 to 50°	0.5°				
6609-12	35 to 70°	0.5°	12 / 305			
6609-13	0 to 50°	1.0°				
6609-14	0 to 70°	1.0°				
6609-15	0 to 90°	1.0°				

6609-1







# **Brix Hydrometers**

Brix Hydrometers are calibrated to show the percentage of sucrose by weight at 20°C/68°F

The Brix scale is based on the following: 1°Brix = 1% sucrose by weight at the specified temperature (20°C/68°F)

### **Brix Hydrometers**

### **Brix Hydrometers**

Item No.	Brix Range	Subdivisions	Approx. Length (in/mm)	Item No.	Brix Range	Subdivisions	Approx. Length (in/mm)
6601-1	0 to 12°	0.1°	13 / 330	6601-10	0 to 35°	0.5°	13 / 330
6601-2	9 to 21°	0.1°	13 / 330	6601-11	35 to 70°	0.5°	13 / 330
6601-3	19 to 31°	0.1°	13 / 330	6601-13	-5 to 5°	0.1°	13 / 330
6601-4	29 to 41°	0.1°	13 / 330	6601-14	5 to 15°	0.1°	13 / 330
6601-5	39 to 51°	0.1°	13 / 330	6601-15	15 to 25°	0.1°	13 / 330

6601-15

# **Brix Hydrometers with °C Thermometers**

Brix Hydrometers are calibrated to show the percentage of sucrose by weight at 20°C/68°F

The Brix scale is based on the following: 1°Brix = 1% sucrose by weight at the specified temperature (20°C/68°F)

- The thermometer includes a correction scale printed opposite the temperature scale for conversion of the reading to 20°C/68°F

### Brix Hydrometers with °C Thermometers

Item No.	Brix Range	Subdivisions	Approx. Length (in/mm)	Thermometer Range	Thermometer Subdivision	Thermometer Liquid Fill
6601TS-1	0 to 12°	0.1°	14.5 / 370			SafetyBLUE
6601TS-2	9 to 21°	0.1°	14.5 / 370	_		
6601TS-3	19 to 31°	0.1°	14.5 / 370	_		
6601TS-4	29 to 41°	0.1°	14.5 / 370			
6601TS-5	39 to 51°	0.1°	14.5 / 370			
6601TS-6	49 to 61°	0.1°	14.5 / 370	0 to 50°C	100	
6601TS-8	69 to 81°	0.1°	14.5 / 370	- 0 to 50°C	1°C	
6601TS-10	0 to 35°	0.5°	14.5 / 370	_		
6601TS-12	0 to 70°	1.0°	14.5 / 370			
6601TS-13	-5 to 5°	0.1°	14.5 / 370			
6601TS-14	5 to 15°	0.1°	14.5 / 370			
6601TS-15	15 to 25°	0.1°	14.5 / 370			

6601TS-2

# **Brix Hydrometers with °F Thermometers**

Brix Hydrometers are calibrated to show the percentage of sucrose by weight at 20°C/68°F

The Brix scale is based on the following: 1°Brix = 1% sucrose by weight at the specified temperature (20°C/68°F)

- The thermometer includes a correction scale printed opposite the temperature scale for conversion of the reading to 20°C/68°F

### **Brix Hydrometers With °F Thermometers**

Item No.	Brix Range	Subdivisions	Approx. Length (in/mm)	Thermometer Range	Thermometer Subdivision	Thermometer Liquid Fill		
6601TS-1F	0 to 12°	0.1°	14.5 / 370					
6601TS-2F	9 to 21°	0.1°	14.5 / 370	_				
6601TS-3F	19 to 31°	0.1°	14.5 / 370	30 to 130°F	2°F	SafetyBLUE		
6601TS-10F	0 to 35°	0.5°	14.5 / 370	_				
6601TS-13F	-5 to 5°	0.1°	14.5 / 370					

6601TS-1F



# Calcium Chloride (CaCl2) Specific Gravity/Freezing Point Hydrometer

This hydrometer measures the specific gravity and freezing point of calcium chloride in solution

- Calibrated for use at 15°C/60°F

Calcium Chloride (CaCl2) Specific Gravity/Freezing Point Hydrometer

Item No.	Range	Subdivisions	Approx. Length (in/mm)
6611.0	1.000 to 1.280 Specific Gravity	0.002	12 / 305
6611-3	+30 to -40°F Freezing Point	5°	12 / 305

# **Sodium Chloride Hydrometers**

Measure the percentage of sodium chloride in solution with these hydrometers

Model 6611-1 measures the % saturation of sodium chloride in solution; Model 6611-2 measures the % by weight of sodium chloride in solution.

- Calibrated for use at 15°C/60°F

# **Sodium Chloride Hydrometers**

Item No.	Range	Subdivisions	Approx. Length (in/mm)	
6611-1	0 to 100% saturation	1.0%	12 / 305	
6611-2	0 to 26.5% by weight	0.5%	12 / 305	



**Specific Gravity,** also known as relative density, is the ratio of the mass of a liquid to the mass of an equal volume of distilled water.

# **Ultra-Precision Specific Gravity Hydrometers**

Featuring smaller ranges and higher resolution, with subdivisions of 0.0005

- Calibrated for use at 15°C/60°F

### Ultra-Precision Specific Gravity Hydrometers

olita-i recision opecine dravity riyurometers							
Item No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)				
6602-2	0.700 to 0.770	0.0005	13 / 330				
6602-3	0.760 to 0.830	0.0005	13 / 330				
6602-4	0.820 to 0.890	0.0005	13 / 330				
6602-5	0.880 to 0.950	0.0005	13 / 330				
6602-6	0.940 to 1.010	0.0005	13 / 330				
6602-7	1.000 to 1.070	0.0005	13 / 330				
6602-8	1.060 to 1.130	0.0005	13 / 330				
6602-9	1.120 to 1.190	0.0005	13 / 330				
6602-10	1.180 to 1.250	0.0005	13 / 330				
6602-11	1.240 to 1.310	0.0005	13 / 330				
6602-12	1.300 to 1.370	0.0005	13 / 330				
6602-13	1.360 to 1.430	0.0005	13 / 330				
6602-14	1.420 to 1.490	0.0005	13 / 330				
6602-15	1.480 to 1.550	0.0005	13 / 330				
6602-19	1.720 to 1.790	0.0005	13 / 330				



# **Precision Specific Gravity Hydrometers**

Precision Specific Gravity Hydrometers feature subdivisions from 0.001 to 0.01

- Calibrated for use at 15°C/60°F

# **Precision Specific Gravity Hydrometers**

Item No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)
6603-10	0.700 to 0.810	0.001	12 / 305
6603-11	0.800 to 0.910	0.001	12 / 305
6603-12	0.900 to 1.000	0.001	12 / 305
6603-1	1.000 to 1.220	0.002	12 / 305
6603-2	1.200 to 1.420	0.002	12 / 305
6603-3	1.400 to 1.620	0.002	12 / 305
6603-4	1.600 to 1.820	0.002	12 / 305
6603-5	1.800 to 2.020	0.002	12 / 305
6603-21	2.200 to 2.420	0.002	12 / 305
6603-22	2.400 to 2.620	0.002	12 / 305
6603-13	0.650 to 1.000	0.005	12 / 305
6603-6	1.000 to 1.600	0.005	12 / 305
6603-7	1.000 to 2.000	0.010	12 / 305
6603-8	2.000 to 3.000	0.010	12 / 305
6603-9	3.000 to 4.000	0.010	12 / 305



# **Precision Specific Gravity Hydrometers with Thermometer**

No need for a separate thermometer to check liquid temperature

- These Hydrometers feature an internal thermometer for calculating temperature corrections
- Calibrated for use at 15°C/60°F

### Precision Specific Gravity Hydrometers with Thermometer

Item No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)	Thermometer Range	Thermometer Subdivision	Thermometer Liquid Fill
6602TS-7	1.000 to 1.070	0.001	15 / 385			
6602TS-8	1.060 to 1.130	0.001	15 / 385	20 to 150°F	2°F	SafetyBLUE
6603TS-1	1.000 to 1.220	0.002	12 / 305			

6602TS-7

# **Precision Short-Form Specific Gravity Hydrometers**

Short-Form Hydrometers use less volume of the test liquid, while still offering 0.001 resolution

- Calibrated for use at 15°C/60°F

# **Precision Short-Form Specific Gravity Hydrometers**

Item No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)
6602-4S	0.820 to 0.890	0.001	6.5 / 165
6602-5S	0.880 to 0.950	0.001	6.5 / 165
6602-6S	0.940 to 1.010	0.001	6.5 / 165
6602-7S	1.000 to 1.070	0.001	6.5 / 165
6602-8S	1.060 to 1.130	0.001	6.5 / 165
6602-9S	1.120 to 1.190	0.001	6.5 / 165
6602-10S	1.180 to 1.250	0.001	6.5 / 165

6602-7S



# **ASTM Specific Gravity Hydrometers**

Manufactured to the specifications of the American Society for Testing and Materials (ASTM)

Select the hydrometer that you need by matching the ASTM number shown in the table.

- Calibrated for use at 15°C/60°F

# **ASTM Specific Gravity Hydrometers**

# **ASTM Specific Gravity Hydrometers**

	No.	Specific Gravity Range	Subdivi- sions	Approx. Length (in/mm)	Item No.	ASTM No.	Specific Gravity Range	Subdivisions	Approx (in/mm)
67111H	111H	1.000 to 1.050	0.0005	13 / 330	67125H	125H	1.000 to 1.050	0.001	12 / 26
67112H	112H	1.050 to 1.100	0.0005	13 / 330	67126H	126H	1.050 to 1.100	0.001	12 / 260
67113H	113H	1.100 to 1.150	0.0005	13 / 330	67127H	127H	1.100 to 1.150	0.001	12 / 260
67114H	114H	1.150 to 1.200	0.0005	13 / 330	67128H	128H	1.150 to 1.200	0.001	12 / 260
67115H	115H	1.200 to 1.250	0.0005	13 / 330	67129H	129H	1.200 to 1.250	0.001	12 / 260
67116H	116H	1.250 to 1.300	0.0005	13 / 330	67130H	130H	1.250 to 1.300	0.001	12 / 260
67117H	117H	1.300 to 1.350	0.0005	13 / 330	67131H	131H	1.300 to 1.350	0.001	12 / 260
67118H	118H	1.350 to 1.400	0.0005	13 / 330	67132H	132H	1.350 to 1.400	0.001	12 / 260
67119H	119H	1.400 to 1.450	0.0005	13 / 330	67133H	133H	1.400 to 1.450	0.001	12 / 260
67120H	120H	1.450 to 1.500	0.0005	13 / 330	67134H	134H	1.450 to 1.500	0.001	12 / 260

# 67125H

# **Wide-Range Specific Gravity Hydrometers**

Ideal when you aren't sure of the approximate value of your liquid

- Also use for multiple liquids with different specific gravity within the range of the hydrometer, allowing multiple uses with one unit
- Calibrated for use at 15°C/60°F

# Wide-Range Specific Gravity Hydrometers

Item No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)
6604-1	1.000 to 1.250	0.005	6.5 / 165
6604-2	1.200 to 1.450	0.005	6.5 / 165
6604-3	1.400 to 1.650	0.005	6.5 / 165
6604-4	1.600 to 1.850	0.005	6.5 / 165
6604-5	1.800 to 2.050	0.005	6.5 / 165

6604-1



# **Universal Specific Gravity Hydrometer**

Ideal when you aren't sure of the approximate value of your liquid

This dual-scale model measures specific gravity and Baume. Also use for multiple liquids with different readings, allowing multiple uses with one hydrometer.

- Calibrated for use at 15°C/60°F

# **Universal Specific Gravity Hydrometer**

Item No.	Specific Gravity Range	Subdivisions	Baume Range	Subdivisions	Approx. Length (in/mm)
6605-1	0.700 to 2.000	0.005 (from 0.7 to 1.0); 0.01 (from 1.0 to 2.0)	0 to 72° (Heavy); 70 to 10° (Light)	1°	15 / 380

6605-1

# **Dual Scale Specific Gravity / Baume Hydrometers**

Dual-Scale Hydrometers have two scales on the stem, covering different applications in one unit

The Baume scale has two different scales: one for liquids heavier than water and one for liquids lighter than water. See table below for Baume ranges.

- Calibrated for use at 15°C/60°F

# **Dual Scale Specific Gravity / Baume Hydrometers**

Dual Could Opcome dravity / Dualite Hydromotors						
Item No.	Specific Gravity Range	Baume Range	Subdivisions	Approx. Length (in/mm)		
6603DS-1	1.000 to 1.220	0 to 26° (Heavy)	0.002 / 0.2°	12 / 305		
6603DS-1S	1.000 to 1.225	0 to 26° (Heavy)	0.005 / 1°	12 / 305		
6603DS-2	1.200 to 1.420	24 to 42° (Heavy)	0.002 / 0.2°	12 / 305		
6603DS-2S	1.200 to 1.425	24 to 42° (Heavy)	0.005 / 1°	12 / 305		
6603DS-3	1.400 to 1.620	42 to 55° (Heavy)	0.002 / 0.2°	12 / 305		
6603DS-4	1.600 to 1.820	55 to 65° (Heavy)	0.002 / 0.2°	12 / 305		
6603DS-5	1.800 to 2.020	64 to 72° (Heavy)	0.002 / 0.2°	12 / 305		
6603DS-6	1.000 to 1.450	0 to 45° (Heavy)	0.005 / 0.5°	12 / 305		
6603DS-6S	1.000 to 1.400	0 to 41° (Heavy)	0.010 / 1°	12 / 305		
6603DS-7	1.000 to 2.000	0 to 72° (Heavy)	0.010 / 1°	12 / 305		
6603DS-13	0.600 to 1.000	10 to 100° (Light)	0.005 / 1°	12 / 305		
6603DS-04	0.70 to 1.00	10 to 70° (Light)	0.01 / 1°	6.5 / 165		

6603DS-1



# Liquid Petroleum Gas (LPG) Specific Gravity Hydrometer with Thermometer

LPG Hydrometers measure the specific gravity of LPG such as propane, butane, and other light hydrocarbons includes an internal thermometer for temperature measurement.

- Calibrated for use at 15°C/60°F

# Liquid Petroleum Gas (LPG) Specific Gravity Hydrometer with Thermometer

Item No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)	Thermometer Range	Thermometer Subdivision	Thermometer Liquid Fill
67101HTS	0.500 to 0.650	0.001	14.4 / 365	30 to 90°F	1°F	SafetyBLUE

67101HTS

# **Soil ASTM Specific Gravity Hydrometer**

Measure the specific gravity of soil samples in suspension

Often used with a sedimentation cylinder for soil analysis.

- Calibrated for use at 20°C/68°F

# Soil ASTM Specific Gravity Hydrometer

Item No.	ASTM No.	Specific Gravity Range	Subdivisions	Approx. Length (in/mm)
67151H	151H	0.995 to 1.038	0.001	11 / 280





# **Soil ASTM Soil Colloids Hydrometer**

Measure the amount of soil colloids in grams of soil samples in suspension Often used with a sedimentation cylinder for soil analysis.

- Calibrated for use at 20°C/68°F

# Soil ASTM Soil Colloids Hydrometer

Item No.	ASTM No.	Soil Colloids Range	Subdivisions	Approx. Length (in/mm)
67152H	152H	-5 to 60g	1.0	11 / 280





### **Hydrometer Cases**

**PVC Hydrometer Cases help protect your hydrometer from damage** 

Ideal for storing your hydrometer when not in use, these cases feature a foam or felt padded insert and padded bottom and cap.

### **Hydrometer Cases**

Item No.	Description	Approx. Length (in/mm)
66CS-8	PVC Case	8.5 / 216
66CS-8F	PVC Case, short	8.5 / 216
66CS-13	PVC Case	13.5 / 343
66CS-13J	PVC Case, wide-body	13.5 / 343
66CS-16	PVC Case	16.5 / 420



Made of borosilicate glass, these jars feature a beaded rim and a circular base

Sized properly to eliminate surface tension effects that lead to inaccurate readings.

- Model 21050-1000-L and 21000-1000-L have a graduation line at 1000 mL for precise filling and are used as a sedimentation cylinder

**Hydrometer Jars** 

Item No.	Description	Maximum Hydrometer Length (mm)
21050-100	100 mL Hydrometer Jar	185
21050-250	250 mL Hydrometer Jar	280
21050-500	500 mL Hydrometer Jar	360
21050-600	600 mL Hydrometer Jar	390
21050-1000	1000 mL Hydrometer Jar	400
21050-1000-L	1000 mL Hydrometer Jar with Graduation Line	400

21050-1000-L



Hydrometer Series	Description	Approx. Length (mm)	Uses Case	Uses Jar
6601-Series	Brix	330	66CS-13	21050-500
6601TS-Series	Brix w/Thermometer	370	66CS-16	21050-600 or -1000
6602, 6603-Series	Specific Gravity	305	66CS-13	21050-500
6603DS-Series	Specific Gravity & Baume	305	66CS-13	21050-500
6604-Series	Specific Gravity	165	66CS-8F	21050-100
6605-1	Specific Gravity & Baume: Universal	380	66CS-16	21050-600 or -1000
6609-Series	Baume: Heavy	305	66CS-13	21050-500
6611-Series	Salt	305	66CS-13	21050-500
6612-1	Alcohol: Tralle and Proof	305	66CS-13	21050-500
6612-2TS	Alcohol: Tralle and Proof w/Thermometer	355	66CS-16	21050-600 or -1000
6613-Series	Alcohol: Proof, IRS	230-305	66CS-13J	21050-250 or -500
6614TS-Series	Plato w/Thermometer	360-380	66CS-16	21050-600 or -1000
675xHTS-Series	API, ASTM w/Thermometer	380	66CS-16	21050-600 or -1000
67111H-67134H	Specific Gravity: ASTM	260-330	66CS-13	21050-250 or -500
67101HTS	LPG w/Thermometer	365	66CS-16	21050-600 or -1000
6715xH	Specific Gravity or Soil Colloids: ASTM	280	66CS-13J	21050-250 or -500







### **RF3000®**

Packed with features for comfortable and accurate pipetting

The powerful pump offers rapid, yet sensitive aspiration and dispensing, has adjustable valves and is microprocessor controlled for constant reliability and accuracy

- Long-life li-ion battery, with quick charge-up to 20 hours use
- Toggle operation by pressing the large mode button
- An adapter can be fitted to the nose cone to stabilize small volume and thinner pipettes
- New USB connecter ensures the appropriate charger is used with the unit
- Low-battery indicator
- Can be recharged when in use



Nose cone assembly only



Storage Stand doubles as a charging station

Item No.	volume Range mL	L X W X H in	L X W X H cm
HS3000	1–100	9 x 8.5 x 5	22.9 x 21.6 x 12.7
Replacement p	arts		
HS5000RFFIL	RF3000® Replacement Filters, <i>Pack of 5</i>		
120669	RF3000® Universal Power Adaptor 5VDC Micro USB Plug		
HS300502	RF3000® Li-lon Rechargeable Battery, One year warranty on battery included		
120487	RF3000® Nose Cones Assorted* ■ ■ ■ Set of	of 5	
HS5000RFRFA	RF3000® Replacement Filter Adapter		
HS5000RFRP	RF3000® Replacement Pipette Silicone Adapter		

\*Assorted colors (Black, Blue, Grey, Red and Yellow)

### **Easy Identification**

5 different colored nose cones help differentiate department, task or user; as well as a labeling area on the side of the unit



Magnetic & Screw Wall-Mount



### 3 Modes

Bright LEDs indicate 3 operation modes: high, low, & gravity



Specifications	
Power	Micro-USB connector Rechargeable li-lon battery (replaceable)
Filter	PTFE/25mm dia/0.45 µm filter
Pipettes	1-100 mL, glass and plastic
Pipette Nozzle	Polypropylene









Package	Includes	

RF3000® Pipette Controller

Charging Stand and Wall-Mountable/Magnetic Holder

Nozzle Adapter for use with Small Pipettes

Li-Ion Rechargeable Battery

Five Nose Cones: Black, Yellow, Red, Grey and Blue

5VDC Micro USB Plug

6 ft (183 cm) Cable

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs









### **Pipette Controller**

### Precise fluid control at your finger tips

Comfortable to hold and easy to use, the manual inline pipette controller was designed with an ergonomic grip and bulb for simple and precise operation

A T-bar thumb lever controls the aspirating and dispensing speeds from serological pipettes from 1 to 100 mL in size

- Conveniently accommodates either glass or plastic serological or volumetric pipettes ranging from 1-100 mL
- Includes a replaceable .45 µm hydrophobic filter to protect your samples and internal components



Nose cone assembly only

Item No.	Color	Volume Range mL	L x W x H in	L x W x H
120479	Blue	1–100	9 x 2.1 x 2.1	23 x 5.4 x 5.4

### Replacement Parts

Item No.	Description	Color
120604	Three Noses Cone with Adapter	Green, White, Blue 🔳 🗌 🔳

### **Operation Is Simple:**

Squeeze the bulb

to create a vacuum to draw up fluid

Push the control lever backward (up)

to aspirate liquid through a pipette

Push the control lever forward (down) to dispense the liquid

Press the central bulb on the unit

to blow out any remaining fluid if necessary





### **RF1000**

Lightweight and easy to use

The RF1000 manual pipette controller is uniquely designed for one-handed operation

Simply squeeze the main bulb and push the toggle forward to aspirate or back to dispense. Push the blowout bulb to expel any liquid remaining in the pipette

- Tapered silicone rubber pipette adapter allows a large range of pipette sizes ranging 1-100 mL and is replaceable
- Includes a replaceable .45  $\mu m$  hydrophobic filter to protect your samples and internal components




Nose cone assembly only

Item No.	Color	Volume Range mL	L x W x H in	L x W x H cm
HS1000RFU	Blue	1–100	7.1 x 1.8 x 4.7	18 x 4.5 x 12

### Replacement Parts

Item No.	Description Color		
HS5000RFFIL	. Replacement Filters, <i>Pack of 5</i>		
120604	Three Noses Cones with Adapter Blue		

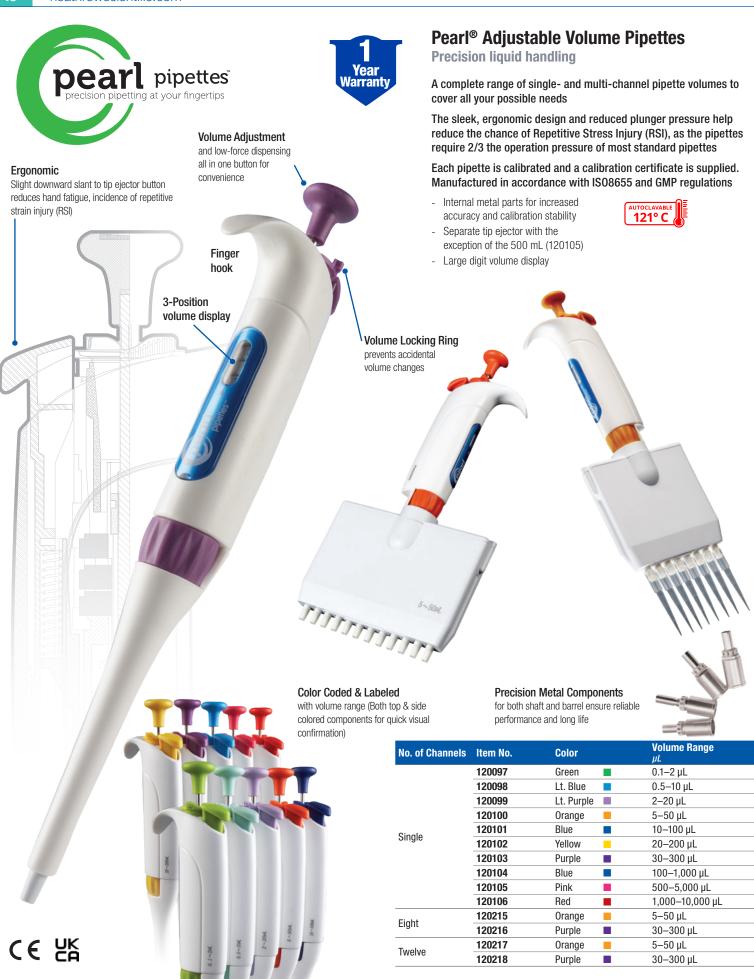



### Pipette Starter Kit Item No. 120554

**Universal Carousel Pipette Stand** For more information see page 78

Upper Holder Color	Lower Holder Color	Base Width in	Top W x H in	Base Width cm	Top W x H cm
Assorted/Clear	Clear	6.2	9 x 16.1	15.8	22.9 x 40.8

\*Assorted colors (Purple, Blue, Green, Grey, Red Orange, Yellow and Clear)

Pearl® Pipettes - Single Channel For more information see page 40

Color	Volume Range	
	μL	
Lt. Purple	2-20 μL	
Yellow	20-200 μL	
Blue	100-1,000 μL	

Reagent Reservoirs For more information see page 84

Color	<b>Volume</b> mL	Material	Sterility	L x W x H in	L x W x H cm
White	50	Polystyrene	Sterile	4 x 1.1 x 1.2	10.2 x 2.8 x 3 Pack of 5





### Rapid, accurate colony counter

Combines a multi-function electronic counter with a Sharpie® Pen for marking to prevent missing or double-counting colonies

- Light, ergonomic body is balanced for easy manipulation
- Select between an audible beep or an LED flash to acknowledge count
- Store up to 32 separate counts and display the total of all counts
- Count up/down verification
- Quick and easy to change between different colored markers to check different colony counts in the same view

### Package Includes:

- eCount Colony Counter
- Black Sharpie® Pen (Page 48)
- Long-life button battery included
- Hand tally Adapter
- Weighted Stand





US Patent No. D564,385 EP Registered Design No. 705751-0001

		Er riogiotoroa Booigir No. 1 00101 0001
Item No.	LxW	LxW
	in	ст
120000	7.4 x 1.2*	17 x 3*
120228	CR1225 Battery	

\*Measurements shown are with Sharpie  $^{\footnotesize @}$  Pen inserted



eCount with Mini-Light Box with Mini-Magnifier, page 43



Weighted Stand for bench top use







### **Mini-Light Box and Mini-Magnifier**

Perfect accessories for the eCount or use as stand-alone items

### Mini-Light Box



- LED Light source. Color temperature 5,300K +/- 270K
- Viewing area (125 x100 mm) is perfect for quick screening of 96-well plates, petri dishes, electrophoresis gels, slides, and transparencies
- Powered by six 1.5V AAA batteries (included), not included in warranty

### Mini-Magnifier

- Swiveling molded lens 1.75x magnification
- Enlarges objects for easier viewing

Item No.	Description	L x W x H in	L x W x H cm
120661	Mini-Light Box*	6.2 x 5.9 x 0.6	15.8 x 15 x 1.4
HS15092D	Mini-Magnifier	6.3 x 5.3 x 8.1	16.1 x 13.5 x 20.7

\*Batteries included



### **Hand Tally Counter**

Use for counting blood cells, bacterial colonies, drops of liquid or any repetitive event

- Tallies up to 9999
- Features a quick-reset knob
- Lens for easy reading, and finger ring for a secure grip

Item No.	Dia. x W	Dia. x W	
	in	ст	
HS6594	1.8 x 1.6	4.6 x 4.1	



### **Benchtop Counters**

Mechanical, benchtop counters are perfect for blood cells (label included), bacterial colonies or repetitive counting

- Each key records up to 999
- Blood cell key labels
- Total window keeps track of total number of strokes made on all other keys; bell rings when total reaches 100

Item No.	No. of Keys	L x W x H in	L x W x H cm
HS5360A	6-Key	9.8 x 3 x 2.2	25 x 7.6 x 5.6
HS5360B	9-Key	12.5 x 3 x 2.2	32 x 7.6 x 5.6

### **Workstation Essentials**

Products necessary for workflow efficiency and keeping your work space clean and organized



Lab Mats
Notebooks & Timers
Drying Racks
Organization

### Curiosity drives us beyond yesterday's way of thinking:

We are an inquisitive bunch focused on understanding the lab science world and how to make a difference.

We change, we adapt and we expect creative thinking and solutions from everyone we align with.

We don't just make lab products, we make the work lab professionals do easier, more efficient and more fun.





### **Workstation Lab Mat**

Safely facilitate workflow processes

- Protect counter tops from hot items up to 200°C (392°F)
- Lip around the edge can hold approx. 250 mL of liquid
- Small footprint, 12 x 14.9 x 0.1 in (30.5 x 37.9 x 0.2 cm)





### **Gridded Surface**

for tube staging, preventing movement of objects

### Reversible

Side One: combination gridded and smooth surface Side Two: smooth surface only





Item No.	Color	LxWxH in	L x W x H cm
120747	Blue	12 x 14.9 x 0.1	30.5 x 37.9 x 0.2
120748	Purple	12 X 14.9 X U.1	30.3 X 37.9 X 0.2

Large Working Surface accommodates benchtop equipment,

weigh boats, stir bars, scales, etc.

### Silicone Lab Mat

Large working space

The lab mat can be used with items such as centrifuges, vortex mixers, hot plates, stirrers, pipettes and so much more

- Mat features designated magnetic area to keep stir bars in place
- Protect countertops from hot items up to 200°C (392°F)



### Raised Edges

enable liquid retention of approximately 300 mL of fluid







Item No.	Side #1 Color	Side #2 Color	L x W x H in	L x W x H cm
120506	Yellow	Grey/Blue		
120507	Purple	Purple/Grey	14 x 23.5 x 0.1	35 x 60 x 0.2
120599	Blue	White/Blue		

### **Laboratory Notebooks**

Keep information safe and well organized with Laboratory Notebooks

Record research data, experimental procedures, design and development ideas and more

- Quality Smyth Sewn binding holds pages securely and lies flat when opened
- Made of acid-free paper, pages are numbered and include space for title, date, project number, book number, date and witness signature
- User data, notebook instructions, table of contents, periodic table, and conversion tables are included in every book




Keep Information Safe with a durable hardbound.

with a durable nardbound water-resistant cover

Case of 12

Marker Ribbon to keep your place

					0400 01 12
Item No.	Color	Pages	Page Format.	L x W in	L x W cm
HS8610A	Black	100	Grid		
HS8610B	Blue	200	dilu	0.6 v 11.1	21.8 x 28.1
HS8610D	Red	200	Lined	0.0 X 11.1	21.0 X 20.1
HS8610E	Black	100	Lilleu		

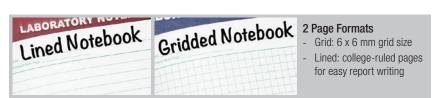





### **Twin Tip & Fine Tip Sharpie® Permanent Pens**

Ink dries quickly and resists both fading and water

- Pen has a fine and an ultra fine tip and comes in a bright selection of colors
- Ink dries quickly and resists both fading and water
- AP certified

Item No.	Description	Color	
HS15094	Twin Tip Marker Set, Pack of 8	Assorted	
HS15092ZP	Fine Point Marker, Pack of 12	Black	

\*Assorted colors (Black, Turquoise, Navy, Lime, Orange, Purple, Sky Blue, Berry) Sharpie® is a registered trademark of Sanford, A Newell Rubbermaid Company

### Lab Alert® Clocks & Timers





### Pocket Timer/ Stopwatch/Clock

- Four independent timing channels with visual and audible alarms
- Use as a stopwatch, a count up/down timer, or a clock
- Includes spring clip and magnetic back





### **Triple Timer/Clock**

- Three channels to count up/down with
- Simultaneous display, time-out, and 60 second audible alarm
- Channels may be started independently or simultaneously
- Timer automatically counts up after counting down to zero
- Memory recall





### 4-way Timer with Whiteboard/Pen

- Clock feature count up or down to 24 hours 59 minutes
- 4 timer channels
- Stands up or attaches to a magnetic surface (magnets on back of unit)
- Whiteboard and marking pen included
- Sounds when time expires

	Item No.	No. of Channels	Timing Capacity	Count Direction	Display Digits, H mm	Battery Included	L x W x H in	L x W x H cm
а	HS24490D	Four	100 hr.	Up/Down	6	One LR44 1.5V	2.4 x 2.8 x 1	6.1 x 7.1 x 2.5
b	HS24490M	Three	20 hr.	Up/Down	10.5	One AAA	3.4 x 2.6 x 0.9	8.6 x 6.6 x 2.2
С	120365	Four	24 hr.	Up/Down	10	Two AAA	4.8 x 3.4 x 0.9	12.3 x 8.6 x 2.1
d	HS24645	One	24 hr.	Up	8.5	One LR44 1.5V	3.4 x 3.1 x 0.8	8.7 x 7.8 x 2.1
е	HS24670	Four	99 hr.	Down	14	Two LR44 1.5V	2.9 x 2.3 x 0.5	7.3 x 5.9 x 1.3





### Multi-function Stopwatch/Clock

- Quick start, quick stop timer features time, split, hourly chime, calendar, alarm
- 100 cm neck cord





### **Four Channel Timer/Clock**

- Independently or simultaneously set time on four channels
- Audible and visual alarms signify elapsed time
- Easily toggle between channels and clock with top buttons
- Includes spring clip with magnetic back

### **Lab Drying Rack**

### Bringing color to a mundane task

### Most Versatile Rack in the market


- Channels designed for more useful configurations
- Large and Small pegs for all glassware varieties
- One of a kind Drying Basket to collect smaller items

### Designed to protect against breakage

- Peg angle allows for better hold especially large, heavy pieces
- Patent Pending SureGrip Hook securely grips glassware into place

### Included:

- 72 Pegs ~ 36 Small (Green) ~ 36 Large (Blue)
- 2 Drying Baskets
- Drainage Tubing
- Mounting Hardware

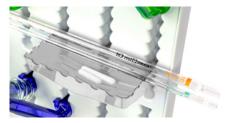

**Drying Basket** 

to keep small items and slotted to hold pipets safe while drying



### SureGrip Hook

Patent-pending design securely grips glassware into place



### Use on counter or wall

Stands alone on integrated hinged stand. Mounts easily on wall with included guide and hardware

Can hold up to 77 lbs

(34.9 kg)

Rack made of durable polystyrene



### Lab Drying Rack

U.S. Patent No. D923,275

Item No.	L x W x H in	L x W x H	
120791	17.4 x 26 x 5.5	44.1 x 66 x 13.9	

### **Replacement Parts**

Item No.	Description
120793	Peg 10 Piece Kit, (5) Small, (5) Large
120792	Drying Basket Kit, 2 ea
120813	Replacement Draining Set, Tube and Elbow

### **Wire Drying Rack**

Saves glassware from chips and cracks

Drain and dry test tubes, beakers, flasks, and other labware

- Place self-standing rack on a table or near a sink, or use hooks to hang rack on the wall
- Holds up to 90 pieces of labware on 18 x 127 mm pins and  $72 \times 76$  mm pins
- Constructed of spot-welded, HDPE-coated steel wire for durability with polystyrene tray and back
- Easily accessible and removable bottom drainage pan



Item No.	L x W x H in	L x W x H cm
HS23243A	18.2 x 7.4 x 20.7	46.2 x 18.2 x 52.5



Workstation Lab Mat, Page 47

### **Workstation Lab Mat**

Keep glassware in place while drying

### **Peg Racks**

**Perfect for drying small items** 

Fiberglass-reinforced polypropylene racks are ideal for bench work or tube storage

- Use in sink or on countertop
- Sturdy pegs hold items in place





Open tear drop wells allow for water to drain





Pack of 2

Item No.	Color	Tubes mL	Wells	Well Dia. mm	L x W x H in	L x W x H cm
HS24311A	White	10	50	14–17	7.4 x 4.1 x 2.7	18.8 x 10.4 x 6.9
HS24311B	Blue	10	30	14-17	7.4 X 4.1 X 2.7	18.8 x 10.4 x 6.9
HS24312A	White	7	06	10–13	3.2 x 5 x 2.7	8.2 x 12.6 x 6.9
HS24312B	Blue	1	96	10-13	3.2 X 3 X 2.1	0.2 X 12.0 X 0.9



Item No.	Description	Color		L x W x H	L x W x H cm
HS23500	MagLab Series, Pack of 4	Assorted	*	_	_
HS23502	MagWipe	Green		5.1 x 3.5 x 3.8	13 x 9.1 x 9.8
HS23503	MagPette	Purple		5.1 x 4.4 x 2.9	13.1 x 11.2 x 7.6
HS23504	MagBox	Orange		5.1 x 3 x 4	13 x 7.8 x 10.2

\*Assorted colors: Teal, Green, Purple and Orange

plates that adhere to any surface

### **Workspace Organization**

Studies\* show that by clearing clutter from your work environment, you are better able to focus and process information, resulting in increased productivity.

\*Princeton University Neuroscience Institute



### **Workstation Organizer**

Keeps supplies organized and easily accessible

For easy access to items on counter, shelf or inside drawer

- Protect fragile items and contain small items

- Flip-top covers on generous openings make filling and dispensing easy and

keep products dust-free

Item No.	Color	L x W x H in	L x W x H
HS2538A	White $\square$	17.9 x 12.1 x 2.5	45.5 x 30.7 x 6.3

Compartments	L x W	L x W cm
2 Small Compartments	17 x 3.1	43 x 8
1 Large Compartment	17 x 4.3	43 x 11





Item No. Color Compartments LxWxH LxWxH HS23401 7.2 x 7.6 x 9.1 18.4 x 19.4 x 23.2 HS23402 2 5.9 x 6.9 x 14.4 15 x 17.7 x 36.8 Clear HS23403 3 8.5 x 6.6 x 8.5 21.7 x 17 x 21.7 HS23404 4 28.7 x 17 x 21.7 11.3 x 6.6 x 8.5

### **Tubby® Storage Container** Sturdy storage tub with lid and dividers

Our handy Tubby® stackable storage container is designed to help organize a laboratory


- The removable dividers can create as big or as small a compartment as needed
- Perfect for keeping tubes, tips, pipettes, gloves, and other small items organized and out of the way
- The lid snaps in place to secure the contents and for safe transport

Pack of 5

Item No.	L x W x H in	L x W x H cm
HS23453	13 x 7.9 x 4.5	33 x 20 x 11.5



### More Dispensers for organization and to keep safety PPE products at the ready

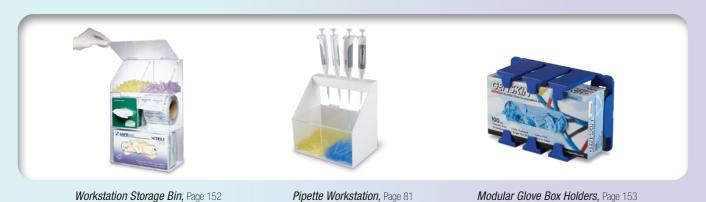

### **Test Tube Basket**

Ideal for washing and storing test tubes

Strong basket for use with mechanical washers and autoclaves

- Basket for use in washers, autoclaves or for storage of lab consumables
- Made from aluminum



Item No.	L x W x H in	L x W x H cm
HS20341B	10 x 6 x 6	25.6 x 15.2 x 15.2



### Glassware

High-quality glassware made from heavy-duty borosilicate glass and meets ASTM standards



Beakers
Separatory Funnels
Graduated Cylinders
Volumetric Flasks

### **Glass Beakers**

These low-form Griffin beakers feature easy-pour spouts and permanent white graduations

- Uniform wall thickness ideal for heating applications
- Made from high-quality borosilicate tubing with a low coefficient of expansion
- Graduation accuracy ±5% of full capacity

### **SIBATIA Beakers**

Item No.	Capacity (mL)	Graduation Range (mL)	Graduation Interval (mL)	Approx. O.D. x H	Approx. O.D. x H	Pack Qty.
10020-50A	50	20 to 50	10	1.8 x 2.4	46 x 60	10
10020-100A	100	20 to 100	10	2.2 x 2.8	55 x 70	10
10020-200A	200	50 to 200	25	2.6 x 3.5	67 x 89	10
10020-300A	300	50 to 300	50	3 x 4	78 x 103	10
10020-500A	500	100 to 500	50	3.5 x 4.3	90 x 120	10
10020-1000A	1000	200 to 1000	100	4.3 x 5.9	110 x 150	6
10020-2000	2000	250 to 2000	250	5.3 x 7.9	135 x 200	1
10020-3000	3000	500 to 3000	250	6 x 8.9	153 x 225	1

- Beakers are double-graduated metric

Meet ASTM E-960 Type I Specifications for Classification, Design, Capacity & Dimensions and Markings **VEE GEE Beakers** 

Item No.	Capacity (mL)		Graduation Interval (mL)	Approx. O.D. x H	Approx. O.D. x H		Case Qty.
20229-20	20	5 to 15	5	1.7 x 1.8	32 x 45	12	48
20229-3000	3000	250 to 2500	125	5.9 x 8.3	149 x 212	1	6



### **Glass Burets**

Used to accurately dispense a liquid, for example in a manual titration

- Designed "To Deliver"
- Made from high-quality borosilicate tubing with a low coefficient of expansion
- Feature durable, permanent markings; fine, sharp lines; and large easy-to-read numbers
- Supplied with PTFE stopcocks to eliminate freezing and lubricant contamination
- Meet ASTM E-287 Class A Specifications

### SIBATIA Burets

Item No.	Capacity (mL)	Graduation Interval (mL)	Tolerance (±mL)	Stopcock Bore (in/mL)	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Case Qty.
21310-10	10	0.05	0.02	0.2 / 4	0.4 x 20	9.1 x 570	1	20
21310-25	25	0.10	0.03	0.2 / 4	0.4 x 25.2	11 x 640	1	20
21310-50	50	0.10	0.05	0.2 / 4	0.5 x 33	12.6 x 840	1	25

### **VEE GEE Buret**

Item No.		Graduation Interval (mL)				Approx. O.D. x H mm	Pack Qty.	
20114-100	100	0.20	0.10	0.08 / 2	0.7 x 30.5	17 x 775	1	24

### **Separatory Funnels**

Feature smooth tooled tips, large white block letters, and large marking squares

- Supplied with PTFE stopcocks to eliminate freezing and lubricant contamination
- Made from high-quality borosilicate tubing with a low coefficient of expansion

### SIBATA Separatory Funnels with Stopper - Squibb Style

Item No.	Capacity	Stopper	Stopcock	Approx.	Approx.	Approx. Total		Case
142CT 200	(mL)	No. 24 / 40	Plug Size	0. <b>D.</b> x <b>L</b> in	<b>0.D. x L</b> <i>mm</i>	Length (in/mm) 15.41 x 390	Qty.	Qty.
1426T-300	300		4	******			- 1	0
1426T-500	500	24 / 40	4	0.4 x 2.8	10 x 70	17 x 430	<u> </u>	<u>b</u>
1426T-1000		29 / 42	6	0.5 x 3.1	12 x 80	19.3 x 490	1	4
1426T-2000	2000	29 / 42	6	0.5 x 3.1	12 x 80	20 x 510	1	2



Durable, permanent markings



PTFE stopcock



### **VEE GEE Separatory Funnels with Stopper - Squibb Style**

- All sizes supplied with a polyethylene standard taper stopper
- Meet ASTM E-1096 Type 4 Specifications for Classification, Design, and Markings

### VEE GEE Separatory Funnels with Stopper - Squibb Style

Item No.	Capacity (mL)	Stopper No.	Stopcock Plug Size	Approx. O.D. x L in	Approx. O.D. x L mm	Approx. Total Length (in/mm)	Pack Qty.	Case Qty.
20149-250	250	22	2	0.4 x 2.6	10 x 65	13.6 x 345	1	24
20149-500	500	27	4	0.4 x 2.6	10 x 65	15.4 x 392	1	242

### **FILTRATION** Filter Holders

Glass Microanalysis, 47 mm

SIBATA 6168-Series Filter Holders are offered with a wide variety of components for different applications. Available with fritted glass support base for general and biological applications, or with removable stainless steel screen support base and PTFE gasket for ultraclean filtrate or proteinaceous solutions. The borosilicate glass funnels are available in 100 and 300 mL.

### **Specifications**

- Effective Filtration Area: 9.6 cm2 (filter dia. 47 mm, not included)
- Maximum Capacity: 100 mL Funnel: 150 mL; 300 mL Funnel: 350 mL
- Graduated Volume: 100 mL Funnel: 100 mL; 300 mL Funnel: 300 mL
- Connection: Silicone Stopper (No. 8)
- Sterilization: Autoclavable, except PTFE Gasket (at 121°C / 250 °F for 20 minutes)

### **Filtration Kits**

Cat No.	Kit Includes	
	6168-4711	Funnel, 300 mL
6168-0001	6168-4722	Fritted Glass Support Base
0100-0001	6168-4751	Silicone Stopper No. 8
	6168-4762	Aluminum Anodized Clamp
	6168-4711	Funnel, 300 mL
	6168-4721	Glass Support Base
6168-0002	6168-4731	Stainless Steel Screen
0100-0002	6168-4741	PTFE Gaskets (10 pk)
	6168-4751	Silicone Stopper No. 8
	6168-4762	Aluminum Anodized Clamp
	6168-47100	Funnel, 100 mL
6168-0003	6168-4722	Fritted Glass Support Base
0100-0003	6168-4751	Silicone Stopper No. 8
	6168-4762	Aluminum Anodized Clamp
	6168-47100	Funnel, 100 mL
	6168-4721	Glass Support Base
6168-0004	6168-4731	Stainless Steel Screen
0100-0004	6168-4741	PTFE Gaskets (10 pk)
	6168-4751	Silicone Stopper No. 8
	6168-4762	Aluminum Anodized Clamp





**6168-47100**Funnel
100 ml



6168-4711 Funnel 300 mL



6168-4762 Aluminum Anodized Clamp



6168-4731 Stainless Steel Screen



**6168-4741**PTFE Gaskets (10 pk)



Base

C1C0 4751

**6168-4751**Silicone
Stopper No. 8



6168-4722 Fritted Glass Support Base

### **SIBATA Erlenmeyer Flasks**

- Made from high-quality borosilicate tubing with a low coefficient of expansion
- Uniform wall thickness ideal for heating applications
- Graduation accuracy is  $\pm 5\%$  of full capacity

### **SIBATIA Erlenmeyer Flasks**

Item No.	Capacity (mL)	Graduation Range (mL)	Graduation Interval (mL)	Approx. O.D. x H in	Approx. O.D. x H mm	Approx. O.D. of Neck (in /mm)	Rubber Stopper No.	Pack Qty.
10530-100A	50	20 to 50	10	2.5 x 4.1	64 x 105	0.9 / 22	2	10
10530-300A	300	20 to 100	25	3.4 x 6.1	87 x 156	1.3 / 34	6	10
10530-500A	500	50 to 200	50	4.1 x 7.1	105 x 180	1.3 / 34	6	10
10530-1000A	1000	50 to 300	50	5.2 x 8.7	131 x 220	1.7 / 42	8	10
10530-2000A	2000	100 to 500	200	6.5 x 11	166 x 280	2.0 / 50	10	1



### **SIBATA Filtering Flasks**

### Constructed of heavy wall borosilicate glass to withstand high vacuum

- Used in conjunction with vacuum filter systems, these filtering flasks feature heavy-duty rims to reduce chipping
- Feature white permanent graduations and large marking square
- Port for tubulation has an outside diameter of 10 mm

### **SIBATIA Filtering Flasks**

Item No.	Capacity (mL)	Approx. O.D. x H in	Approx. O.D. x H mm	Rubber Stopper No.	Pack Qty.	Pack Qty.
1780T-300	300	3.4 x 5.5	86 x 140	8	1	12
1780T-500	500	4.1 x 6.7	105 x 170	8	1	12
1780T-1000	1000	5.2 x 8.7	131 x 220	8	1	8
1780T-2000	2000	6.5 x 10.6	166 x 270	8	1	4



### **SIBATA Class A Graduated Cylinders**

Class A, Double-Metric Scale, Graduated, To Contain, Calibration Temperature 20°C/68°F

- Graduated Cylinders feature hexagonal bases for added stability
- All sizes, except the 2000 mL size, come with a protective, polyethylene hexagonal bumper guard
- Constructed of borosilicate glass with durable white enamel markings
- Exceeds ASTM E-1272 Class A, Style I Specifications for Classifications and Tolerances

### **SIBATA Class A Graduated Cylinders**

Item No.	Capacity (mL)	Graduation Range (mL)	Graduation Interval (mL)	Tolerance (±mL)	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Case Qty.
2351A-5	5	0.5 to 5.0	0.1	0.05	0.5 x 5	12 x 126	4	48
2351A-10	10	0.7 to 10.0	0.1	0.05	0.5 x 7.2	13 x 182	4	48
2351A-20	20	1.0 to 20.0	0.2	0.05	0.6 x 7.9	16 x 200	4	48
2351A-25	25	1.4 to 25.0	0.2	0.15	0.7 x 8.3	17 x 210	4	48
2351A-50	50	3.5 to 50.0	0.5	0.25	0.9 x 8.6	24 x 218	4	48
2351A-100	100	5.0 to 100.0	1	0.40	1.2 x 9.3	30 x 236	4	48
2351A-250	250	10.0 to 250.0	2	1.00	1.6 x 11.9	41 x 301	4	24
2351A-500	500	25.0 to 500.0	5	2.00	2.1 x 14.3	53 x 362	2	12
2351A-1000	1000	50.0 to 1000.0	10	3.00	2.8 x 16.6	70 x 422	2	8
2351A-2000	2000	100.0 to 2000.0	20	6.00	3.5 x 18.9	90 x 478	1	6

### **SIBATA Class B Graduated Cylinders**

Class B, Double-Metric Scale, Graduated, To Contain, Calibration Temperature 20°C/68°F

- Graduated Cylinders feature hexagonal bases for added stability
- All sizes, except the 2000 mL size, come with a protective, polyethylene hexagonal bumper guard
- Constructed of borosilicate glass with durable white enamel markings
- Exceeds ASTM E-1272 Class B, Style I Specifications for Classifications and Tolerances

### SIBATA Class B Graduated Cylinders

Item No.	Capacity (mL)	Graduation Range (mL)	Graduation Interval (mL)	Tolerance (±mL)	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Case Qty.
2351-10	10	0.7 to 10.0	0.1	0.10	0.5 x 7.2	13 x 182	4	48
2351-25	25	1.4 to 25.0	0.2	0.30	0.7 x 8.3	17 x 210	4	48
2351-50	50	3.5 to 50.0	0.5	0.40	0.9 x 8.6	24 x 218	4	48
2351-100	100	5.0 to 100.0	1	0.60	1.2 x 9.3	30 x 236	4	48
2351-250	250	10.0 to 250.0	2	1.40	1.6 x 11.9	41 x 301	4	24
2351-500	500	25.0 to 500.0	5	2.60	2.1 x 14.3	53 x 362	2	12
2351-1000	1000	50.0 to 1000.0	10	5.00	2.8 x 16.6	70 x 422	2	8
2351-2000	2000	100.0 to 2000.0	20	10.00	3.5 x 18.8	90 x 478	1	6

### **SIBATA Student Graduated Cylinders**

Plastic Base, Double-Metric Scale, Graduated, To Contain, Calibration Temperature 20°C/68°F

- These economical graduated cylinders are made of two parts: a strong, borosilicate glass graduate and a detachable, sturdy plastic base which absorbs impacts and reduces breakage
- All cylinders come with a protective, polyethylene, hexagonal bumper guard
- All markings are in durable white enamel
- Exceeds ASTM E-1272 Class B, Style I Specifications for Classifications and Tolerances

### **SIBATA Student Graduated Cylinders**

Item No.	Capacity (mL)	Graduation Range (mL)	Graduation Interval (mL)	Tolerance (±mL)	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Case Qty.
2355-10	10	0.4 to 10.0	0.1	0.10	0.5 x 7	13 x 177	6	24
2355-25	25	1.4 to 25.0	0.2	0.30	0.6 x 8.3	16 x 210	6	24
2355-50	50	2.0 to 50.0	1.0	0.40	0.9 x 8.3	24 x 212	6	24
2355-100	100	5.0 to 100.0	1.0	0.60	1.2 x 9.1	30 x 232	6	24





### **SIBATA Class A Volumetric Flasks with Snap-Cap**

Feature machine-blown bodies and heavy-beaded, heavy-tubing necks that are tooled for snap-cap fit


- Designed "To Contain" with a temperature of calibration at 20°C/68°F
- Constructed of borosilicate glass with a low coefficient of expansion
- Graduation line is sharply defined and of uniform width
- Each flask has a large marking square and includes a polypropylene snap-cap
- Meet ASTM E-288 Class A Specifications & USP Standards for Volumetric Flasks

### SIBATA Class A Volumetric Flasks with Snap-Cap

Item No.	Capacity (mL)	Tolerance (±mL)	Snap-Cap No.	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Pack Qty.
2303A-5	5	0.02	1	0.8 x 2.3	21 x 59	6	72
2303A-10	10	0.02	1	1 x 3.3	26 x 85	6	72
2303A-25	25	0.03	2	1.5 x 4.3	37 x 110	6	72
2303A-50	50	0.05	3	1.9 x 5.4	49 x 136	4	48
2303A-100	100	0.08	3	2.3 x 6.5	58 x 164	4	48
2303A-200	200	0.10	4	2.9 x 8.1	74 x 205	4	24
2303A-250	250	0.12	4	3.1 x 8.9	78 x 225	4	24
2303A-500	500	0.20	5	3.8 x 10.2	96 x 260	2	12
2303A-1000	1000	0.30	6	4.8 x 11.8	122 x 300	2	8

### Snap-Cap for SIBATA Class A Volumetric Flask

Item No.	Description	Fits Class A Volumetric Flask	Color	Each Qty.
SC#1	Snap Cap Size #1	2303A-5, 2303A-10	Yellow	1
SC#2	Snap Cap Size #2	2303A-25	Orange	1
SC#3	Snap Cap Size #3	2303A-50, 2303A-100	Red	1
SC#4	Snap Cap Size #4	2303A-200, 2303A-250	Green	1
SC#5	Snap Cap Size #5	2303A-500	Blue	1
SC#6	Snap Cap Size #6	2303A-1000	Yellow	1

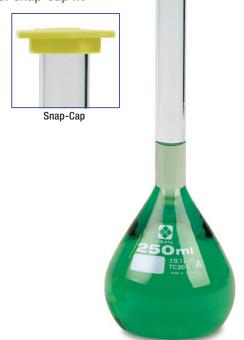

### **SIBATA Class A Volumetric Flasks with Ground Glass Stopper**

Feature machine-blown bodies and heavy-beaded, heavy-tubing necks that are tooled for standard taper glass stoppers

- Designed "To Contain" with a temperature of calibration at 20°C/68°F
- Constructed of borosilicate glass with a low coefficient of expansion
- Graduation line is sharply defined and of uniform width
- Each flask has a large marking square and includes a glass stopper
- Meet ASTM E-288 Class A Specifications & USP Standards for Volumetric Flasks

### SIBATA Class A Volumetric Flasks with Ground Glass Stopper

Item No.	Capacity (mL)	Tolerance (±mL)	Snap-Cap No.	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Pack Qty.
2306A-1	1	0.01	9	0.4 x 3.1	10 x 79	6	96
2306A-2	2	1.015	9	0.5 x 3.5	12 x 90	6	96
2306A-5	5	0.02	9	0.8 x 2.8	21 x 72	6	72
2306A-10	10	0.02	9	1 x 3.8	26 x 96	6	72
2306A-25	25	0.03	9	1.5 x 4.8	37 x 123	6	72
2306A-50	50	0.05	13	1.9 x 6	49 x 152	4	48
2306A-100	100	0.08	13	2.3 x 7.1	58 x 181	4	48
2306A-200	200	0.10	16	2.9 x 8.7	74 x 220	4	24
2306A-250	250	0.12	16	3.1 x 9.4	78 x 239	4	24
2306A-500	500	0.20	19	3.8 x 11	96 x 279	2	12
2306A-1000	1000	0.30	22	4.8 x 12.9	122 x 327	2	8
2306A-2000	2000	0.50	27	6.2 x 15.6	157 x 396	2	8



### **SIBATA Amber Class A Volumetric Flasks**

Constructed of low-actinic amber glass, providing protection for materials that are sensitive to lig

- Feature machine-blown bodies and heavy-beaded, heavy-tubing necks that are tooled for standard taper glass stoppers
- Constructed of borosilicate glass with a low coefficient of expansion
- Designed "To Contain" with a temperature of calibration at 20°C/68°F
- Each flask has a large marking square and includes a glass stopper
- Meet ASTM E-288 Class A Specifications & USP Standards for Volumetric Flasks

### SIBATA Amber Class A Volumetric Flasks

Item No.	Capacity (mL)	Tolerance (±mL)	Snap-Cap No.	Approx. O.D. x H in	Approx. O.D. x H mm	Pack Qty.	Pack Qty.
2307A-10	10	0.02	9	1.1 x 3.7	28 x 95	6	72
2307A-20	20	0.03	9	1.4 x 4.5	36 x 115	6	72
2307A-25	25	0.03	9	1.5 x 4.8	39 x 122	6	72
2307A-50	50	0.05	13	1.9 x 5.9	49 x 150	4	48
2307A-100	100	0.08	13	2.4 x 7.1	60 x 180	4	48
2307A-200	200	0.10	16	3.0 x 8.7	75 x 220	4	24
2307A-250	250	0.12	16	3.1 x 9.4	78 x 239	4	24
2307A-500	500	0.15	19	3.9 x 11	100 x 280	2	12
2307A-1000	1000	0.30	22	5 x 12.8	127 x 325	2	8



### **SIBATA Glass Stoppers**

Stoppers are lightweight yet very strong to reduce chipping and breakage

- Flat top design allows the stopper to stand on its head, freeing the user's hands
- Fully interchangeable with other common flask stoppers and fit all standard labware of comparable standard taper size
- Meet ASTM E-675 Specifications

### **SIBATA Glass Stoppers**

Item No.	Standard Taper Stopper No.	Tinting	Approx. Length of Ground Zone (in/mm)	Height Above of Ground Joint (in/mm)	Fit Volumetric Flasks	Pack Qty.
03081	9	Clear	0.6 / 14	0.4 / 10	2306A-1, -2, -5, -10, -25	6
03082	13	Clear	0.6 / 14	0.5 / 13	2306A-50, -100	6
03083	16	Clear	0.6 / 16.5	0.7 / 18	2306A-200, -250	6
03084	19	Clear	0.7 / 19	0.7 / 18	2306A-500	6
03085	22	Clear	0.8 / 20.5	0.8 / 20	2306A-1000	6
03087	27	Clear	0.7 / 21	0.8 / 21	2306A-2000	6
03081-A	9	Amber	0.6 / 14	0.4 / 10	2307A-5, -10, -20, -25	6
03082-A	13	Amber	0.6 / 14	0.5 / 13	2307A-50, -100	6
03083-A	16	Amber	0.6 / 16.5	0.7 / 18	2307A-200, -250	6
03084-A	19	Amber	0.7 / 19	0.7 / 18	2307A-500	6
03085-A	22	Amber	0.8 / 20.5	0.8 / 20	2307A-1000	6





### Sample Collection & Transport

Make everyday easier and safer with products for sample and blood collection



Portable Trays
Totes & Transports
Clinical Tube Racks
Sample Trays & Racks



## Ratcheting Handle folds out of the way for easy access, even over tall sharps containers Supply Drawer Side Pockets

### Tray Inserts

come in two configurations and include an extra-deep compartment for a sharps container



Disposable Insert A for 13 mm tube rack



Disposable Insert B for 16 mm tube rack

US Patent No. D544703

### **Droplet Blood Collection Tray**

**Designed with phlebotomists in mind** 

Keeping a clean and organized tray has never been easier. Wash the ABS tray as needed, and replace the insert periodically for a fresh start

- Multiple compartments keep a variety of supplies neatly at hand, and organized
- Built in side pockets offer additional storage space
- Handle includes soft rubber grip to reduce hand and finger strain
- Each tray includes either a 13 or 16 mm Heathrow OneRack® tube rack (page 115) and two disposable polystyrene inserts

**Droplet Configurations** 

D. 0p.0t 00	ropiot comigurations							
Item No.	Tray Insert	Tube Rack	L x W x H	L x W x H	Handle in	Height cm		
HS2200A	HS2201A	13	— 16.8 x 10.4 x 6.7	40.0 × 06.5 × 17.0	-11	20		
HS2200B	HS2201B	16	10.0 X 10.4 X 0.7	42.9 x 26.5 x 17.2	11	28		

Droplet Inserts Pack of 4

Item No.	Description	L x W in	L x W cm	
HS2201A	Tray insert for 13 mm tube rack	— 14.2 x 8.2	36 x 21	
HS2201B	Tray insert for 16 mm tube rack	14.2 X 0.2	30 X Z I	



Keeps contents clean and secure

### **Droplet Lid**

Cover for the Droplet blood collection tray

Protect the contents of the Droplet Tray with this Droplet lid. Custom designed sturdy lid is a snug fit, protects contents and provides access to the tray handle in its vertical position

- Made of polystyrene
- Can be secured with cable ties for extra security
  - Stackable

Item No.	Color	L x W x H	L x W x H	
HS2202	White	16.5 x 9.4 x 8	42 x 24 x 20.5	

### X-tourn® Tourniquets

Hi-viz orange alert for safer practice

### High-quality disposable tourniquets eliminate roping or curling when used

- Made from supple and smooth polyisoprene
- Latex free
- Less abrasive than other brands
- No rubber proteins, resins or fatty acids that can produce hypersensitive reactions
- Available in bulk or on a perforated roll

Item No.	Color	Packaging	L x W in	L x W cm
120108	Green	100/Bulk	18 x 1	45.7 x 2.5
120110	Orongo	250/Bulk	10 v 1	45.7 × 0.5
120111	− Orange ■	25 Rolls	– 18 x 1	45.7 x 2.5











### **DuraPorter®**

**Specimen or Sample Transport Box** 

Ideal for transporting biological and clinical samples, as well as instruments and products that you want to keep clean and dry under testing environmental conditions

The Duraporter® transparent design is OSHA 29 CFR 1910.1031 and UN3373 compliant

- Polycarbonate construction and available in three colors for user/ owner identification
- Accommodates 13 mm, 16 mm and 20 mm tubes with our Half and Full Size OneRack®
- Use the separators provided to create your own compartment sizes




Item No. Color		L x W x H	L x W x H
120052	Clear with Blue Handle		
120077	Red with Red Handle	15 x 7.7 x 6.3	38 x 19.6 x 16
120078	Yellow with Yellow Handle		









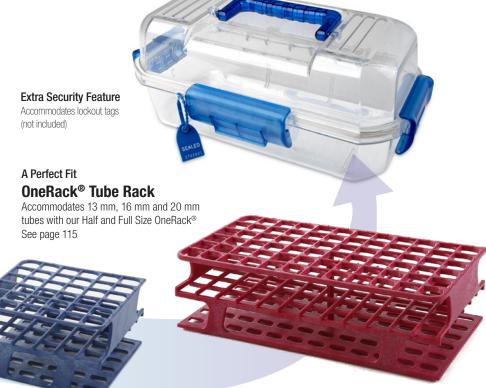

### For field sample processing, learn more about our Portable Centrifuge Kit

Includes Mini-Centrifuge and 12V Car Adapter for processing outside the laboratory



### **DuraPorter® XL**

### **High-Capacity Sample Transport Tote**

Your reliable solution for the safe and convenient transportation of laboratory samples, instruments, and valuable items. Crafted to meet the stringent demands of clinical and research laboratories, this oversized transport container ensures the protection and preservation of your cargo

### Built to Last and Designed with the User in Mind

- Sturdy and spacious, ideal for transporting a large volume of samples
- Transparent design helps comply with OSHA 29 CFR 1910.1030 standards without having to re-label containers and enabling the user to visibly check sample integrity and meets UN3373 secondary packaging transport requirements
- Equipped with articulated handles featuring silicone grips for comfortable and secure handling



### **Emphasizing Security and Safety**

- A silicone seal ensures a tight enclosure, safeguarding samples from most environmental elements
- Dual latches with security tag feature guarantee a secure closure, providing tamper-evidence and ensuring the integrity and authenticity of your samples



### **Exceptional Versatility**

- Includes a MultiTube Rack Insert to accommodate various tube sizes.
   Additionally, two other insert options for blood tubes up to 10 mL and general research tubes up to 15 mL are available separately
- The bottom compartment is generously sized, allowing for carrying most tube racks, instruments or supplies according to your specific requirements



# Soft Grip Handle for comfortable handling Removable Lid for ease of loading and unloading 4-Point Handle Security ensures durable transport Silicone Seal for secure content enclosure Dual Latches with tag out feature for added security

### DuraPorter® XL

Image shown with 120961 Small Tube Rack Insert (not included)

### Includes One 120962 MultiTube Rack Insert

Item No.	Color	Dimensions (not including handle) in	Dimensions (not including handle) cm	UOM
120926	Blue	16 x 16.25 x 7	40.64 x 41.28 x 17.78	1 ea
120927	Red	16 x 16.25 x 7	40.64 x 41.28 x 17.78	1 ea
120928	Yellow	16 x 16.25 x 7	40.64 x 41.28 x 17.78	1 ea

### **DuraPorter® XL**

**Tube Rack Inserts Accessories** 



### **120960**Duraporter® XL Large Tube Insert


Holds 100 tubes with diameter up to 18 mm Most 10 mL blood tubes & 15 mL conical tubes Array 10 x 10 Pack of 3



### 120961

Duraporter® XL Small Tube Insert

Holds 168 tubes with diameter up to 13 mm Most 3.5 - 7.5 mL blood tubes Array 12 x 14 Pack of 3

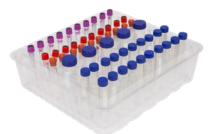

### 120962

Duraporter® XL MultiTube Rack Insert

Holds tube with various diameters: 13, 18 & 30 mm
Most 15 mL & 50 mL conical tubes
and 3.5 mL - 10 mL blood tubes
Pack of 3

### **DuraPorter® XL Tube Rack Inserts**

Item No.	Tubes	Cut-out Diameter	UOM
120960	15 mL, Large blood tube	18 mm	Pack of 3
120961	Small blood tubes	13 mm	Pack of 3
120962	50 mL, 15 mL, Small & Large blood tubes	30 mm (center), 18 mm (large), 13 mm (small)	Pack of 3

### **Blood Collection Tray**

Light weight, heavy on organization

Organize your phlebotomy, I.V. start or catheter insertion supplies with this spacious, robust and easy to clean carry tray

- Lightweight easy-to-clean ABS construction
- Customize the compartments to suit your requirements with the dividers provided
- Comfort grip on durable acetal plastic handle tested to hold up to 10 kg

Item No.	L x W x H in	L x W x H	Handle Height in	Handle Height cm
120261	16.4 x 12.2 x 3.3	41.9 x 31.1 x 8.5	9.4	24.1





### For safe disposal of syringes, needles and other small sharps see our Sharps Containers



Sharps Containers, Page 150



Item No.	Color	L x W x H in	L x W x H
120021	Blue	17 x 11 x 11	43.2 x 27.9 x 27.9

### **Totemax® Blood Collection Tray**

Keep off the floor, work top or bed

Unique phlebotomy tray design and patented hook feature allows the tray to be on or over most bed rails, foot-boards and handrails found in a typical patient room

- Made from food grade HDPE resins for a durable non-breakable design
- Deep tray design helps reduce spills






### **Clinical 4-Way Tube Rack**

The only tube rack specifically designed to hold round, flat & conical blood collection and culture tubes

- Accommodates all brands of 13 and 17 mm round, flat and conical bottom blood collection & culture tubes
- Customize the rack to fit the number of wells required, infinite units can be linked together or use one standalone
- Molded alphanumeric identification for quick location and retrieval of tubes
- Robust, made from polypropylene that is chemical resistant and easy to clean
- Racks can be used in water baths and placed in freezers
- Patent Pending

### Customizable

Racks can slide together to create many combinations

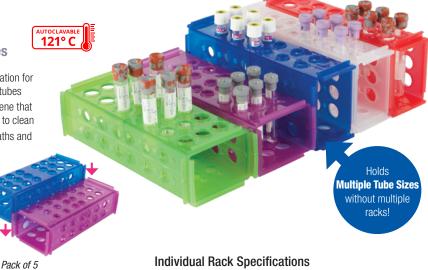

US Patent No. D900337

Item No.	Color	Item No.	Color	
120684	Blue	120687	Natural	
120685	Purple	120683	Assorted*	
120686	Green			

\*Assorted colors (Blue, Purple, Green, Red, Natural)


	ide	Wells	Well Dia.	L x W x H	L x W x H cm	
3. G		12	13			
ALCO DE			1.0			

В		mm	in	ст
	12	13		
A	18	13	01 x 2 0 x 2 0	23.2 x 7.0 x 9.5
D	12	17	9.1 X 2.0 X 3.0	23.2 X 7.0 X 9.3
The second second	18	17		

### Mega Rack®

### **High-capacity storage for your tubes**

Fully assembled, polypropylene mega racks are great for high throughput labs requiring large sample storage

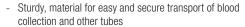

- Alphanumeric identification for quick location and retrieval of tubes
- Place different test groupings into one rack eliminating the need for multiple racks



### Mega Racks® 10-13 mm - Holds 5-7 mL sample tubes

Item No.	Color	Wells	Well Dia. mm	L x W x H in	L x W x H cm	
HS24322A HS24322E	Blue White	<b>— 216</b>	13	13.3 x 9.3 x 3.6	33.7 x 23.5 x 9.1	
120112 120113	White  Blue	<b>-</b> 432	13	26.5 x 9.3 x 3.5	67.3 x 23.5 x 9.0	

### Mega Racks® 13-16 mm - Holds 10 mL sample tubes

Item No.	Color	Wells	Well Dia. mm	LxWxH in	L x W x H cm
HS24323A	Blue	<b>– 120</b>	16	13.3 x 10.5 x 3.6	33.7 x 26.7 x 9.1
HS24323E	White	- 120	10	13.3 X 10.3 X 3.0	33.7 X 20.7 X 9.1
120114	White	<b>- 240</b>	16	26.5 x 10.5 x 3.5	67.3 x 26.7 x 9.0
120115	Blue	- 240	10	20.3 x 10.3 X 3.3	07.3 X 20.7 X 3.0

### Mega Racks® 18-20 mm - Holds 15 mL conical/swab tubes and 20 mm round bottom tubes

Item No.	Color	Wells	Well Dia. mm	L x W x H in	L x W x H cm	
120834	White	— 120	20	14.9 x 11.9 x 3.6	37.7 x 30.1 x 9.1	
120835	Blue	120	20	14.9 X 11.9 X 3.0	37.7 X 30.1 X 9.1	



### **Tube Roller Racks**

### Efficient and practical solution for tube roller workflows

### Expand capacity of your workflow

Each Tube Roller Rack accommodates 10 tubes, doubling or tripling the capacity of your tube roller equipment

- Saves Process Time
   Large batch processing reduces the number of roller cycles required for lab procedures
- Saves Money
  Handles more without the need to invest
  - in multiple machines to process the same amount of tubes
- Saves Valuable Space
   Additional capacity especially helpful with limited incubator capacity



- Organize samples pre and post mixing
- Use as a storage tube rack for upright storage or submerse in an ice bath
- Polypropylene construction is wearresistant and offers reliable sterilization through autoclaving



Silicone rings for smooth rotation and gentle contact with roller

Accommodates 10 tubes

Tubes snap securely into holders, ensuring a stable hold and aligning with the roller's rotational axis







Fits most standard and larger size tube roller equipment models.

### Available in two sizes

- Large: Fits flat and conical bottom blood collection & culture tubes (10 mL) and standard 15mL tubes up to 18-20 mm
- Small: Fits flat and conical bottom blood collection & culture tubes 10 17 mm (3.5 7 mL)

Item No.	Tube Size	Unit Diameter in (cm)	Holding Prongs Diameter (at rest) in (cm)	Unit Length in (cm)	UOM
120956	Large	2.54 in (6.5 cm)	$0.59 \pm 0.01$ in (1.5 $\pm 0.02$ cm)	11 in (27.9 cm)	Pack of 3
120957	Small	2.54 in (6.5 cm)	0.468 ±0.008 in (1.2 ±0.02 cm)	9.4 in (23.9 cm)	Pack of 3

### **Clinical Tube Racks**

Learn about more of our Tube Racks that hold 10-17 mm blood collection tubes





### **Droplet Corrugated Sample Storage Tray**

Affordable sample organization and storage

Corrugated storage trays can be used with the included dividers for specimen containers or remove them for general lab product storage or transport

- Made from corrugated polypropylene

the chance of mold

- Stackable with other trays
- Moisture resistant, significantly reducing
- Pack includes: 10 trays and 10 dividers

Item No.	Color	Wells	Well Dia.	L x W x H	L x W x H
120375	White $\square$	24	58	16 x 10.2 x 3.5	40.6 x 26 x 8.9



### Removable Inserts

allow outer shell to be used for storage or transport of laboratory products



Lightweight and Sturdy

to hold up to 24 full specimen containers



### **Droplet Sample Storage Tray**

Sample organization for the busy lab

Hard wearing storage trays that can be used for specimen containers with the insert provided, or for general lab products without the insert

- Bin and tray are made from polystyrene
- Lightweight and sturdy to hold up to 24 full specimen containers (Specimen containers not included)
- Stackable with other trays to support as many as four high
- Can be easily cleaned with common disinfectants and germicides



### Removable Insert

allows outer shell to be used for storage or transportation of laboratory products

Well Dia.

**Label Holder** 





### **Coated Wire Rack - Urine Containers**

**Durable racks for transport or sample handling** 

- Chemically-resistant, HDPE-coated steel
- Fits bottles with 55-58 mm diameters

- Mesh bottom support



Item No.	Color	Tube Array	Wells	Well Dia.	LxWxH	LxWxH
				mm	in	ст
120091	Blue	6 x 4	24	58	15.5 x 10.3 x 2	39.3 x 26.2 x 5

### **Liquid Sample Handling**

Whether it be dispensing or aspirating microliters or milliliters, we offer you precision at your fingertips




Pipettors, Fillers & Bulbs
Glass & Transfer Pipettes
Reagent Reservoirs & Chillers
Funnels & Lead Rings
Pipette Racks & Accessories

# **High-quality Pipettors offering precision at your fingertips**





# **High-quality Pipettors offering precision at your fingertips**

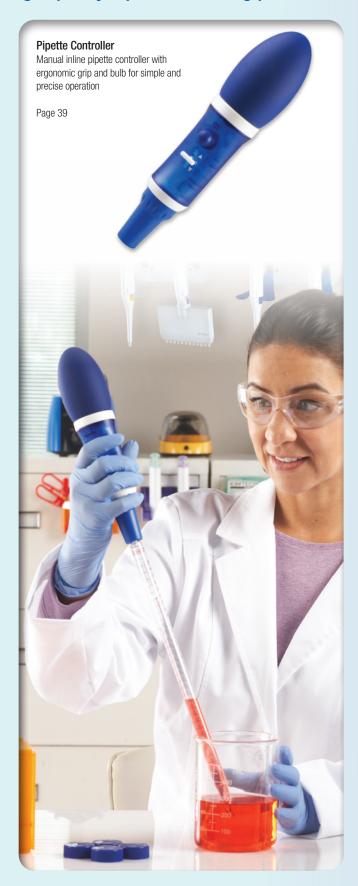



## **Fillers and Bulbs**

- 1 Latex Dropper Bulbs
  - Designed for use with Pasteur pipettes, small pipettes, and medicine droppers
- 2 Dropper Bulbs
  - Designed for use with Pasteur pipettes, small pipettes, and medicine droppers
- Rubber Bulbs
  - Contoured shape and heavy-wall neck increase grip for easier use with pipettes and syringes
- Vacuum/Pressure Bulb
  - Ideal for gas sampling or for pumping gas through apparatus
  - Two rubber valves with PVC connections (tubing not supplied)



- without twisting

   Pipette does not contact hard plastic, which reduces the risk of cracking and chipping
- Elastic band for hanging on reagent bottle or wrist
- Filler fits standard size pipettes

# Pipette Filler

- Efficient one-handed operation
- Control air-evacuation, liquid uptake, and liquid dispensing with three valves
- Pinch release-valve for smooth liquid flow
- Red pipette filler is safe for pipetting of caustics and does not contain carbon black and will not smudge



Latex Dropper Bulbs (Natural Rubber)

Pack of 72

Item No.	Color	<b>Volume</b> mL	Length in	Length mm	Outer Diameter in	Outer Diameter mm	Hole Diameter in	Hole Diameter mm
HS20622A	— Yellow	1	1.7	43	0.51	13	0.20	5
HS20622B		2	2.0	52	0.51	13	0.24	6

Dropper Vial Bulbs (Silicone Rubber)

Pack of 24

Item No.	Color	Volume mL	Length in	Length mm	Outer Diameter in	Outer Diameter	Hole Diameter in	Hole Diameter
HS20629A	— Blue	1	1.1	28.7	0.47	12	0.24	6.0
HS20629B	Diue	3	1.6	41	0.70	17.8	0.26	6.5

#### Rubber Bulbs (Natural Rubber)

Item No.	Color	Volume mL	Length in	Length mm	Outer Diameter in	Outer Diameter	Hole Diameter	Hole Diameter mm
HS20630B		15	1.9	50	1.3	34	0.16	4
HS20630C	- Pod -	30	2.6	66	1.8	46	0.24	6
HS20630D	— Red ■	60	3.1	80	2.1	54	0.24	6
HS20630E	_	150	4.1	103	2.9	73	0.24	6

#### Vacuum/Pressure Bulb (Natural Rubber)

Item No.	Color	Volume	Length	Length	Outer Diameter	Outer Diameter	Hole Diameter	Hole Diameter
		mL	III	mm	in	mm	III	mm
HS20631B	Black	85	4.2	106	2.1	52	0.29	7.4

#### Safety Bulb Pipette Filler (Natural Rubber)

Pack of 3

Item No.	Color	Volume mL	Length in	Length mm	Outer Diameter in	Outer Diameter mm	Hole Diameter in	Hole Diameter mm
HS20623	Red	15	4.2	106	1.6	41	0.39	10

#### Pipette Filler (Natural Rubber)

Item No.	Color	<b>Volume</b> mL	Length in	Length mm	Outer Diameter in	Outer Diameter mm	Hole Diameter in	Hole Diameter
HS20632A	Black	_						
HS20632B	Red	35	5.1	130	2.1	54	0.24	7
HS20632C	Blue							



# **Resuable Glass Pipets**

Our SIBATA reusable glass pipets are meticulously manufactured for tolerances and accuracy.

Glass is a superior material for measuring samples with a wide variety of chemical compounds.

Compared to plastic, reusable glass pipets are autoclavable and sustainable with less waste.

These quality pipets offer markings and color coding that are permanently enameled to the glass.

## **Volumetric Pipets** – To Deliver

Accurately transfer a specific volume of solution for which it is calibrated.

- Meet ASTM E-969 Class A Specifications and USP Standards for volumetric glassware

#### **VOLUMETRIC PIPETS**

Item. No.	Capacity (mL)	Tolerance (±mL)	Approx. Length (in/mm)	Color Code	Pack Qty.	Case Qty.
2040A-05	0.5	0.005	11.8 / 300	2-Black	12	72
2040A-1	1.0	0.006	13.5 / 345	Blue	12	72
2040A-2	2.0	0.006	13.9 / 355	Orange	12	72
2040A-3	3.0	0.01	14.1 / 360	Black	12	72
2040A-4	4.0	0.01	14.3 / 365	2-Red	12	72
2040A-5	5.0	0.01	14.7 / 375	White	12	72
2040A-6	6.0	0.01	15.3 / 390	2-Orange	12	72
2040A-7	7.0	0.01	16.1 / 410	2-Green	12	72
2040A-8	8.0	0.02	16.9 / 430	Blue	12	72
2040A-10	10.0	0.02	17.3 / 440	Red	12	72
2040A-13	13.0	0.03	18.1 / 460	White	6	36
2040A-15	15.0	0.03	18.1 / 460	Green	6	36
2040A-20	20.0	0.03	18.7 / 475	Yellow	6	36
2040A-25	25.0	0.03	18.8 / 480	Blue	6	36
2040A-50	50.0	0.05	21.2 / 540	Red	6	36
2040A-100	100.0	0.08	23.6 / 600	Yellow	6	36

## **Volumetric Pipets - Wide-Tip Opening** – To Deliver

Ideal for water treatment facilities, winemakers, tissue culture work and essential oil testing

- Meet ASTM E-969 Class B Specifications for tolerance
- Wide-tip opening facilitates rapid fill and drain of viscous liquids, slurries and suspensions
- No drawn tips on 1, 2 & 5 mL sizes. Standard drawn tips on 10-100 mL sizes

#### **VOLUMETRIC PIPETS - WIDE-TIP OPENING**

	-					
Item No.	Capacity (mL)	Tolerance (±mL)	Approx. Length (in/mm)	Approx. Tip I.D. (in/mm)	Pack Qty.	Case Qty.
2040-1-1.9MM	1.0	0.01	12.6 / 320	0.08 / 1.9	12	72
2040-2-2.4MM	2.0	0.01	13 / 330	0.09 / 2.4	12	72
2040-5-3MM	5.0	0.02	13.6 / 346	0.1 / 3	12	72
2040-10-3MM	10.0	0.02	17.3 / 440	0.1 / 3	12	72
2040-20-3MM	20.0	0.03	18.8 / 478	0.1 / 3	6	36
2040-25-3MM	25.0	0.03	18.7 / 475	0.1 / 3	6	36
2040-50-3MM	50.0	0.05	21.2 / 538	0.1 / 3	6	36
2040-100-3MM	100.0	0.10	23.8 / 605	0.1 / 3	6	36





Minimum IO mIS SERTA

## **Serological Pipets**

Graduated with calibration marks for delivering variable volumes of solution as needed

## **Mohr Measuring Pipets** – To Deliver, DRAIN OUT

- Meet ASTM E-1293 Style 1, Class A Specifications and USP Standards for volumetric glassware

#### MOHR MEASURING PIPETS

Item No.	Capacity (mL)	Tolerance (±mL)	Graduation Interval (mL)	Approx. O.D. x Height (in)	Approx. O.D. x Height (mm)	Color Code	Pack Qty.	Case Qty.
2010A-5	5.0	0.02	0.1	0.3 x 14.4	8 x 365	Blue	12	72
2010A-10	10.0	0.03	0.1	0.4 x 14.6	10 x 370	Orange	12	72
2010A-25	25.0	0.05	0.1	0.5 x 17.9	13 x 455	White	12	72

## **Serological Pipets** – To Deliver, BLOW OUT

- Meet ASTM E-1044 Style 1, Class A Specifications

#### **SEROLOGICAL PIPETS**

Item No.	Capacity (mL)	Tolerance (±mL)	Graduation Interval (mL)	Approx. O.D. x Height (in)	Approx. O.D. x Height (mm)	Color Code	Pack Qty.	Case Qty.
2011A-01-C	0.1	0.0025	0.01	0.2 x 11.8	6 x 300	White	12	72
2011A-1-C	1.0	0.01	0.01	0.3 x 14.2	7 x 360	Yellow	12	72
2011A-2-D	2.0	0.01	0.1	0.3 x 14.2	7 x 360	Green	12	72
2011A-5	5.0	0.02	0.1	0.3 x 14.4	8 x 365	Blue	12	72
2011A-10	10.0	0.03	0.1	0.4 x 14.6	10 x 370	Orange	12	72
2011A-25	25.0	0.05	0.1	0.5 x 17.9	13 x 455	White	12	72

## **Serological Pipets-Wide-Tip Opening** – To Deliver, BLOW OUT

**Ideal for biochemical applications** 

#### **SEROLOGICAL PIPETS - WIDE-TIP**

Item. No.	Capacity (mL)	Tolerance (±mL)	Graduation Interval (mL)	Hole Diameter x Length (in)	Hole Diameter x Length (mm)	Color Code	Pack Qty.	Case Qty.
20200-11A	1.0	0.02	0.2	0.06 x 13.8	1.5 x 350	Red	12	72
20200-51A	5.0	0.1	0.5	0.08 x 13.8	2.0 x 350	Blue	12	72
20200-101A	10.0	0.2	0.5	0.1 x 13.8	3.0 x 350	Orange	12	72

#### SEROLOGICAL PIPETS - GRADUATED TO THE TIP

Item. No.	Capacity (mL)	Tolerance (±mL)	Graduation Interval (mL)	Hole Diameter x Length (in)	Hole Diameter x Length (mm)	Color Code	Pack Qty.	Case Qty.
20140-10A	10	0.05	0.1	0.05 x 14.4	1.4 x 366		12	72
20140-20A	20	0.1	0.1	0.06 x 17.1	1.5 x 435		12	72

## **Transfer Pipettes**

One-piece, unbreakable polyethylene plastic pipettes

Well-suited for blood banking, urinalysis, hematology, wet chemistry, microbiology, and serology

- Bulk packed versions come in handy dispenser box
- Sterile versions available single wrapped or 20 wrapped at a time
- Lot-to-lot traceability on sterile versions

	Item No.	Vol. mL	Working Vol. mL	Length mm	Graduations mL	Packaging 500 per box	Sterile
	HS206371A					Bulk	
a	HS206371B	5	1	150	0.25	20/bag	✓
	HS206371C					1/single wrap	✓
	HS206372A					Bulk	
b	HS206372B	5	2	155	0.5	20/bag	✓
	HS206372C					1/single wrap	✓
	HS206373A					Bulk	
С	HS206373B	7	3	155	0.5	20/bag	✓
	HS206373C					1/single wrap	✓





## **Universal Carousel Pipette Stand**

Pipettes never looked so good, on your bench top


Stylish pipette stand holds and stores up to 8-single channel, 4-multi channel, or a combination of both pipettes

#### Versatile

- Patent pending design enables this stand to securely hold virtually every known brand of pipette in the marketplace
- Modular upper holder inserts accommodate the fit of different pipette widths
- Assorted color inserts allow for color coding configurations

#### **Durable Utility**

- Stable weighted base allows for easy one handed removal and replacement of pipettes
- Rotating head makes pipette selection fast and easy
- Moving the stand while keeping pipettes in place is simple with a quick grip handle allowing a user to grab and go

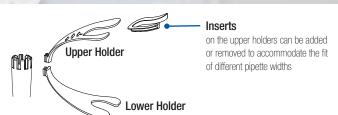

Item No.	Upper Holder Color	Quantity	Lower Holder Color	Quantity
120480	Blue/Green	8 ea/color	Blue	8 ea
120481	Purple/Clear	8 ea/color	Purple <	8 ea
120499	Assorted/Clear*	1 ea/color	Clear	8 ea
		8 clear		

\*Assorted colors (Purple, Blue, Green, Grey, Red, Orange, Yellow and Clear)

Base Width in	Top W x H	Base Width	Top W x H
	in	cm	cm
6.2	9 x 16.1	15.8	22.9 x 40.8

## **Universal Pipette Holder**

**Wall mounted storage solution** 

The pipette holder comes with two wall mounting brackets that each hold two pipettes

#### Secure

- Upper and lower segment inserts easily slide in place to securely hold virtually every known brand of pipette
- Segments unique design includes traction grips that center pipettes allowing for easy one handed removal and replacement

#### Adaptable

Item No.

120575


**Lower Holder** 

2 Blue

2 Green

2 Clear

 Inserts come in two separate colors allowing for color coding, quick identification or just to add a pop of color to your workspace

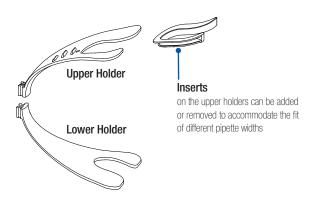

**Upper Holder** 

2 Green

1 Purple 1 Yellow 1

2 Blue

LxH

6.9 x 1.8





## Clamp-ette® Pipette Holder

Create an instant pipette holding point within your lab

Secures to tabletops, benches and shelves with a twist

Chemical-resistant polypropylene and fits most major pipette brands






Secures Clamping for tabletops, benches and shelves

#### Tri Clamp-ette®

Item No.	Color	L x W x H in	L x W x H cm
HS20162A	Blue	6.2 x 4.7 x 3.7	15.7 x 12 x 9.5

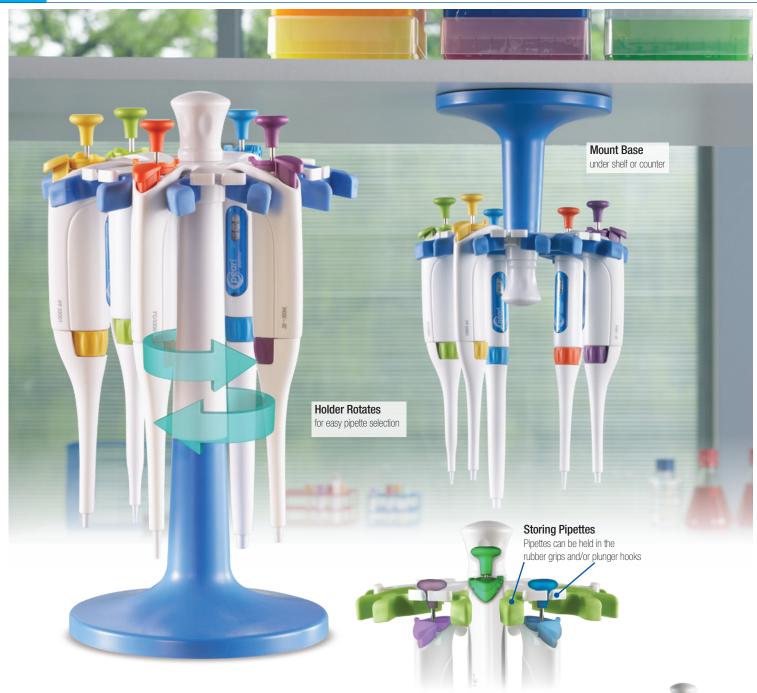

# Flip & Grip® Pipette Holder Giving labs a choice

Unique pipette holder that gives you the option of using on the bench top or hanging under shelf to help maximize your working space  ${\bf p}$ 

- Will hold up to 12 single-channel pipettes, or a combo of 3 single- and 3 multichannel pipettes
- Holds most common brands
- Rubber grips prevent scratching pipette
- Weighted base
- Installation hardware included

US Design Patent No. D754,872 & D763,464

Item No.	Color	Bench Top Style W x H <i>in</i>	Bench Top Style W x H cm	Under Shelf Style W x H in	Under Shelf Style W x H cm
120026	Blue	— 13.1 x 7.6	33.2 x 19.3	7.4 x 7.6	18.7 x 19.3
120056	Lime	13.1 X 1.0	33.2 X 19.3	1.4 X 1.0	10.1 X 19.3





## **Pipette Workstation**

Organize your pipettes and tips

#### Keep pipette and disposable tips together in one convenient stand

- Acrylic workstation is opaque on sides with clear front and bin lids for easy content identification
- Holds up to four pipettes on top shelf
- Dual bins hold multiple size tips
- Self-closing lid provides dust-free storage
- Oversized lip allows for easy access when wearing gloves

Item No.	Color	L x W x H in	L x W x H cm
HS234631	Clear	9.6 x 6.5 x 9.7	24.4 x 16.5 x 24.6



## **Pipette Tip Disposal Box**

Quick flip of port cover makes tip ejection easy

Convenient for one-handed disposal of used pipette tips, eliminates any need to handle pipette tips used with hazardous fluids

- The 0.375 in (10 mm) thick acrylic provides protection from Beta radiation
- The large, hinged lid allows easy disposal of contents
- Use in conjunction with small biohazard disposal bags, if hazardous fluids are not acrylic friendly

Item No.	Lid Opening	L x W x H	L x W x H
	mm	in	cm
HS234634	30	6.1 x 6.7 x 8.1	15.4 x 17.1 x 20.6



## **Pipette Stands**

**Convenient and safe storage** 

3-, 4-, and 6-Place Pipette Stands are designed to hold most major pipette brands

- Acrylic stands are transparent/clear; ABS stands are solid opaque colors

Item No.	Color	Material	Pipette Capacity	Slot Upper in	Slot Upper mm	Slot Lower in	Slot Lower mm
HS206203			Three				
HS206204	Clear	Acrylic	Four		0.5	0.0	10
HS20620			Six	ı	25	0.8	19
HS20613C	Blue	Acrylic	Six				
HS20613E	Blue	ABS	Six	1	25	0.8	20

## **Cargo Pipette Rack**

Meets the need of any lab


Racks securely lock together horizontally or vertically to desired configuration or can be used individually


Can be used on top of benchtop or conveniently hung on any metallic surface with strong magnets that are designed into the racks

- Angled shelves store 0.1 to 50 mL pipettes
- Protective lid encloses the unit to prevent dust build-up and provides additional shelf space

Pack of 3

Item No.	Color	L x W x H in	L x W x H cm
HS20615M	Blue	11.9 x 3.9 x 5	30.1 x 10 x 12.8

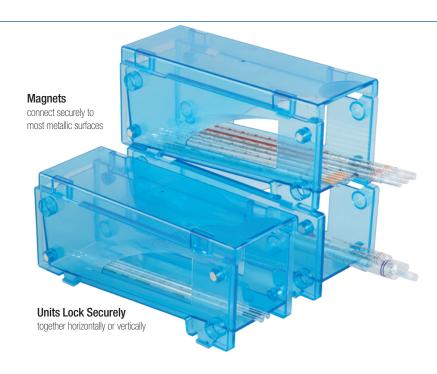

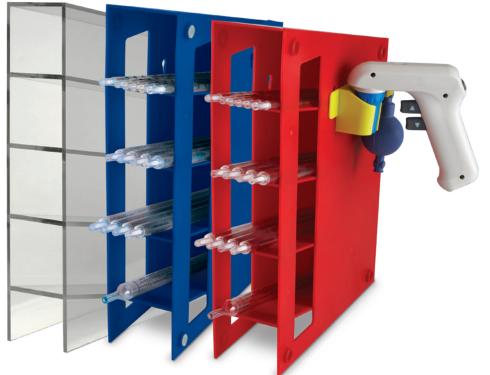

## **Pipette Racks**

Protects and organizes glass and plastic pipettes

Four angled-shelved compartments keep your 0.1 to 50 mL pipettes or cryo vial canes within easy reach while preventing damage and clutter

- Heavy-duty acrylic version provides crystal-clear viewing from any angle
- ABS versions have cut-out viewing window and come with two removable attachments that hold pipette bulbs, fillers, and pumps
- HS20612M magnetized ABS version has strong magnets on both sides and mounts securely to steel surfaces

Item No.	Color	Material	L x W x H in	L x W x H cm
HS20611	Clear	Acrylic		
HS20612	Blue	ABS	11 4 4 2 5 4 15 0	201 v 0 0 v 40 1
HS20613	Red	ABS	- 11.4 x 3.5 x 15.8	29.1 x 8.8 x 40.1
HS20612M	Blue	ABS (mounting magnets)		



## **Liquid Weighing Pour Boats**

Deep vessel with sturdy handle designed for convenient handling and pouring of liquids

- Flat bottom ensures stability and prevents tipping
- Notched pour spout design controls liquid pour
- Operating temperature range 0 to 40°C (104°F)
- Made from FDA food grade polystyrene in an ISO Class 8 clean room
- Material biologically inert and resistant to dilute and weak acids, aqueous solutions, alcohols and bases
- Recyclable, marked
- Available in 3 Sizes
  - 50mL 100mL 200 mL

Pack of 200

Item No.	Color	Size	Volume mL	L x W x H in	L x W x H cm
120890	White	Small	50	3.6 x 1.6 x 1.7	9.2 x 4.0 x 4.3
120891	White	Medium	100	4.2 x 3.5 x 1.9	10.7 x 9 x 4.8
120892	White	Large	200	5.8 x 4.6 x 2	14.8 x 11.7 x 5.1



#### **Reservoir Chiller**

Ice-free, non-toxic way to extend working period and preserve integrity for temperature sensitive reagents

- Reagents are kept at less than 5°C (41°F) for approximately 1 hour
- Engineered custom fit design for Heathrow Scientific 50 or 100 mL reservoirs to ensure a maximum contact for longer uniform cooling
- Eliminates the risk of contamination from messy ice bath
- Not for use with liquid nitrogen

Item No.	Color	Volume mL	L x W x H in	L x W x H cm
120730	Purple			
120731	Blue	50	6.1 x 2.9 x 1.7	15.5 x 7.4 x 4.2
120732	Assorted*			
120733	Purple			
120734	Blue	100	6.2 x 9 x 1.7	15.7 x 9 x 4.2
120735	Assorted*	_		
				*Assorted colors (Purple and Blue)

Pack of 2



**Symmetrical cooler design**Allows reagent reservoirs to be mounted in either direction



**Dual Solution Reservoir** 

Innovative shape guides and protects pipette tips

Multi-functional reagent reservoir has a single compartment 50 mL basin on one side and 12 slotted 5 mL basins on the other side

- Angled wells guide pipette tips toward a center position to help avoid cross contamination and waste
- Individual 5 mL basin wells numerical imprinted to identify samples
- Cover included to prevent contamination





Pack of 25



Item No.	Volume mL	LxWxH in	L x W x H
HS20821A	5/50	4.9 x 2 x 1.5	12.5 x 5 x 3.7

## **Reagent Reservoirs**

Reagent Reservoirs by Heathrow Scientific offer a fresh and modern design with enhanced features to streamline pipetting procedures

- Superior durability with enhanced walls and crush-resistant "turrets" at the corners
- Easy handling with tabs to prevent stacked reservoirs from sticking together

#### **Sterile Versions**

- -Manufactured and gamma irradiated in an ISO certified facility
- -Certified DNase, RNase and pyrogen free
- -Packaged individually in an easy open HDPE package
- -Traceability Certificate available upon request

STERILE R



allows pipetting of two different reagents or maximizes recovery if using only single channel and small amount of reagent



on the top edge and inner walls allows measurements of liquid to be easily identified and reduces wastage of expensive reagents



conveniently located on all four corners enables easy spill-free pouring and reduces chance of dripping



Item No.	Volume mL	Graduation Marks mL	Material	Sterile	Packaging	Interior Dimensions L x W x H <i>in</i>	L x W x H
120637	5	5	Polystyrene	Yes	50/Single Bag	2.8 x 0.8 x 0.7	7.2 x 2.1 x 1.8
120638	10	5 and 10	Polystyrene	No	100/Bulk	2.8 x 0.8 x 0.7	7.2 x 2.1 x 1.8
120639	10	5 and 10	Polystyrene	Yes	100/Single Bag	2.0 x 0.0 x 0.7	1.2 A 2.1 A 1.0
120640	25	15 and 25	Polystyrene	Yes	200/40 x 5pk		
120641	25	15 and 25	Polystyrene	Yes	100/Single Bag	4.4 x 0.9 x 1.1	11.3 x 2.4 x 2.7
120642	25	15 and 25	Polystyrene	No	100/Bulk	4.4 X U.3 X 1.1	11.3 X 2.4 X 2.1
120884	25	15 and 25	Polystyrene	Yes	100/Single Bag		
120643	50	25 and 50	Polystyrene	Yes	200/40 x 5pk		
120644	50	25 and 50	Polystyrene	Yes	100/Single Bag		
120645	50	25 and 50	Polystyrene	No	100/Bulk	4.5 x 1.1 x 1.4	11.5 x 2.9 x 3.5
120646	50	25 and 50	PVC	Yes	100/20 x 5pk		
120690	50	25 and 50	PVC	No	100/Bulk		
120647	100	50 and 100	Polystyrene	Yes	200/40 x 5pk		
120648	100	50 and 100	Polystyrene	Yes	100/Single Bag	4.6 x 1.9 x 1.0	11.7 x 4.8 x 2.5
120649	100	50 and 100	Polystyrene	No	100/Bulk	4.0 X 1.9 X 1.0	11./ X 4.0 X 2.3
120650	100	50 and 100	PVC	No	100/Bulk		
120651	25 (10 and 15 mL divided)	10/15	Polystyrene	Yes	200/40 x 5pk	4 x 1 x 0.8	10.2 x 2.5 x 2.1
120652	25 (10 and 15 mL divided)	10/15	Polystyrene	Yes	100/Single Bag	4 1 1 1 0.0	10.2 x 2.3 x 2.1

# **Eco smartFunnel™ Paper Funnel**

A clean transfer every time

- Supplied flat-pack to minimize storage
- With a minimal amount of finger pressure the funnel takes its working shape
- Maximum flexibility cut the stem to any size or shape to meet with your needs
- Fits a wide variety of openings and can be used with virtually any liquid\* or powder



						Pa	ack of 25
Item No.				Bottom Dia.	Bottom Dia.	Height in	Height cm
120219	186	2.4	6	0.3	0.6	5	12.8

Eco smartFunnel™ is a trademark of LevGo Inc.

## **Lead Rings**

Stabilize vessels and labware on the bench or in water baths

- Vinyl-coated lead rings will not mark or scratch benchtops and labware
- C-shape format allows some manual size variation for a customized fit



Item No.	Shape	Color	Flask Sizes mL	ID in	ID cm	<b>OD</b> in	OD cm	Weight g
HS8881A	С	Blue	125-500	1.9	4.2	3	7.5	215
HS8882A	0	blue =	125-500	1.9	4.8	3	7.7	210
HS8881B	С	Blue	250-1000	2	5.1	3.5	9	465
HS8882B	0	blue =		2	5.1	3.5	9	490
HS8881C	С	Blue	500-2000	2.1	5.4	3.9	9.8	595
HS8882C	0	blue =	300-2000	2.2	5.7	3.9	9.8	610
HS8881D	С	Dlue	1000-4000	2.6	6.6	4.4	11.3	890
HS8882D	0	Blue	1000-4000	2.8	7	4.5	11.5	935

\*Weights may vary by +/- 3%

## **Safety Bottle Carriers**

Secure and convenient bottle transport

- Ergonomic handle prevents hand fatigue and related accidents when carrying heavier loads
- Handle swings down to either side for easier loading and unloading
- Centering ring stabilizes bottle in the center of carrier for safer, steadier transport
- Carriers come in four colors with designated labeling areas for easy identification
- Made of chemical resistant thermoplastic elastomer PP blend



Available in 4 Colors for easy identification and organization

#### Safety Bottle Carrier, 4.5 L

li	tem No.	Color	Dimensions (Dia. x H) in	Dimensions (Dia. x H)		
1	20911	Red				
1	20912	Black	7.2 x 10.5	18.3 x 26.7		
1	20913	Blue	(16 in height, handle extended)	(40.6 cm height, handle extended)		
1	20914	Green	<del></del>			

#### Safety Bottle Carrier, 1.5 L

ou.ory Do	ttio ouiii	o.,o =				
Item No. Color			Dimensions (Dia. x H) in	Dimensions (Dia. x H)		
120915	Red					
120916	Black		4.5 x 6.5	11.4 x 16.5		
120917	Blue		(10 in height, handle extended)	(25.4 cm height, handle extended)		
120918	Green		<del></del>			



**Movable Handle** swings down to allow clear access to vessel



# **Sample Handling**

Complete range of products and tools for gathering, weighing, and working with various types of samples



Pill Counting Trays
Weigh Boats, Dishes, Papers
Sampling Scoops

**Utensils: Spatulas, Forceps, Cell Spreaders, Loops & Picks** 

## **Pill Counting Trays**

Counting trays that offer safety and convenience when counting pills and capsules


Intended for single use to prevent cross-contamination for batch sampling or distribution of pills and capsules  $\frac{1}{2} \left( \frac{1}{2} \right) = \frac{1}{2} \left( \frac{1}{2} \right) \left( \frac{1}{2}$ 

#### Safety

- Disposable to ensure no chemical cross-contamination for patient safety
- Made from FDA food grade polystyrene in an ISO Class 8 clean room
- Meets quality control and GMP requirements for manufacturing batch sampling
- Durable and shatterproof trays constructed of inert, BSE/TSE free material
- Lot traceable

#### Convenience

- Pour spout in corner for dispensing control
- Ergonomic spatula (length 6.9 in/17.6 cm)
- Available in both, sterile\* and non-sterile for intended usage

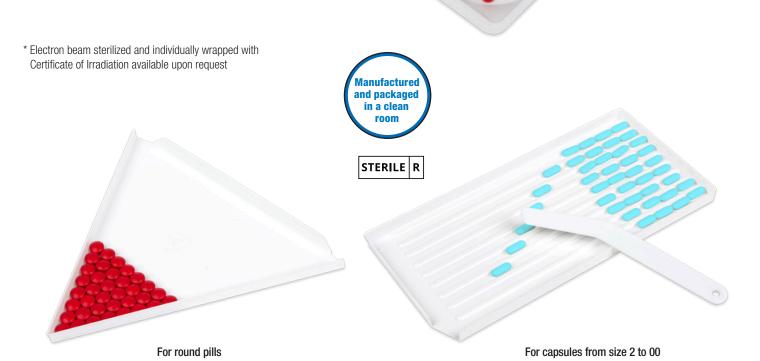

Item No.	Description	L X W X H in	L x W x H cm	UOM
120845	Pill Counting Tray Rectangular Sterile	10.3 x 5.2 x 0.5	26.2 x 13.1 x 1.2	Pack of 10 Case 40; 4 x 10/Single Bag
120846	Pill Counting Tray Rectangular Non-Sterile	10.3 x 5.2 x 0.5	26.2 x 13.1 x 1.2	Pack of 10 Case 40; 4 x 10/Bulk
120847	Pill Counting Tray Triangular Sterile	9.1 x 9.1 x 0.5	23.0 x 23.0 x 1.3	Pack of 10 Case 30; 3 x 10/Single Bag
120848	Pill Counting Tray Triangular Non-Sterile	9.1 x 9.1 x 0.5	23.0 x 23.0 x 1.3	Pack of 10 Case 30; 3 x 10/Bulk

## **Aluminum Weighing Boats & Dishes**

Dishes for general-purpose weighing, dispensing, storage, food processing, evaporating or drying

#### **Aluminum Weighing Boats**

 Sturdy yet lightweight boats feature a smooth, flat bottom for sample handling









#### **Aluminum Weighing Boats**

Item No.	Volume mL	Dimensions in	Dimensions mm	Packaging	Tab	Wall Finish
120937	140 mL	4.1 x 0.8	103 x 20	100		Smooth
120936	80 mL	2.9 x 0.7	73 x 18	100		Smooth
120930	70 mL	2.8 x 0.7	71 x 17	100		Crinkle
HS14521C	50 mL	2.8 x 0.6	70 x 16	100	Х	Crinkle
HS14521B	39 mL	2.2 x 0.6	57 x 14	100	Х	Crinkle
120931	20 mL	2 x 0.5	52 x 13	100		Crinkle
HS14521A	20 mL	1.7 x 0.5	43 x 12	100	Х	Crinkle
120932	8 mL	1.1 x 0.4	28 x 10	100	Х	Crinkle
120933	1.12 mL	0.98 x 0.3	25 x 7	50	,	Smooth

#### **Aluminum Weighing Dish**

 Designed for moisture determination balances with slightly tapered walls for stacking







HS14522




#### **Aluminum Weighing Dish**

Item No.	Volume mL	<b>Dimensions</b> in	Dimensions mm	Packaging	Tab
120934		2.7 x 0.2	70 x 6	50	
120935	50 mL	2.2 x 0.6	57 x 16	100	
HS14522	30 mL	4 x 0.3	101 x 8	50	
120940	0.35 mL	0.5 x 0.1	13 x 3.5	100	Х
120938	39 mL	3.5 x 0.3	90 x 7	50	

## smartPans™

Ideal for microbalances, organic micro elementary and thermal analysis, and pharmaceutical use

- Made of Ultra High Vacuum (UVH) Foil "A" wettable and completely free from oils normally found on other grades of aluminum foil
- Foil meets ASTM B-479, specification 3.1.4 and 10.3.1
- "Cahn style" accessory pans are compatible with most commercially available analyzers and microbalances

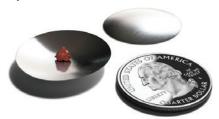



Pack of 20

Item No.	Dia. x H in	Dia. x H mm	
120627	1 x 0.1	25 x 3.6	

smartPans™ is a trademark of LevGo Inc.

## **Weighing Papers**

Smooth surface for good sample recovery

- Nitrogen-free
- Strong enough to support dense material





Pack of 500

Item No.	LxH	LxH	
	in	ст	
120116	3 x 3	7.6 x 7.6	
120117	4 x 4	10.2 x 10.2	
120118	6 x 6	15.2 x 15.2	

## **Weighing Boats**

Flexible polystyrene boats feature smooth surfaces to provide accurate pour-outs with minimal sample loss

- Select white, black, or blue for best color contrast to ensure sample is completely removed
- Biologically inert and resistant to dilute and weak acids, aqueous solutions, alcohols, and bases
- Temperature range -10 to 70°C (14 to 158°F)
- Antistatic versions available making the boats nonreactive with most substances and suitable for weighing staticaffected samples
- Wide flat bottoms ensure stability and prevent tipping
- Rounded corners allow simple flexing for controlled pouring of liquid, powdered or granular samples
- Can be used with microwave moisture analyzers
- Made from FDA food grade polystyrene





**Standard Weighing Boats** 

#### **Sterile Weighing Boats**

Ideal for applications where there is concern of cross contamination and 'clean' is a necessity in life science research, food, biopharmaceutical, environmental and molecular research environments

- Electron beam sterilized
- Lot number traceability for quality assurance
- Certificate of Irradiation available upon request

STERILE R



**Square Weighing Boats** 



**Hexagonal Weighing Boats** 



**Diamond Weighing Boats** 



**Pour-boat Weighing Dishes** 

## **New Liquid Weighing Pour Boats**

Flat bottom ensures stability and prevents tipping

- Notched pour spout design controls liquid pour
- Available in 3 Sizes
  - 50 mL
  - 100 mL
  - 200 mL



Weigh Boats							Pack of 500
Shape	Regular Material Item No.	Antistatic Material Item No.	Color	Size	Volume mL	L x W x H	L x W x H
	HS1420A	HS1420AA	White				
	HS1421A	HS120222	Blue	Small	10	1.8 x 1.8 x .3	4.6 x 4.6 x .8
	HS1423A	HS1423AA	Black <b></b>				
	HS1420B	HS1420BB	White				
Standard	HS1421B	HS120223	Blue	Medium	85	3.4 x 3.4 x .9	8.5 x 8.5 x 2.4
	HS1423B	HS1423BB	Black <b></b>				
	HS1420C	HS1420CC	White				
	HS1421C	HS120224	Blue	Large	280	5.5 x 5.5 x .9	14 x 14 x 2.2
	HS1423C	HS1423CC	Black ■				
	_	HS14251A		Small	6	1.5 x 1 x .4	3.8 x 2.5 x .9
Hexagonal	_	HS14251B	— — White □	Medium	46	2.5 x 2 x .7	6.4 x 5.1 x 1.9
пехауонан	_	HS14251C	— White $\square$	Large	194	4 x 3 x .9	10.2 x 7.6 x 2.2
	_	HS14251D		X-Large	364	5 x 3.5 x 1.2	12.7 x 8.9 x 3.1
	HS1419A	HS120225		Small	14	1.7 x 2.3 x .5	4.3 x 5.8 x 1.3
Pour-boat	HS1419B	HS120226	White	Medium	137	3.3 x 5.2 x 1	8.3 x 13.2 x 2.6
	HS1419C	HS120227		Large	259	4.3 x 7.2 x 1	10.8 x 18.3 x 2.6
Square	HS1422	_	White $\square$	Medium	51	3.5 x 3.5 x .3	9 x 9 x .7
	HS1424AA	HS1425AA	White	——— Small	5	2.1 x 1.3 x .2	5.3 x 3.3 x .6
	HS1427AA	HS1426AA	Black ■	Jillali	J	2.1 X 1.3 X .2	J.J X J.J X .U
	HS1424A	HS1425A	White	Medium	25	2.6 x 1.9 x .6	6.7 x 4.7 x 1.4
Diamond	HS1427A	HS1426A	Black ■	Mediaiii	23	2.0 % 1.9 % .0	0.7 X 4.7 X 1.4
Diamona	HS1424B	HS1425B	White	Large	30	3.1 x 2.2 x .6	7.8 x 5.6 x 1.4
	HS1427B	HS1426B	Black ■	Large	30	3.1 X 2.2 X .0	7.0 X 3.0 X 1.4
	HS1424C	HS1425C	White $\square$	X-Large	100	4.5 x 3.3 x 0.7	115 x 85 x 17
	HS1427C	HS1426C	Black ■	Λ-Laige	100	4.5 % 5.5 % 0.7	113 x 03 x 17
Sterile Weigh E	Boats						Pack of 200/40 x 5pk
Shape	Regular Material Item No.	Antistatic Material Item No.	Color	Size	<b>Volume</b> mL	L x W x H in	L x W x H cm
Standard	_	120880	White $\square$	Small	10	1.8 x 1.8 x .3	4.6 x 4.6 x .8
Standard	_	120710	White $\square$	Medium	85	3.4 x 3.4 x .9	8.5 x 8.5 x 2.4
Hexagonal	_	120881	White	Small	6	1.5 x 1 x .4	3.8 x 2.5 x .9
Havenenel							

Medium

## **Micro Weighing Boats**

Ideal for microbalance weighing and/or distribution of very small amounts of sampling material

#### Easy to Handle

Hexagonal

- Large tabs make boats easy to use, position and hold, especially with forceps

#### **Superb Sample Recovery**

- Inert material enables the boat and contents to be added to solutions to aid 100% sample recovery, where protocols allow
- Anti-static properties prevent sample retention from static electricity
- Small surface area reduces liquid sample evaporation

#### Micro Weigh Boats

Item No.	Color	Dimensions in	Dimensions mm	UOM
120836	White	.98 x 0.4 x 0.3	25 x 8.9 x 7.5	250/pk


120882



Made of Polystyrene Food Grade Compliant

White

**Easy to Handle Tabs** 

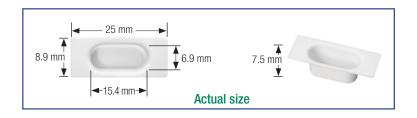


2.5 x 2 x .7

6.4 x 5.1 x 1.9

# Ideal for very small volume pipette calibration and checks and micro sampling



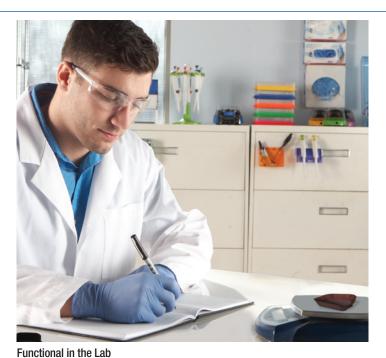

Available in a small and medium in both brown for light colored samples and white for dark samples

## smartBoats™

A truly green solution to sample weighing

Made from FSC® Certified (FSC® Mix Credit) paper (Eco-IN) and when disposed of in the standard waste stream they will naturally biodegrade (Eco-OUT)

- A cost-effective solution to reduce plastics in the lab
- Made of grease-proof paper that meets FDA 21 CFR 176.170 (components of paper in contact with aqueous and fatty foods) and with 2I CFR 176.180 (contact with dry foods)
- Control contents like a polystyrene weigh boat because they are easier to pour
- An alternative to regular weighing papers because they are easier to handle

Item No.	Color	Size	Packaging	L x W x H in	L x W x H cm	
120312	Brown	— Small	Pack of 500	1.4 x 1.4 x 0.6	3.5 x 3.5 x 16	
120314	White	Siliali	rack of 300	1.4 % 1.4 % 0.0	3.3 X 3.3 X 10	
120313	Brown	Modium	Pack of 500	1.9 x 1.9 x 0.6	4.8 x 4.8 x 14	
120315	White	- Mediuiii	Pack of 500	1.9 X 1.9 X 0.0	4.0 X 4.0 X 14	
120626	Brown	Lorgo	Dook of 250	20 v 20 v 1	9.5 x 9.5 x 2.5	
120195	White	— Large	Pack of 250	3.8 x 3.8 x 1	9.5 X 9.5 X 2.5	

smartBoats™ is a trademark of LevGo Inc.

HS15909













### **Forceps**

Corrugated handles and serrated tips

- ABS plastic

\* Forceps shown at actual size

Pack of 6

Item No.	Color	<b>Length</b> in	Length cm	Tip
HS15998C	Blue	4.4	11.1	Straight





## **Spatulas**

Transfer, lift, and handle powders and granular materials

- Made from nickel-plated, stainless steel wire
- Glazed spatula HS15906 has a Hayman-type end



Item No.	End Type	End L x W in	End L x W	Overall Length in	Overall Length cm	Finish
HS15906	Flat	1.2	30	6	15	Glazed
H313900	Spoon	0.5 x 0.2	12 x 4	O		
HS15907	Round Flat	2 x 0.3	50 x 8	7.7	19.5	Mirror
no 1090 <i>1</i>	Tapered Flat	1.8 x 0.3	45 x 6.5	1.1	19.5	
HS15909	Flat	2	50	8.9	00.5	Minnen
	Spoon	1.5 x 0.6	38 x 15	0.9	22.5	Mirror

## smartSpatula®

Designed to fit perfectly with standard 0.2, 0.5, 1.5/2.0, 15 and 50 mL tubes

- Made of virgin, food-grade polypropylene
- Ideal for dispensing or sample collection of powders or soft solids
- Will not shatter when used with liquid nitrogen frozen samples
- Hollow core minimizes heat transfer—good for hot and cold liquids
- Antistatic, sterile and Eco-Friendly versions available
- Complies with FDA standards for food contact





	Item No.	Description	Color	Suggested Tube mL	Packaging	L x W in	L x W cm	Spoon (Flare) in	Spoon (Flare) mm	Scoop/Pick L x W in	Scoop/Pick L x W mm	Taper
a	120005	Micro Antistatic	Opaque	0.2 and 0.5	300/Bulk	5.5 x .2	14 x .4	_	_	4.3	11 (V-Tip)	15°
b	120008	Eco-Friendly	Opaque	15 and 1.5	300/Bulk	5.5 x .3	14 x .7	9.8 x 5.5	25 x 14	11.8 x 1.2	30 x 3	
С	120003	Standard	Blue		300/Bulk							
	120004	Standard	Opaque	15 and 1.5	300/Bulk	- 00 40	0 01 11 7	21 x .7 9.8 x 5.5	25 x 14	11.8 x 1.2	30 x 3	
d	120007	Std. Et0 Sterile	Opaque	15 and 1.5	100/10 x 10pk	- 8.3 x .3	21 X ./					
	120625	Std. Et0 Sterile*	Opaque		100/10 x 10pk	_						
е	120006	Macro	Green	50	150/Bulk	12.2 x .4	31 x 1	16.5 x 9.3	42 x 23.5	11.8	30	35°

\*Non Toxic/Non Pyrogenic smartSpatula<sup>TM</sup> is a trademark of LevGo Inc.

## **Dual Purpose Scoops**

Just scoop to retrieve, place on scale to weigh

Unique economical solution that combines the function of a scoop and a weigh boat into one, streamlining the process

- Made from anti-static polystyrene that is waterproof and reduces the chance of materials sticking to the scoop
- Ships and stores flat, easy to assemble for use
- Locking tab feature enables a scoop shape to help retain materials
- Coated to stand up to powders, gels, and semi-solids



Case of 100 10 x 10pk



Item No.	Color	Size	L x W in	L x W cm
120656	White	Large	5 x 3.5	12.7 x 9
120657	White	Small	5 x 2.5	12.7 x 6.3



## **Sampling Scoops**

Our line of durable Sampling Scoops designed for easy sample retrieval of granules, powders and pastes

Available in multiple sizes and intended for single use to prevent contamination of bulk sample material

#### **Durable Design**


- Large, wide openings for easy sample retrieval
- Flat bottom to stabilize sample on scale or flat surface
- Lot traceable
- Shatterproof constructed of inert, BSE/TSE free material

#### Safe Sample Recovery

- Disposable to ensure no chemical cross-contamination
- Made from FDA food grade polystyrene in an ISO 8 Class clean room
- Meets quality control and GMP requirements for sampling
- Available in both, sterile\* and non-sterile for intended usage
- \*Electron beam sterilized and individually wrapped

Certificate of Irradiation available upon request



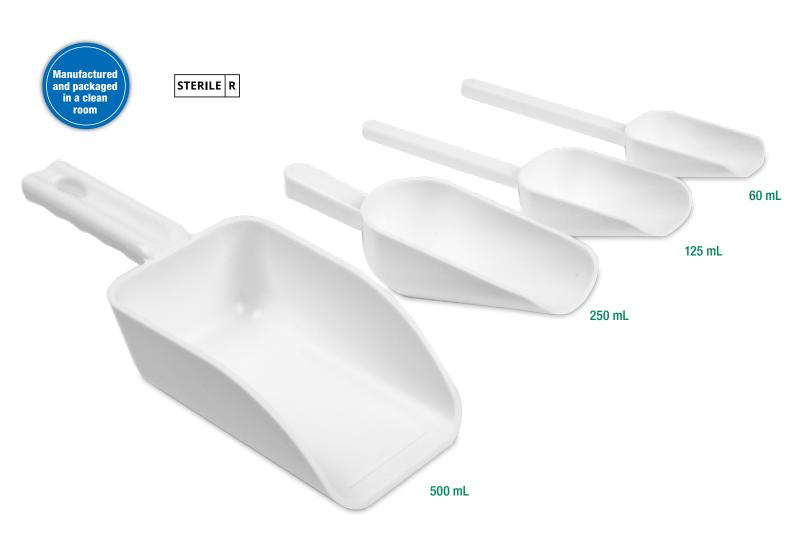





## **Small Sampling Scoop**

Item No.	Description	Volume	L x W x H in	L x W x H mm	Sterile	UOM
120858	Small Sampling Scoop Non-Sterile	60 mL/2 oz.	8.1 x 1.6 x 1.1	205 x 40 x 27		Pack of 10, Case 200; 20 x 10/Bulk
120857	Small Sampling Scoop Sterile	60 mL/2 oz.	8.1 x 1.6 x 1.1	205 x 40 x 27	✓	Pack of 10, Case 200; 20 x 10/Single Bag







## **Medium Sampling Scoop**

Item No.	Description	Volume	L x W x H	L x W x H	Sterile	UOM
120856	Medium Sampling Scoop Non-Sterile	125 mL/4 oz.	8.6 x 1.9 x 1.3	218 x 49 x 38		Pack of 10, Case 200; 20 x 10/Bulk
120855	Medium Sampling Scoop Sterile	125 mL/4 oz.	8.6 x 1.9 x 1.3	218 x 49 x 38	✓	Pack of 10. Case 200: 20 x 10/Single Bag





## Large Sampling Scoop

Item No.	Description	Volume	L x W x H in	L x W x H mm	Sterile	UOM
120854	Large Sampling Scoop Non-Sterile	250 mL/8 oz.	7.7 x 2.3 x 1.7	195 x 58 x 44		Pack of 10, Case 200; 20 x 10/Bulk
120853	Large Sampling Scoop Sterile	250 mL/8 oz.	7.7 x 2.3 x 1.7	195 x 58 x 44	✓	Pack of 10, Case 200; 20 x 10/Single Bag





## **XLarge Sampling Scoop**

Item No.	Description	Volume	L x W x H in	L x W x H	Sterile	UOM
120852	XLarge Sampling Scoop Non-Sterile	500 mL/16 oz.	10.0 x 3.6 x 2.1	255 x 92 x 53		Pack of 10, Case 200; 20 x 10/Bulk
120851	XLarge Sampling Scoop Sterile	500 mL/16 oz.	10.0 x 3.6 x 2.1	255 x 92 x 53	✓	Pack of 10, Case 200; 20 x 10/Single Bag

## smartPicks® Inoculation Picks

The 'smart' way to isolate individual colonies

Available in three lengths and two needle sizes to suit a wide range of application needs

- Naturally biodegradable, made from 100% U.S.A. grown white birch wood
- Manufactured in an FDA-registered factory and complies with FDA standards for food contact
- Available in semi or needle point styles for delicate operations





Mini Sized SmartPicks fit neatly into a 1.5 mL microfuge tube allowing the lid to be closed for full sample transfer/recovery

Drop it in & Close it up!



Item No.	Size	Point	Tip	Suggested Tube mL	Packaging	Length in	Length cm	
120010	Mini	Needle Point		Fits 1.5 tube (closed)	2,000	1.3	3.3	
120011	Standard	Semi Point	Cinalo Tin		1.600	2.5	6.4	
120012	Standard	Needle Point	Single Tip		1,000	2.3	6.4	

smartPicks® is a trademark of LevGo Inc.

## **Inoculating Needle/Loops**


Combination provides choice to spread or pick colonies and cultures

# Transfer bacterial cultures without risk of contamination

- Use loop end to uniformly streak a gel surface
- Use needle end for removal of colonies
- Choose between flexible polypropylene or the more rigid ABS versions
- Lot-to-lot traceability





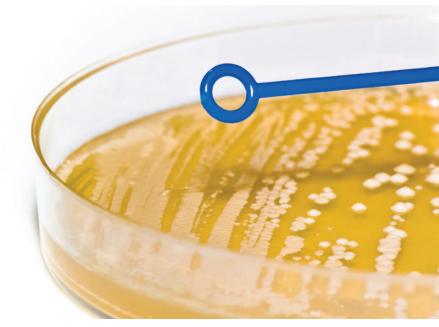

Item No.	Color	Material	<b>Loop</b> μL	L x D in	L x D cm	Packaging 1,000 per Box	Sterile
HS81121A	— Blue ■	Flexible PP	1	7.9 x .2	20 x .4	50 x 20pk	./
HS81121B	Diue	Rigid ABS	Į.	7.9 X .Z	20 X .4	ου χ Ζυμκ	•
HS81121C	— Blue ■	Flexible PP	10	7.9 x .2	20 x .4	50 x 20pk	./
HS81121D	Diue	Rigid ARS	10	7.9 X .2	20 X .4	50 X 20pk	V

## **Cell Spreaders**

#### Delta™

- Made from ABS
- Flexible rippled handle and domed spreading bar to help reduce agar tears

STERILE R

Domed Headpiece

promotes smooth spreading

#### L-Shape

- Made from ABS
- Lot-to-lot traceability
- 10 Pack comes in a zip bag

STERILE R

#### Stainless Steel

- Flame-sterilizable
- Reusable

AUTOCLAVABLE 121° C



Item No.	Shape	L x W x H in	L x W x H cm	Packaging	Sterile
HS8151	Delta™	6.8 x 1.5 x .2	17.2 x 3.7 x .4	500/Single Wrap	✓
HS8171A	I Chana	E 7 v 1 E v 0	14.6 x 3.8 x .4	500/Single Wrap	✓
HS8171B	L-Shape	5.7 x 1.5 x .2	14.0 X 3.0 X .4	500, 50 x 10 Pack	✓
HS86655	Stainless Steel	8 x 1.3 x .2	20.3 x 3.3 x .4	1/Single Wrap	

## **Petri Dish Dispenser**

- Sturdy acrylic holds up to 30 x 100 mm dishes in two easy-access partitions
- Set on counter or wall mount with included hardware

Item No.	LxWxH	LxWxH
	in	cm
HS23452	8.7 x 4.2 x 9.3	22.2 x 10.7 x 23.6





## Need help counting cells?

For accurate cell counts, check out our Digital and Manual cell counting solutions



# **Tubes & Vials**

Reliable sample storage



Micro & Conical Centrifuge Tubes
Screw-top Tubes
Cryogenic Vials

## **Micro & Conical Centrifuge Tubes**

Available in natural and pure black for light sensitive samples

# Microcentrifuge Tubes: 0.2, 0.5 and 1.5 mL (a-d)

- Attached flat caps for easy one-handed use
- 0.2 mL is suitable for thermal cyclers
- 0.5 and 1.5 mL versions have a frosted writing area and a syringe pierceable cap
- 100% polypropylene





Certified DNase, RNase & Pyrogen Free (a, b, c only)



#### Microcentrifuge Tubes

	Item No.	Color	Vol. mL	Outer Dimension mm	Closed Height mm	Graduations mL	Max RCF Rating	Packaging
a	HS4423	Natural	0.2	6.9	22.2	None	15,000 x g	1,000/Bulk
b	HS4422	Natural	0.5	10.1	31.7	0.1, 0.3, 0.5	15,000 x g	1,000/Bulk
С	HS4323	Natural	1.5	10.8	39.0	0.1, 0.5, 1	16,400 x g	500/Bulk
С	HS4325	Natural	1.5	10.8	39.0	0.1, 0.5, 1	16,400 x g	1,000/Bulk
d	HS4323K	Black	1.5	10.2	40.6	0.1, 0.5, 1, 1.5	18,000 x g	500/Bulk

## Conical Centrifuge Tubes:

15 and 50 mL (e-i)

- Lot-to-lot traceability with easy tear packaging
- Writing area and flat polyethylene lids for easy labelling
- Imprinted graduations and a self-standing version (only in clear bulk pack)
- Recyclable chipboard used for the racked versions, making them more environmentally friendly







Skirted

Certified DNase, RNase & Pyrogen Free (f & h)

#### **Conical Centrifuge Tubes**

	Item No.	Color	<b>Vol.</b> mL	Outer Dimension mm	Closed Height mm	Graduations mL	Max RCF Rating	Packaging 500 per Case
0	HS4426B	Natural	15	17	120	0.5	11,000 x g	Bulk/10 x 50pk
е	HS4426R	Natural	15	17	120	0.5	11,000 x g	Rack/10 x 50pk
f	HS4428	Black ■	15	17	120	0.5	9,400 x g	Bulk/10 x 50pk
~	HS4427B	Natural	50	29	114	5.0	9,500 x g	Bulk/20 x 25pk
g	HS4427R	Natural	50	29	114	5.0	9,500 x g	Rack/20 x 25pk
h	HS4429	Black ■	50	29	114	5.0	9,400 x g	Bulk/20 x 25pk
i	HS4427S	Natural	50	29 - skirted	114	5.0	9,500 x g	Bulk/20 x 25pk

## **Screw-top Tubes with 0-ring Cap**

#### Safer long-term storage

- Screw-on polypropylene lids with TPE O-ring (white)
- Clear polypropylene tubes provide excellent view of samples
- Skirted, self-standing base makes filling and retrieval easy
- Non-sterile







						1111		
Item No.	Color	Vol.	Outer Dimensions	Closed Height	Graduations	Max RCF Rating	Packaging	
		mL	mm	mm	mL			
HS10060	Natural	2.0	13.3	47.1	0.5, 1, 1.5	18,000 x q	500/Bulk	

## **Cryogenic Vials**

#### Innovative design saves space in your freezer

- Puncture/rupture-resistant polypropylene
- Silicone lid seal & patented lid with ergonomic
   1.5 twist to open and close
- Flush cap saves space during storage
- Printed writing area and graduations
- Recommended for mechanical freezers and vapor phase LN<sub>a</sub> (-196°)



Patented easy-grip faceted lid

Writing Area & Graduations



Traceable Sterility Certificate
Sterile, non-pyrogenic, non-cytotoxic, DNase/RNase-free vials

Slender Cryo Gripper Design

Ste

LOW TEMP	ΗĒ
_106° C	l∎E
-130 C	<b>)</b>
	_





\*Assorted colors: Blue, Red, Green, Purple, Grey

## **Cryogenic & Vaccine Vial Grippers**

#### **The Perfect Grip for Optimal Precision**

These two styles of hand-operated Grippers are designed to pick up cryogenic and crimped 10R-size vaccine/serum vials, respectively.

Both are durably made of reinforced polypropylene which is cold-tolerant and resistant to everyday wear. The designated squeeze target ensures maximum precision and grip with every use.

Each style is available in a pack of 5 assorted colors for easy identification and organization. A combo pack with one of each style is also available.



Ridged Vaccine Gripper Head

f each style is also available.

AUTOCLAVABLE
121° C

PK 5 / CS 100

Available in 5 colors

for easy identification & organization

Item No.	Gripper Type	Color	L x W x H	L x W x H
120908	Cryogenic	Assorted* ■ ■ ■ ■	5.6 x 0.6 x 0.9	14.3 x 1.5 x 2.3
120909	Vaccine	Assorted* ■ ■ ■ ■	6 x 0.8 x 1.4	15.2 x 2 x 3.5

Assorted colors: red, blue green, purple, yellow.

PK 2 / CS 200

Item No.	Gripper Type	Color	L x W x H in	L x W x H
120910	Combination	Blue ■	6 x 0.8 x 1.4 (Vaccine) 5.6 x 0.6 x 0.9 (Cryo)	15.2 x 2 x 3.5 (Vaccine) 14.3 x 1.5 x 2.3 (Cryo)



gripper head

Marked thumb area for optimal ergonomics and control

# **Tube Racks & Storage**

Classic, adaptable, interlocking racks and everything in between



PCR & Microtube

Adaptable & Multi-Tube

Classic Single Racks

## We like to push boundaries







#### 32 Place PCR Rack

Designed for pre- and post-PCR sampling, streamline your PCR workflow with the 32-Place PCR Rack. This versatile rack accommodates 0.2 mL tubes both individually and in strips of 8, as well as 32-place assay plates (8x4)

#### **Space Saving Design**

- Small and compact design is stackable, saving valuable space in the freezer
- Right sized when processing smaller batches
- Lid hinge design allows lid to sit at 180 degrees, clear of sample processing, or can be removed entirely

#### **Reusable Convenience**

- Transparent lid protects samples and allows for quick assessment of contents
- Embossed numbers and letters on the wells provide easy identification and tracking of samples
- Light-weighted construction made durable for reusability

Pack of 10

Item No.	Tubes mL	No of Wells	L x W x H in	L x W x H cm
120889	2.0	32	3.8 x 6.9 x 1.2	9.8 x 17.5 x 3.1

Made from 100% recyclable Polypropylene

Stores down to -80°C

#### **PCR Tube Racks 96-Well**

Perfect for pre- or post-PCR sampling

- Made of polypropylene
- Accommodates individual tubes or strips
- Removable hinged lid lies flat
- Alphanumeric identification for quick location and retrieval of tubes


Pack of 5

Item No.	Color		Item No.	Color
HS23461A	Assorted*		HS23461C	Blue
HS23461B	Natural		HS23461D	Green
		*Assorted colors (Bi	lue, Green, Pink,	Yellow, Orange)

Tubes	Wells	Well Dia.	L x W x H	L x W x H
mL		mm	in	cm
0.2	96	5.2	5 x 3.9 x 1.2	12.7 x 9.9 x 3.0





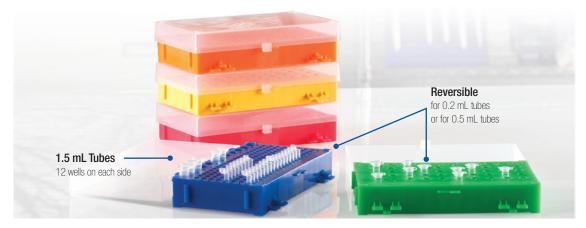

## **PCR Tube Racks Reversible**

Removable hinged lid fits both sides



- Made of polypropylene
- PCR side of the rack holds 8 or 12 x 0.2 mL tube strips or individual tubes
- Reverse side holds 0.5 mL tubes
- Both sides of the rack have 12 wells for 1.5 mL tubes

Pack of 5

lor Tubes Wells Well Dia. L x W x H

ml in (cm)

11011111101	00.0.	mL		mm	in (cm)
HS2344A		0.2	168	5.5	— 8.1 x 5.2 x 2.1
	*	0.5	40	7.6	— (20.6 x 13.1 x 5.3)
		1.5/2.0	12	11.4	(20.0 x 13.1 x 3.3)

\*Assorted colors (Blue, Green, Pink, Yellow, Orange)

## PCR Cooler, 96-Well

**Ice-Free Cooling System for PCR Preparation and Storage** 

Cools industry standard format PCR tubes, strips and single format as well as most standard and low profile 96-Well plates

#### Consistent & Reliable

- Modern design cools faster in freezers and stays cold longer on the bench top
- Maintains the entire 96-Well plate at a temperature of 0 °C for over an hour
- Maintains sample temperature of less than 5°C for up to three hours (with lid off)
- When the temperature ceiling of 6°C has been exceeded, users are signaled by a contrasting color change of the cooler or a single well



for quick location and retrieval of tubes



COLD

Lid doubles as a stand

enables the rack to stay cold longer as there is no conductive heat transfer from surface beneath





provide a secure hold, reducing the chance of dropping the PCR rack, and makes handling and transportation easy, especially when using gloves




WARM TO

#### Two Color Contrasting **Options Available**

purple to pink or dark blue to light blue allowing easy monitoring of samples

				US Patent Pending, (	Community Design No: 0041	01756-0001-0002	Pack of 2
Item No.	Warm Color	Cold Color	Wells	Tube Sizes	LxWxH	LxWxH	
				mL	in	ст	
120727	Pink	Purple <		0.2			
120728	Light Blue	Dark Blue	96		4.3 x 5.8 x 1.9	11 x 14.8 x 4.8	В
120729	Assorted*			Assay Plates			

\*Assorted colors (Light Blue/Dark Blue & Pink/Purple

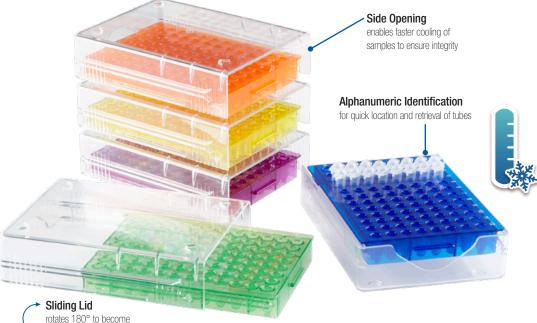

## Low-Temp PCR Rack: 96-Well

Slide, set, work & store PCR samples

Multi-purpose rack converts from compact workstation to lidded storage solution accommodating individual and strip 0.2 mL tubes

#### Versatile & Robust Solution

- Polycarbonate lid serves as cover and elevated rack base - just slide open land pivot lid under rack for compact workstation
- Lid can be used in right or left-handed orientation and is detachable
- Stackable saving valuable bench-top and freezer
- Multiple colors offered for coding options



US Patent Pending, Community Design No: 004101756-0001-0002 Pack of 5


Item No.	Color	Item No.	Color	
120538	Assorted*	120540	Green	
120539	Blue	120541	Purple	

a convenient stand the rack sits on

\*Assorted colors (Purple, Blue, Green, Orange, Yellow)



Tubes mL	Wells	Well Dia.	L x W x H	L x W x H cm
0.2	96	5.2	4.9 x 3.6 x 1.5	12.5 x 9.2 x 3.7

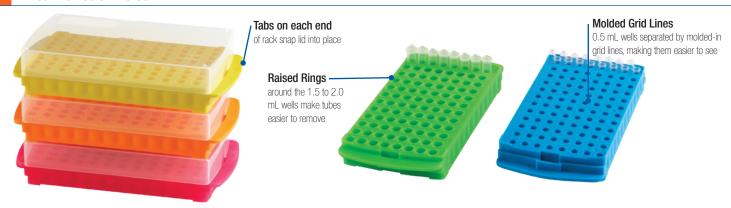

#### **Reversible Tube Racks 96-Well**

Easy-to-grasp tabs for sample handling and storage

96 wells for 0.5 mL tubes on one side and 96 wells for 1.5/2.0 mL tubes on the reverse side

#### **Optimum Function & Identification**

- Wells arranged in standard 8 x 12 array
- Molded in-grid lines separate 0.5 mL wells making them easier to see
- Raised rings around the 1.5 mL/2.0 mL wells make tubes easier to remove
- Alphanumeric identification for quick location and retrieval of tubes

#### Secure lid

- Translucent lid snaps onto rack tabs for secure, visible storage
- Spacious enough to fit most screw-top micro-centrifuge tubes
- Stackable, saving valuable space




Pack of 5

Item No.	Color	Item No.	Color	
HS2345A	Assorted*	HS2345D	Green	
HS2345B	Natural $\Box$	HS2345E	Pink	
HS2345C	Blue	HS2345F	Yellow	

*Assorted colors (Blue, Green, Yello	ow. Orange. Pink)
--------------------------------------	-------------------

Tubes mL	Wells	Well Dia. mm	L x W x H in	L x W x H	
1.5/2.0	— 96	11	9.6 x 4.7 x 2	24.4 x 12 x 5.1	
0.5	90	7.6	9.0 x 4.7 X Z	24.4 X 12 X 3.1	

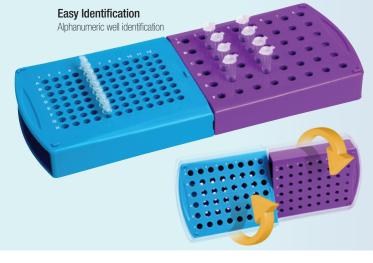

#### Rota-Rack® Duo Tube Rack

**Dilution and concentration studies made easy** 

Transfer solutions between 0.2 mL PCR  $\,$  tubes and 0.5 and 1.5 mL microcentrifuge tubes

- Purple section holds 48 x 0.5 mL and 48 x 1.5/2.0 mL microcentrifuge tubes on the reverse side
- Blue section holds 96 x 0.2 mL tubes and 48 x 1.5/2.0 mL tubes on the reverse side



Item No.	Color		Tubes mL	Wells	Well Dia. mm	L x W x H	L x W x H
	Plue/		0.2	96	5		
HS29041A	Blue/ Purple		0.5	48	8	11.6 x 4.6 x 1.3 29.6	29.6 x 11.8 x 3.4
			1.5/2.0	96	11		

### **Two Tier Microtube Rack 20-Well**


Interlock multiple racks to increase tube storage

#### Improved Experience

- Tiered rack allows for improved staging and visibility for easier, more accurate sample filling and retrieval
- Made of easy to clean polypropylene
- Easily remove tubes using finger slot on front of the rack
- Embossed Index Numbers for quick sample ID and retrieval



						Pack of C
Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
		1.0				
HS4284A	Blue	1.5	20	10	8.3 x 3.5 x 2	21.1 x 9 x 5
		2.0	-			

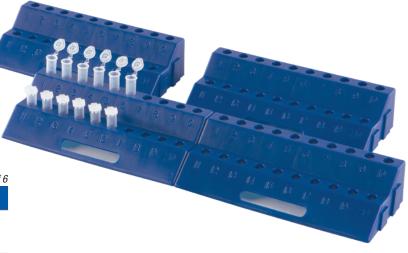

#### **Microtube Racks 80-Well**

Popular fraction-collector-plate style rack

#### Compact 5x16 array makes filling and transfer easy

- Frosted labeling area on bottom of rack
- Embossed index numbers for quick sample ID and retrieval
- Made of polypropylene



Pack of 5

Item No.	Color	Part No.	Color	
HS29025A	Natural	HS29025E	Yellow	
HS29025B	Blue	HS29025F	Orange	
HS29025C	Green <	HS29025G	Assorted*	
HS29025D	Pink	HS29025J	Red	

\*Assorted colors (Blue, Green, Yellow, Orange, Pink)

Tubes mL	Wells	Well Dia.	L x W x H	L x W x H
1.5/2.0	80	11	8.9 x 2.6 x 1	22.5 x 6.7 x 2.8







## **Storage Racks for 80-well Racks**

Strong, clear acrylic construction

These racks help you keep large numbers of microtubes organized on your benchtop.

- Holds up to eight 80-well microtube racks filled with either 1.5 or 2.0 mL microtubes
- Plenty of room to slide the racks in and out

Item No.	Description	Color	LxWxH in	L x W x H cm
HS23263A	Vertical	Clear	6.2 x 9.3 x 9.2	15.7 x 23.6 x 23.4
HS23263B	Horizontal	Clear	12.2 x 9.3 x 4.8	31 x 23.6 x 12.2

## **Microtube Racks**

**Durable, Simple to Use** 

- Made from reinforced nylon offering superior chemical and heat-resistance over most other plastic based racks
- Handles provide easy and safe transportation
- Stable design can be used in water baths
- Alphanumeric identification for quick location and retrieval of tubes




Item No.	Color	Tubes mL	Wells	Well Dia.	Packaging	L x W x H in	L x W x H cm
120241	— Blue ■	1.5/2.0	24	- 11.1	8	5.6 x 3.5 x 2	14.3 x 9 x 5.2
120243	Diue	1.5/2.0	96	- 11.1	4	8.3 x 5.2 x 2.1	21 x 13.1 x 5.4



# Learn more about our convenient multi-functional working rack solutions:





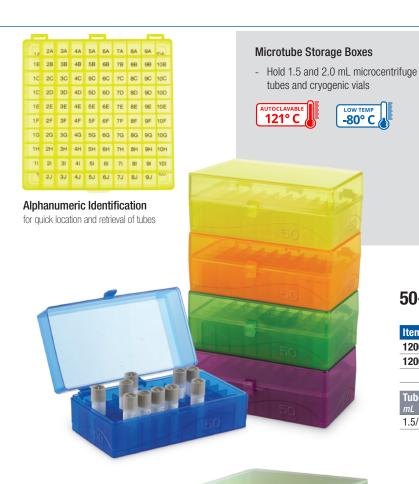

#### Stylish & Identifiable

- Made of translucent, vibrantly colored polypropylene allows for visibility inside box and color coding identification
- Tube capacity is indicated on the outside the box within a distinctive wave design
- Quick visual identification with uniform wells and alphanumeric location identification along with a frosted labeling area

#### Robust

- Easy open and secure lid fitted with either friction fit or snap clasp designs
- Stackable, saving valuable space in the lab

## **50-Well Microtube Storage Boxes**

Pack of 5

Item No.	Color	
120032	Natural	
120033	Assorted*	
		*Assorted colors (Blue, Green, Purple, Yellow, Orange)

Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
1.5/2.0	50	12.3	5.6 x 3.6 x 2.2	14.2 x 9.2 x 5.6

## 81-Well Microtube Storage Boxes

- Black, light blocking version available for light sensitive samples

Pack of 5

Item No.	Color	Item No.	Color	
120034	Natural	120038	Yellow	
120035	Blue	120039	Orange	
120036	Green	120040	Assorted*	
120037	Pink	120041	Black	

\*Assorted colors (Blue, Green, Yellow, Orange, Pink)

Tubes	Wells	Well Dia.	LxWxH	LxWxH
mL		mm		ст
1.5/2.0	81	12.3	5.1 x 5.1 x 1.9	13 x 13 x 4.7



## **100-Well Microtube Storage Boxes**

- Black, light blocking version available for light sensitive samples

Pack of 5

Item No.	Color	Item No.	Color	
120042	Natural	120202	Green	
120043	Black ■	120203	Yellow	
120044	Assorted*	120204	Orange	
120201	Rlug			

\*Assorted colors (Purple, Blue, Green, Yellow, Orange, Pink)

Tubes mL	Wells	Well Dia. mm	L x W x H	L x W x H cm
1.5/2.0	100	12.3	5.6 x 6 x 2.2	14.1 x 15.1 x 5.7

## **Rota-Rack® Tube Racks**

**Get organized with this interlocked rack** 

Each module is interlocked and adjustable with just a twist allowing you to work with multiple sized tubes in the same rack

- Suitable for variety of culture, centrifuge, micro-centrifuge, round, flat or conical bottomed tubes
- Two versions to accommodate small and large tube sizes



Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
	Yellow	50	2	29	8.8 x 4.3 x 4.3	22.4 x 11 x 11
HS29040A	Green	10-15	8	20		
N329040A	Blue	15-20	6	16		
	Purple	10-13	10	13.5		
	Yellow	15	6	16.4		21 x 10 x 10
HS29040B	Green	1.5/2.0	9	11.2	8.2 x 3.9 x 3.9	
по29040В	Blue	0.5	12	8	0.2 X 3.9 X 3.9	
	Purple	0.2	32	6.5		



## **Version available for PCR and Microtubes:**



Rota-Rack, Page 104

## **Cube Tube Rack**

**Interlocking multi-tube racks** 

- Made of polypropylene
- Allows you to work with different sized tubes all in one unit



						T dek of 5
Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
HS29050A	*	50	4	30	4.3 x 4.3 x 4.5	10.8 x 10.8 x 11.3
		15	10	17		
		3-7	12	12.5		
		1.5/2.0	16	11.5	_	
						0/ / // // 0 /

\*Assorted colors (Blue, Green, Pink, Yellow, Orange)





# **Check out our Clinical Version**

Accommodates all brands of 13 and 17 mm round, flat and conical bottom blood collection & culture tubes




# **Transformer Cube Tube Rack**

Low temp storage of large samples. 15 mL & 50 mL tubes

- Quick change over with 2 snap-in inserts that hold 15 and 50 mL tubes
- Rack keyed and numeric wells for tube identification
- Clear Lid included for viewable, stackable storage
- Open space available under insert for extra storage





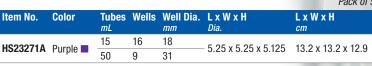



# OneRack® Multi

**Exclusive design, optimal functionality** 

The most flexible tube rack designed for use with a multitude of tube diameters and heights

Holds virtually any brand tube 10-30 mm of varying heights

#### Multiple Configurations for Workflow Efficiency

- Interchangeable tube rack inserts hold multiple tube sizes firmly in place when alternating between sample vessels
- Inserts can be pivoted in either an upper or lower position to securely accommodate the size and height of tube being used

#### Multi-Functional Design

- Open design allows for better filling visibility of tube graduations when filling
- Rack wells have a unique teardrop design that enables better retention of tubes over traditional designs
- Rotate tube rack insert to an angle position to facilitate more comfortable pipetting



Angled position for more comfortable pipetting Alphanumeric grid for quick sample identification and retrieval

**Unique Well Design** Tear drop shape allows

for better retention of

tubes

**Ergonomic handles** 


for ease of transporting the rack




Lower position is ideal height

when working with shorter microcentrifuge tubes


Included: 1 Base + 6 Tube **Rack Inserts** 

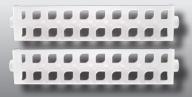

Two 10-13 mm Tube Racks Fits 5 mL, 7.5 mL, 10 mL Tubes

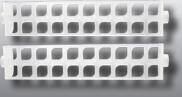

Two 13-17 mm Tube Racks Fits 5 mL, 15 mL Tubes

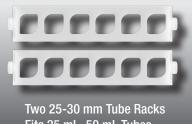

Fits 25 mL, 50 mL Tubes

# OneRack® Multi

US Patent No. 9776188

Item No.	Color	L x W x H in	L x W x H	# Inserts Included	Well Dia./Insert	Wells/Insert
				6 (2 ea/	13 mm	20
120817	Blue/Natural	12.2 x 4.2 x 3	30.9 x 10.7 x 7.7	Diameter)	17 mm	20
				Diameter)	30 mm	6
				6 (2 ea/	13 mm	20
120818	Purple/Natural	12.2 x 4.2 x 3	30.9 x 10.7 x 7.7	•	17 mm	20
				Diameter)	30 mm	6

# Adapt-a-Rack® Adaptable Multi-tube Rack

Each well has fingers that flex that adapt to accommodate a variety of tubes sizes while holding them in a firm upright position

- Made of POM
- Racks can be linked together for additional capacity
- Open sides allow for label viewing
- Flat base to hold tubes upright
- Available in solid or a combination of colors for easy coding identification







Pack of 2

Item No.	Rack #1 Color	Rack #2 Color	
120184	Blue/Green	Green/Blue	
120185	Blue/Yellow	Yellow/Blue	<u> </u>
120186	Pink/White	White/Pink	
120187	Blue/Blue	Blue/Blue	
120188	White/White	White/White	
120189	Pink/Pink	Pink/Pink	
120190	Yellow/Yellow	Yellow/Yellow	
120191	Green/Green	Green/Green	

Tubes	Wells	Well Dia.	L x W x H	L x W x H
mL		mm	in	cm
7–50	4	12-30	7.1 x 2.2 x 3	18.1 x 5.6 x 7.6

# Looking for racks for tube rollers?:



Tube Roller Racks, Page 70

# **One-Well Connecting Tube Racks**

Form a line or use as a single, stand-alone polypropylene rack

- Open Design enables easy identification of samples
- Molded circular base ring keeps conical or round-bottom tubes in place

Pack of 5

Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H
HS23052	*	50	1	29	3.2 x 2.8 x 2.6	8.3 x 7 x 6.7

\*Assorted colors (Blue, Green, Yellow, Orange, Pink)



# One-Well Friction-Fit Tube Rack

- Made from polystyrene
- Open Design enables easy identification of samples
- Holds conical and round bottom tubes
- Slides smoothly across the benchtop

# **One-Well Wire Tube Rack**

- HDPE-coated steel wire
- Holds conical and round bottom tubes
- Open Design enables easy identification of samples
- Base ring keeps tubes in place





#### One-Well Friction-Fit Tube Rack

Pack of 5

one-wen i netion-i it lube nack								
Item No.	Color		Tubes mL	Wells	Well Dia.	W x H in	W x H cm	
HS23053A	Blue		50	1	29	3.6 X 1.2	9.2 x 3	_

#### **One-Well Wire Tube Rack**

Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
HS23051	Yellow	50	1	29	3.5 x 2.6 x 2.6	9 x 6.7 x 6.7





# **Snap-Together Conical Tube Rack**

**Space Saver - shipped flat and easy-to-assemble** 

#### Conveniently holds 15 and 50 mL tubes simultaneously

- ABS racks snap together securely and are designed for stable use in water baths
- Alphanumeric identification for quick location and retrieval

Item No.	Color	Tubes mL	Wells	Well Dia.	LxWxH in (cm)
HS24306	Orongo	15	30	17	6.7 x 8.1 x 2.3
	Orange _	50	20	28	(17.1 x 20.6 x 5.9)





Item No.	Color	Tubes mL	Wells	Well Dia. mm	L x W x H in	L x W x H cm
HS24330A	White/Purple □ ■	15-50	12	16–30	7.6 x 5.5 x 2	19.4 x 14 x 5.1

# **Pop-Up Tube Racks**

Collapse for compact storage and "pop-up" for quick use

#### Designed to save space!

- Easy to open, holds 15 and 50 mL tubes
- Easy to close, stackable for compact storage
- Alphanumeric well identification for quick location and retrieval of samples

Pack of 2

Item No.	Color	Tubes mL	Wells	Well Dia. mm	
HS24320B	Blue	15. 50	21, 12	17. 30	
HS24320P	Pink	15, 50	21, 12	17, 30	
HS24318	Purple	15	45	17	
HS24319	Green	50	18	30	

L x W x H	L x W x H	L x W x H	L x W x H
(Expanded) in	(Collapsed) in	(Expanded) cm	(Collapsed) cm
10 x 5.4 x 2.8	10 x 5.4 x 0.8	25.5 x 13.7 x 7.2	25.5 x 13.7 x 2.1



Rack easily "Pops" open for easy use







# **Fold and Snap Tube Rack**

#### **Economical durable rack**


- Polypropylene racks ship flat in one piece
- Easy-to-assemble, just fold and snap together
- Stable design can be used in water baths
- Open design allows samples to be incubated at the same temperature
- Color code for easy identification










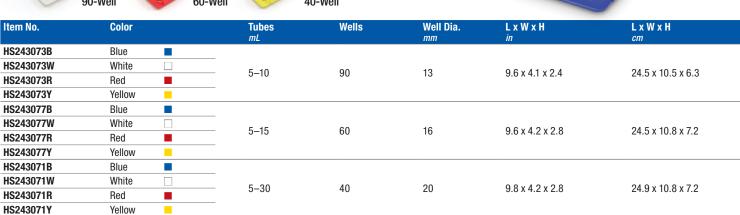

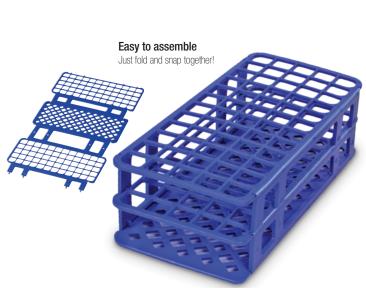



#### Available in two materials:

## Fiberglass-Reinforced Polyoxymethylene (POM)

- Available in half and full rack sizes for 5 different tube diameters
- Offers the highest chemical and heat resistance (-40 to +121°C)

#### Polypropylene

- Available in full rack sizes for 2 different tube diameters
- Offers good chemical and heat resistance (0 to 110°C)

#### OneRack® Full-Size POM





Ulichack	i uii-Size	UIV		_		
Item No.	Color		Tubes mL	Wells	Well Dia.	LxWxH in (cm)
HS27511A	White					
HS27511B	Blue					7.9 x 4 x 2.3
HS27511C	Green		5-10	72	13	(20.2 x 10.4
HS27511D	Magenta					x 5.9)
120205	Yellow					
HS27512A	White					
HS27512B	Blue					9.8 x 5 x 2.7
HS27512C	Green		5–10	72	16	
HS27512D	Magenta					(25 x 12.7 x 7)
120206	Yellow					
HS27513A	White					9.9 x 3.9 x 3.2
HS27513B	Blue		10–18	40	20	(25.2 x 10 x
HS27513C	Green		10-10	40	20	8.3)
HS27513D	Magenta					0.3)
HS27514A	White					
HS27514B	Blue		10–18	40	25	11.8 x 4.7 x 3.6
HS27514C	Green		10-10	40	23	(30 x 12 x 9.2)
HS27514D	Magenta					
HS27515A	White					
HS27515B	Blue					11.1 x 4.3 x 3.3
HS27515C	Green		15-50	24	30	(28.2 x 11 x
HS27515D	Magenta					8.5)
120207	Yellow					

# OneRack® Half-Size POM





Ulichack	Hall-Size	I ON	1		<u></u>	
Item No.	Color		Tubes mL	Wells	Well Dia.	LxWxH in (cm)
HS27501A	White					
HS27501B	Blue		5–10	36	13	4.1 x 4.1 x 2.3
HS27501C	Green		J=10	30	13	(10.4 x 10.4 x 5.9)
HS27501D	Magenta					
HS27502A	White					
HS27502B	Blue			36	16	5 x 5 x 2.7
HS27502C	Green		5-10			(12.7 x 12.7 x 7)
HS27502D	Magenta					(12.7 x 12.7 x 1)
120399	Yellow					
HS27503A	White					4 x 5 x 3.2
HS27503B	Blue		10–18	20	20	(10 x 12.7 x 8.3)
HS27503C	Green					(10 x 12.7 x 0.3)
HS27504A	White		10–18	16	25	4.7 x 4.8 x 3.6 (12 x 12.2 x 9.2)
HS27505A	White					
HS27505B	Blue		15–50	9	30	4.3 x 4.3 x 3.3
HS27505C	Green		15-50	J	30	(11 x 11 x 8.5)
HS27505D	Magenta					

OneRack® Full-Size Polypropylene

Item No.	Color	Tubes mL	Wells	Well Dia.	LxWxH in (cm)
HS27551A	White				
HS27551B	Blue	5–10	70	13	7.9 x 4 x 5.9
HS27551C	Green	5-10	12	13	(20.2 x 10.4 x 5.9)
HS27551D	Magenta	-			
HS27552A	White				
HS27552B	Blue	- 5–10	70	16	9.8 x 5 x 2.7
HS27552C	Green	5-10	12	10	(25 x 12.7 x 7)
HS27552D	Magenta	-			


# **Coated Wire Racks**

# **Chemically resistant HDPE Coated Wire Racks**

- Holds all brands and types of conical, round or flat bottom, glass or plastic, tubes and/or bottles
- Square wells and three deck levels help keep tubes secure





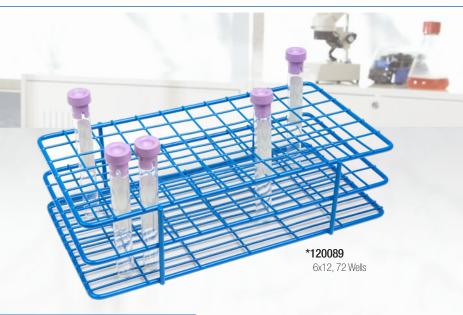

Fits 10-13 mm Tubes

Item No.	Color	Tube Array	Wells	Well Dia. mm	L x W x H in	L x W x H cm
120089*	Blue	6x12	72	13	8.1 x 4.1 x 2.6	20.5 x 10.5 x 6.5
120759	Blue	9x12	108	13	8 x 6.4 x 2.6	20.4 x 16.1 x 6.5



Fits 13-16 mm Tubes

Item No.	Color	Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
120760	Blue	4x10	40	16	7.9 x 3.3 x 2.4	20 x 8.3 x 5.9
120761*	Blue	4x12	48	16	9.4 x 3.3 x 2.4	23.7 x 12.7 x 5.9
HS23072	Blue	6x12	72	16	9.5 x 5 x 2.6	24.1 x 18.8 x 6.5
120763	Blue	9x12	108	16	9.3 x 5 x 2.6	23.7 x 12.7 x 6.5
HS23108	Red	9x12	108	16	9.3 x 5 x 2.6	23.7 x 12.7 x 6.5

Fits 16-18 mm Tubes - 1 5mL


Item No.	Color		Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
HS23215*	Blue		3x5	15	20	5 x 3 x 3	12.7 x 7.6 x 7.6

Fits 16-20 mm Tubes - 15 mL

Item No.	Color	Tube Array	Tube Array Wells		LxWxH in	L x W x H cm
120764	Blue	4x10	40	20	9.5 x 3.9 x 3.1	24 x 9 x 8
120765	Blue	8x10	80	20	9.5 x 7.6 x 3.1	24 x 19 x 8



\*HS23215 3x5, 15 Wells

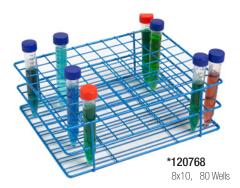

#### Fits 20-24 mm Tubes

Item No.	Color	Tube Array	Wells	Well Dia.	LxWxH in	L x W x H cm
120766	Blue	4x10	40	24	10.2 x 4.2 x 3.2	26 x10.7 x 8

#### Fits 20-25 mm Tubes

Item No.	Color	Array	Wells	mm	LxWxH in	L x W x H cm
120769	Blue	3x8	24	25	10.7 x 4.1 x 3.2	27.3 x 10.5 x 8
120767	Blue	4x10	40	25	11.4 x 4.7 x 3.1	29 x 11.8 x 8
120768*	Blue	8x10	80	25	11.4 x 9.1 x3.1	29 x 23 x 8

# **Coated Wire Racks**

# **Chemically resistant HDPE Coated Wire Racks**



Fits 25-30 mm Tubes - 50 mL

Item No.	Color	Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
HS232461*	Blue	2x4	8	30	5.8 x 3.2 x 3.5	14.7 x 8 x 8.9
120770	Blue	6x8	48	30	10.9 x 8.3 x 3.9	27.7 x 20.9 x 9.8

Fits 30-40 mm Tubes

Item No.	Color		Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
120771*	Blue		3x8	24	40	14.1 x 5.5 x 3.9	35.7 x 13.9 x 9.8
120772	Blue		6x8	48	40	14 x 10.6 x 3.8	35 x 26 x 9.8





Fits 23-25 mm Bottles - Bijou Bottles

Item No.	Color	Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
120085	Blue	2x6	12	25	6.4 x 2.5 x 1.3	16.4 x 6.4 x 3.2
120086*	Blue	5x10	50	25	10.7 x 5.7 x 1.3	27.1 x 14.5 x 3.2

Fits 30-33 mm Bottles - Universal & McCartney Bottles

Item No.	Color	Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
120087*	Blue	2x6	12	33	8.3 x 3.1 x 2.4	21.2 x 8 x 6
120088	Blue	5x10	50	33	8.7 x 3.2 x 2.4	21.8 x 8.2 x 6





Fits 55-58 mm Bottles - Urine Containers

Item No.	Color	Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
120091*	Blue	4x6	24	58	15.5 x 10.3 x 2	39.3 x 26.2 x 5

## Fits 55-60 mm Bottles - Centrifuge Bottles

Item No.	Color	Tube Array	Wells	Well Dia.	L x W x H in	L x W x H cm
HS23224*	Blue	2x3	6	60	8.1 x 5.3 x 4.1	20.5 x 13.5 x 10.4



# Mega Rack®

## High-capacity storage for your tubes

Fully assembled, polypropylene mega racks are great for high throughput labs requiring large sample storage

- Sturdy, material for easy and secure transport of blood collection and other tubes
- Alphanumeric identification for quick location and retrieval of tubes
- Place different test groupings into one rack eliminating the need for multiple racks

#### Mega Racks® 10-13 mm - Holds 5-7 mL sample tubes

Item No.	Color	Wells	Well Dia. mm	LxWxH in	L x W x H
HS24322A HS24322E	Blue White	- 216	13	13.3 x 9.3 x 3.6	33.7 x 23.5 x 9.1
120112 120113	White □ Blue ■	<b>- 432</b>	13	26.5 x 9.3 x 3.5	67.3 x 23.5 x 9.0

AUTOCLAVABLE 121° C

#### Mega Racks® 13-16 mm - Holds 10 mL sample tubes

Item No.	Color	Wells	Well Dia. mm	L x W x H in	L x W x H
HS24323A HS24323E	Blue  White	- <b>120</b>	16	13.3 x 10.5 x 3.6	33.7 x 26.7 x 9.1
120114	White	- 240	16	26.5 x 10.5 x 3.5	67.3 x 26.7 x 9.0
120115	Blue	240	10	20.0 X 10.0 X 0.0	07.0 X 20.7 X 0.0

#### Mega Racks® 18-20 mm - Holds 15 mL conical/swab tubes and 20 mm round bottom tubes

Item No.	Color	Wells	Well Dia.	L x W x H in	L x W x H
120834	White	- 120	20	14.9 x 11.9 x 3.6	37.7 x 30.1 x 9.1
120835	Blue	120	20	14.9 % 11.9 % 3.0	37.7 x 30.1 x 9.1

# **Peg Tube Racks**

Strong, tough, versatile racks for bench use or tube storage

#### Fiberglass-reinforced polypropylene racks are ideal for bench work or tube storage

- Spacious peg design allows easy viewing of tube contents
- Pegs remain straight even after repeated autoclaving
- Embossed alphanumeric identification for quick location and retrieval of tubes

# Unique tear drop well design

allows for better retention of tubes



Item No.	Color	Tubes mL	Wells	Well Dia.	LxWxH in	L x W x H cm
HS24311A	White	10	50	14–17	7.4 x 4.1 x 2.7	18.8 x 10.4 x 6.9
HS24311B	Blue	10	30	14-17	1.4 X 4.1 X 2.1	10.0 X 10.4 X 0.9
HS24312A	White	7	96	10–13	3.2 x 5 x 2.7	8.2 x 12.6 x 6.9
HS24312B	Blue	1	90	10-13	3.2 X 3 X 2.1	0.2 % 12.0 % 0.9

Pack of 2

# View our complete line of Freezer Storage Boxes



# **Cold Products & Storage**

From cool-looking ice containers to freezer necessities, our innovative range of cold products continues to expand



Cooling Containers

Low Temp Racks & Storage

Freezer Boxes, Racks & Accessories





**Durable containers for chilling temperature-sensitive samples**Suitable for low temperature use. Ideal for ice, ice slurries and CO<sub>2</sub> pellets

- Made from polyurethane with metal inserts for stability

Item No.

HS28714K

HS28721G

- Close-fit lid keeps evaporation to a minimum and features a comfortable recessed knob
- Tapered design for easy stacking
- Precise pour spout and large handles




**Lid nests under bucket** for convenience and to save space

Color

Black

Green



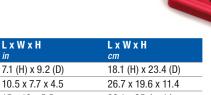

# HS28724B Blue Pan 4 15 x 10 x 5.5 38.1 x 25.4 x 14 HS28729R Red Pan 9 20.5 x 11.5 x 6.4 52.1 x 29.2 x 16.3

Style

Bucket

Pan

# Learn more about our open rack systems that work perfectly with our Cooling Containers

Capacity

4

1



Fold and Snap Tube Rack, Page 114

OneRack® Tube Rack, Page 115

# Floating Microtube Rack For use on bench and in ice and

water baths

- Made of polypropylene
- Attach legs and rack stands upright for use on bench top

Color

- Detach legs and samples can be quickly removed by pushing down onto any surface




Item No.

HS2132A



Tubes

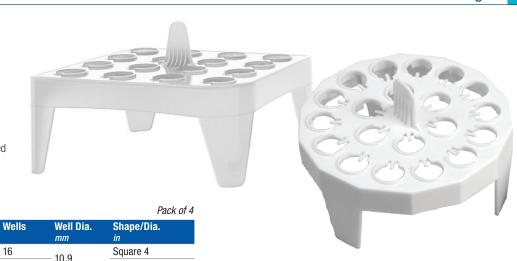

2132 <i>A</i> 2134 <i>A</i>		1.	5/2.0	16 20	10.9	Square 4 Round 3.8				
	Electing	j Tube Ra	acke							
		and flexibl								1 1 2 2
	<ul><li>Brightly color</li><li>polyethylen</li><li>Racks inclu</li></ul>	ating tubes in vored foam rack e that resists wide detachable	ks made of hi vater absorpt	gh-density ion						
	- Temperatur	out of baths re Range 0° to	70° C					1		
	O° C						1111111		4	
						To both tracks	Carlo Martin			
									100	
						- 4			III	

Pack of 5

Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
HS2165A	Green	50	4	29	6.4 x 5.3 x 2.5	16.5 x 13.7 x 6.5
HS2165B	Blue	15	8	16	6.7 x 2.2 x 2.5	17.2 x 5.7 x 6.5
HS2165C	Yellow	1.5/2.0	24	9.3	6 x 3.2 x 2.7	15 x 8.3 x 7
HS2166	Blue	0.2/0.5/1.5/2.0	18	9.3	3.2	8.2 dia.





**Alphanumeric Well Identification** for guick location and retrieval of tubes

#### Lid doubles as a stand

enables the rack to stay cold longer as there is no conductive heat transfer from surface beneath



Vented Lid

enables faster cooling of samples while in freezer



# PCR Cooler, 96-Well

**Ice-Free Cooling System for PCR Preparation and Storage** 

Cools industry standard format PCR tubes, strips and single format as well as most standard and low profile 96-Well plates

#### Consistent & Reliable

- Modern design cools faster in freezers and stays cold longer on the bench top
- Maintains the entire 96-Well plate at a temperature of 0  $^{\circ}\text{C}$  for over an hour
- Maintains sample temperature of less than 5°C for up to three hours (with lid off)
- When the temperature ceiling of 6°C has been exceeded, users are signaled by a contrasting color change of the cooler or a single well



# Two Color Contrasting Options Available

purple to pink or dark blue to light blue allowing easy monitoring of samples

**Ergonomic Textured Finger Grips** provide a secure hold, reducing the chance of

dropping the PCR rack, and makes handling and transportation easy, especially when using gloves







US Patent Pending, Community Design No: 004101756-0001-0002

Pack of 2

Item No.	Warm Color	Cold Color	Wells	Tube Sizes	L x W x H	L x W x H
120727	Pink	Purple <		0.0		
120728	Light Blue	Dark Blue	96	0.2	4.3 x 5.8 x 1.9	11 x 14.8 x 4.8
120729	Assorted*			Assay Plates		


COLD

# Need to keep your reagents cool while working?

Check out our ice-free Reagent Reservoir Chillers with temperature color changing feature



Reservoir Chiller, Page 83

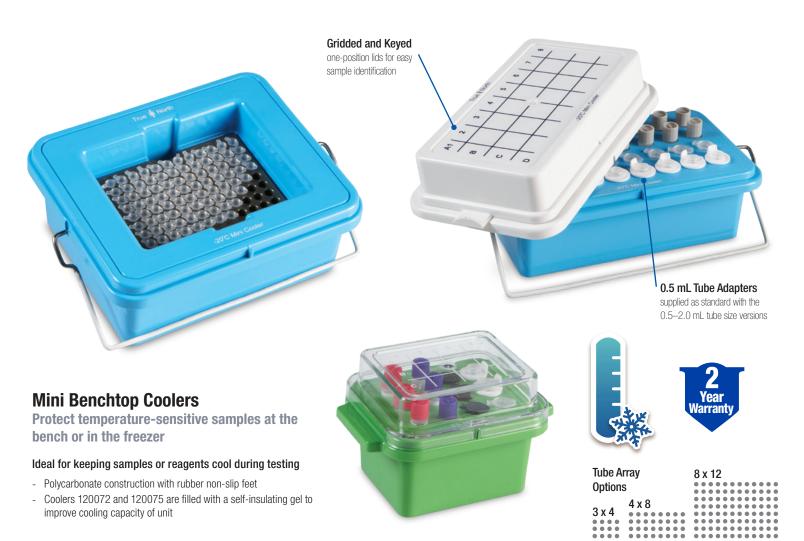

Item No.	Color	Tubes mL	Lid	Tube Array	Wells	L x W x H in	L x W x H cm	Temp.	Temp. Retention
120070	Green		Clear	3 x 4	12	5.7 x 4.1 x 3.9	14.5 x 10.5 x 10	0°C	1°C for 3 hr
120071	dieen	0.5-2.0	Clear	4 x 8	32	9.6 x 4.1 x 3.9	24.5 x 10.5 x 10	0°C	1°C for 2 hr
120072	_		Gel Filled	4 x 8	32	9.6 x 4.1 x 3.9	24.5 x 10.5 x 10	0°C	1°C for 5 hr
120073	Dlue		Clear	3 x 4	12	5.7 x 4.1 x 3.9	14.5 x 10.5 x 10	-20°C	-15°C for 1 hr
120074	Blue	0.5-2.0	Clear	4 x 8	32	9.6 x 4.1 x 3.9	24.5 x 10.5 x 10	-20°C	-15°C for 1 hr
120075	_		Gel Filled	4 x 8	32	9.6 x 4.1 x 3.9	24.5 x 10.5 x 10	-20°C	-15°C for 1 hr
120076	Blue	0.2	Clear	8 x 12	96	8.9 x 6.5 x 4.3	22.5 x 16.5 x 11	-20°C	-15°C for 2 hr



# TracRack® Tube Rack

Take what you need and store the rest

Comprised of a storage box and eight removable tube segments allowing you to easily take what you need - an entire segment or individual sample - and store the rest!

#### **Convenient Rack System**

- Eight polypropylene segments hold 1.5–2.0 mL cryo and microcentrifuge tubes (64 tubes)
- Wells feature a center depression stabilizing conical tubes and a cryogenic vial lock for one-handed cap removal
- Printed index letters on each segment make for quick sample ID and retrieval






# TracRack® Outer Box

- Made from thin film polypropylene material
- Moisture and mold resistant
- 2 in (5.3 cm) boxes fit in standard freezer racks
- Not recommended for autoclaving



#### TracRack® Tube Rack

Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H	L x W x H
120356	Purple	1.5/2.0	64	12.3	5.1 x 5.1 x 1.4	12.9 x 12.9 x 3.5

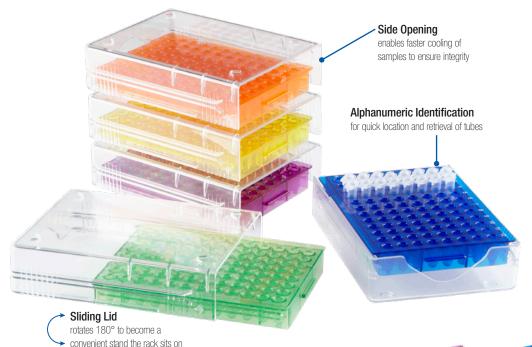

# Low-Temp PCR Rack: 96-Well

Slide, set, work & store PCR samples

Multi-purpose rack converts from compact workstation to lidded storage solution accommodating individual and strip 0.2 mL tubes

#### Versatile & Robust Solution

- Polycarbonate lid serves as cover and elevated rack base - just slide open land pivot lid under rack for compact workstation
- Lid can be used in right or left-handed orientation and is detachable
- Stackable saving valuable bench-top and freezer space
- Multiple colors offered for coding options



US Patent Pending, Community Design No: 004101756-0001-0002

Item No.	Color	Item No.	Color	
120538	Assorted*	120540	Green	
120539	Blue	120541	Purple	

\*Assorted colors (Purple, Blue, Green, Orange, Yellow)

Pack of 5


0.2

96



4.9 x 3.6 x 1.5

12.5 x 9.2 x 3.7

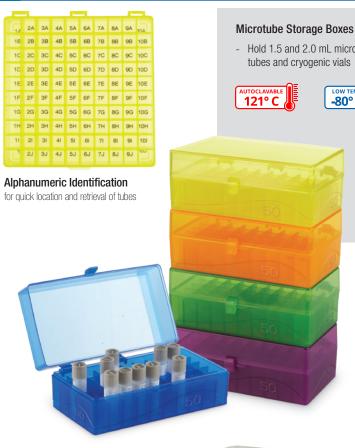

#### Stylish & Identifiable

- Hold 1.5 and 2.0 mL microcentrifuge tubes and cryogenic vials

- Made of translucent, vibrantly colored polypropylene allows for visibility inside box and color coding identification
- Tube capacity is indicated on the outside the box within a distinctive wave design
- Quick visual identification with uniform wells and alphanumeric location identification along with a frosted labeling area

#### Robust

- Easy open and secure lid fitted with either friction fit or snap clasp designs
- Stackable, saving valuable space in the lab

# **50-Well Microtube Storage Boxes**

Pack of 5

Item No.	Color	
120032	Natural	
120033	Assorted*	

\*Assorted colors (Blue, Green, Purple, Yellow, Orange)

Tubes mL	Wells	Well Dia. mm	LxWxH in	L x W x H cm	
1.5/2.0	50	12.3	5.6 x 3.6 x 2.2	14.2 x 9.2 x 5.6	

# **81-Well Microtube Storage Boxes**

- Black, light blocking version available for light sensitive samples

Pack of 5

Item No.	Color	Item No.	Color	
120034	Natural	120038	Yellow	
120035	Blue	120039	Orange	
120036	Green	120040	Assorted*	
120037	Pink	120041	Black	

\*Assorted colors (Blue, Green, Yellow, Orange, Pink)

Tubes	Wells	Well Dia.	L x W x H	L x W x H
mL		mm	in	cm
1.5/2.0	81	12.3	5.1 x 5.1 x 1.9	13 x 13 x 4.7



- Black, light blocking version available for light sensitive samples


Pack of 5

Item No.	Color	Item No.	Color	
120042	Natural	120202	Green	
120043	Black	120203	Yellow	
120044	Assorted* ■■■	120204	Orange	
120201	Blue			

\*Assorted colors (Purple, Blue, Green, Yellow, Orange, Pink)

Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H
1.5/2.0	100	12.3	5.6 x 6 x 2.2	14.1 x 15.1 x 5.7



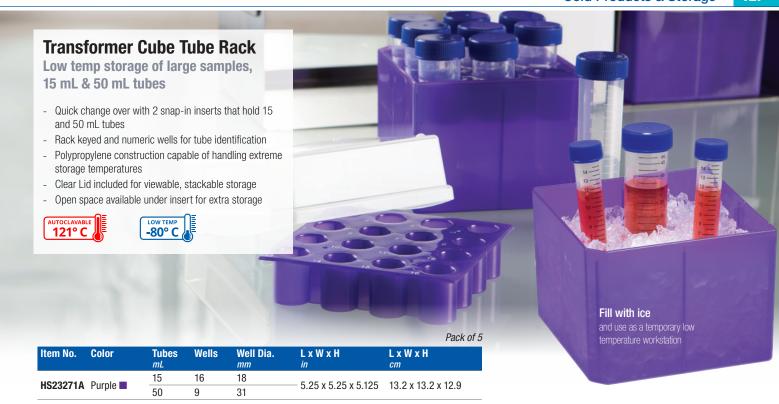



# **Foam Tube Racks**

High storage, minimal footprint

- Two-piece racks are constructed of expanded polystyrene (EPS) foam
- Lightweight microtube racks are ideal for long-term studies or tube storage  $\,$
- Extra spacing between wells allows for easy gripping with fingertips
- Racks stack securely with nesting features on lid and base
- Must use tubes capable of being stored at the required temperature of the user



Pack of 100

Item No.	Color	Tubes mL	Wells	Well Dia.	LxWxH in	L x W x H cm
HS4283P	White 🗌	1.5/2.0	50	11.5	8.3 x 4.3 x 2.7	21 x 11 x 7.1

# Store microslides in low temperatures with our Low Temperature Microscope Slide Boxes



Low Temperature Microscope Slide Boxes, Page 142

# **Cryogenic Vial Rack**

Safely store and transport cryogenic vials

- Made from polycarbonate
- Well tabs allow one-handed opening and closing of vials that have a slot feature
- Embossed index numbers and squared, upper left corner provide easy orientation
- Drainage holes in each well aid in cleaning
- Alphanumeric identification for quick location and retrieval of tubes





Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H	L x W x H
		IIIL		,,,,,,	""	UIII
HS23229B	Blue	2-5	50	12.9	7.8 x 4.1 x 0.9	19.9 x 10.5 x 2.8

provide easy orientation



make handling with gloves easier

# **Cryogenic Vial Canes**

Organize cryogenic samples in dewars

Aluminum canes hold vials for freezer storage or immersion in LN, dewars

- Use in liquid and vapor phase; use approved vials only
- Canes bend to fit all major cryogenic vial brands



Pack of 12

Item No.	Tubes mL	No. of Vials	Length in	Length cm
HS23205	1.5/2.0	5	11.4	29
HS23206	1.5/2.0	6	11.8	30

# **PVC Cryogenic Cane Sleeves**

For extra security and sample safety

Clear plastic sleeves protect vials in aluminum canes during handling and storage

- Safe for cryogenic freezing
- More durable than cardboard sleeves
- Transparency allows sample identification without removing sleeve



Pack of 100

Item No.	Diameter in	Length in	Diameter cm	Length cm
HS23210	0.6	10.8	1.5	27.5

# Don't miss our line of Cryogenic Vials

Our innovative flush design saves space in your freezer





Cryogenic Vials, Page 99

Pack of 12

# **Cryogenic Vial Boxes & Partitions**

For storage of small samples and tubes

- Made from cardboard
- Colors allow for easy identification and sorting purposes
- All box footprints are 13.3 x 13.3 cm, select from two different heights



## **Cardboard Cryogenic Vial Boxes**

Color Tubes LxWxH Item No. LxWxH White 1.5/2.0 5.2 x 5.2 x 2 13.3 x 13.3 x 5 HS2860A HS2860B White 1.5/2.0 5.2 x 5.2 x 3 13.3 x 13.3 x 7.5 HS2860C Assorted Blue 1.5/2.0 5.2 x 5.2 x 2 13.3 x 13.3 x 5 HS2860CB HS2860CR Red



<b>Partitions</b>	Partitions m	Pack of 12		
Item No.	Array	Tubes mL	Wells	Well Dia. mm
HS2860E	10 x 10	1.5/2.0	100	11.5
HS2860H	9 x 9	1.5/2.0	81	12
HS2860F	8 x 8	1.5/2.0	64	15
HS2860G	7 x 7	1.5/2.0	49	17.5

# **Microtube Storage Boxes**

Safe microtube storage

Low profile boxes are available in two sizes to accommodate standard 0.5 mL and 0.2 mL microtubes and vessels of a similar size  $\,$ 

- Partitions included in the boxes
- Will fit into most standard freezer racking



Pack of 10

Item No.	Color	Tubes mL	Wells	Well Dia. mm	LxWxH in	L x W x H cm
120220	White	0.2	169	9.5	5.2 x 5.2 x 1.1	13.3 x 13.3 x 2.9
120221	White	□ 0.5	100	12	5.2 x 5.2 x 1.6	13.3 x 13.3 x 4.1

# of 10

# 15/50 mL Tube Storage Boxes & Partitions

- Write-on surface for descriptions or self-numbering systems
- Partitions available in two sizes





Pack of 10

Item No. Color		L x W x H in	L x W x H cm
HS2860M	White $\square$	5.7 x 5.7 x 4.8	14.5 x 14.5 x 12.2

<b>Partitions</b>	Partition	ns must be ordered	separately	Pack of 10
Item No.	Array	Tubes mL	Wells	Well Dia. mm
HS2860T	6 x 6	15	36	23
HS2860P	4 x 4	50	16	35

# **Cryogenic Vial Tube Storage Box**

- Two-piece set-up boxes with finger holds
- Boxes made from paper coated cardboard
- 5 x 5 format for tubes up to 14 mm in diameter
- Partitions are included in each box





Case of 200

Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H
120519	Blue	2.0	25	14	3.4 x 3.4 x 1.5	87x87x38



120392 Blue ■ 13 30 19.8

120366 White □ 50 16 31.6

120367 Blue ■ 50 16 31.6

Thin Polypropylene film 1.5–2.0 mL

12036° C

Item No.	Color	Tubes mL	Wells	Well Dia.	L x W x H in	L x W x H cm
120370	White					
120371	Blue	1.5/2.0	81	12.1	5 x 4.9 x 1.9	12.5 x 12.4 x 4.9
120372	It Green					

Pack of 10

LxWxH

13.3 x 13 x 3

13.3 x 13 x 4.2

13.3 x 13 x 5.1

13.3 x 13 x 7.5

LxWxH

14.7 x 14.5 x 12

14.7 x 14.5 x 12

18.8 x 19.8 x 4.1

Pack of 10

Pack of 10



# **Arctic Squares® Polycarbonate Cryostorage Boxes**

Safely store cryogenic vials

Our polycarbonate boxes are designed with unique features to help orient, identify and access vials with ease

- Alphanumeric labeling on lid or base (depending on the unit) for quick location and retrieval of tubes
- Liquid nitrogen vent and drainage holes
- Magnet on top of forceps enables easy storing on metal surfaces
- Stackable, saves valuable space in the lab






US Patent No.: D529,622

EP Registered Design No.: 00551429-0001, 00551429-0002, 00551429-0003

						U	U			
Item No.	Color	Tube Array	Tubes	Wells	Well Dia.	Lid	Packaging	LxWxH	LxWxH	Forceps
			mL		mm			in	ст	
HS23325A	Ruby Red	5 x 5	2.0	25	12.5	Slide On	8	3 x 3 x 2	7.6 x 7.6 x 5.3	
HS23381B	Assorted*	9 x 9	2.0	81	12.5	Hinged	4	5.2 x 5.2 x 2	13.3 x 13.3 x 5.3	✓
HC33300D	Arctic Blue	10 v 10	2 በ	100	12.5	Hinged	1	52 x 5 2 x 2	13 3 v 13 3 v 5 3	<b>✓</b>

\*Assorted colors (Ruby Red, Amethyst Purple, Arctic Blue, Emerald Green )

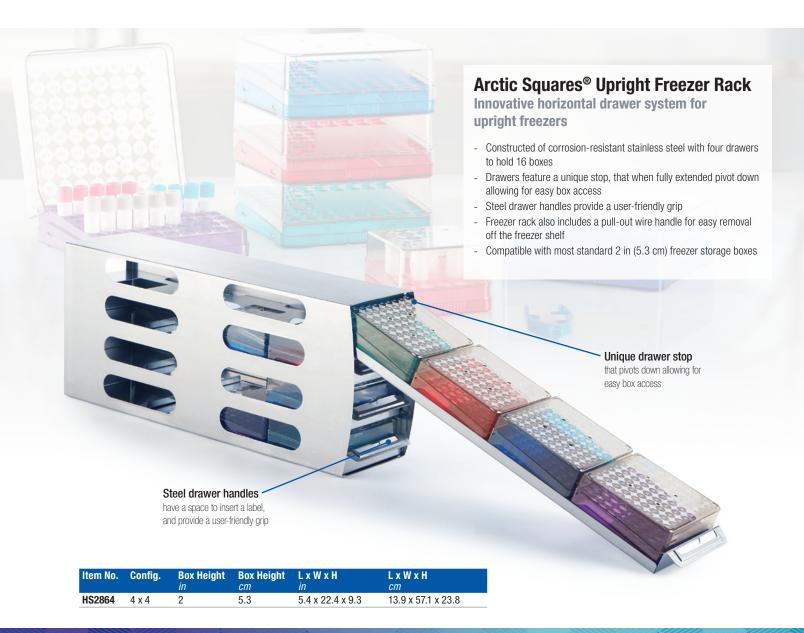



Item No.	Config.	Box Size	Box Size	LxWxH	LxWxH
		in	cm	in	cm
HS28622	2 x 3	6	15.2	18 x 6.2 x 10.6	45.9 x 15.9 x 27.1

15/50 mL Tube Storage Box, Page 129

Transformer Cube Tube Rack, Page 127

# Get a Grip with Cryogenic & Vaccine Vial Grippers

easy use

handles stay securely in place for

The Perfect Grip for Optimal Precision



Cryogenic & Vaccine Vial Grippers, Page 99

#### **Chest Freezer Racks**

Item No.	Config.	Box Size in	Box Size cm	L x W x H in	LxWxH cm
HS2862DB	9 x 1	3	7.6	5.5 x 5.5 x 29	14.1 x 14.1 x 73

# **Upright Freezer Racks: 2 inch Boxes**

Item No.	Config.	Box Size in	Box Size cm	L x W x H in	L x W x H cm
HS2862HA	4 x 3			16.4 x 5.5 x 8.9	41.8 x 14.1 x 22.5
HS2862H	4 x 4	2	5.3	21.8 x 5.5 x 8.9	55.6 x 14.1 x 22.5
HS2862HB	5 x 4			22 x 5.5 x 11	55.6 x 14.1 x 28

# Upright Freezer Racks: 3 inch Boxes

Item No.	Config.	Box Size in	Box Size cm	L x W x H in	L x W x H cm
HS2862KB	4 x 4	2	7.6	22 x 5.5 x 12.8	55.6 x 14.1 x 32.5
HS2862KC	5 x 4	3	7.6	22 x 5.5 x 16	55.6 x 14.1 x 40.6



Microtube Storage Boxes, Page126

#### **Upright Freezer Rack Configurations**











# Handles for easy transfer

Robust mounting brackets assure handles stay securely in place for easy use





HS2862AB



# **Tips for Keeping your Laboratory Freezer Running Efficiently**

#### **Monitor Temperatures**

Installing an alarmed thermometer monitoring system, like Heathrow Scientific's Freezer Thermometer, can alert you early to any fluctuations in temperatures that could indicate the need for cleaning or maintenance.

#### **Organize to Ensure Proper Air Flow**

Our Freezer Rack systems not only use all optimal storage space available; they also allow for air flow through the unit. With materials stacked and organized, it makes it easier to locate samples quickly, greatly reducing the amount of time the door remains open.

# **Remove Ice and Frost Routinely**

An accumulation of ice and frost on racks and doors can interrupt air flow and overwork the freezer. Defrost on a routine basis to keep your unit in proper working order.

# **Clear Debris from Compressor Coils**

Vacuum or use a brush to clean the compressor coils on the outside of the unit to remove any buildup of dust and grime on a regular basis to keep your unit running smoothly.

Also, be sure not to stack any items on or in front of your units venting system. This could cause the compressor to work overtime to maintain temperatures.

#### **Clean Door Gaskets**

To ensure that cold air is not escaping, clean any debris or ice build-up from the door gasket to keep an air-tight seal.

## **Freezer Thermometer**

Protect frozen and chilled samples from potential spoilage

Monitor temperatures inside your freezer or refrigerator with the added convenience of sound and flashing LED lights to alert out-of-range temperatures

- Suction cup attaches the probe inside freezer/refrigerator while the thermometer secures to the outside with hook-and loop tape
- Large digital display shows freezer/refrigerator temperature plus ambient temperature
- Switches toggle for freezer/refrigerator monitoring and °F/°C
- Set freezer range: 14.9° to -9.4°F (-9.5° to -23°C);
   Set refrigerator range: 35.6° to 46.4°F (2° to 8°C)
- Includes one AAA battery



Item No.	LxWxH	LxWxH
	in	ст
HS28630	4.1 x .79 x 0.9	1.4 x 1.9 x 2.8

# Microscopes, Slide Storage & Accessories

Products meeting the demanding standards of clinical, research and educational labs



Microscopes
Boxes & Storage Racks
Folders & Mailers
Holders & Staining

# **VanGuard 1220CM Compound Microscope**

The 1220CM is an economical microscope ideal for use in Education/University, Medical/Veterinary, and Industrial applications

- The drop-down X-Y stage control is positioned close to the focus control for quick and easy adjustments
- Accessories and parts are available; contact VEE GEE for more information

Specifications		1220CM
ptical System		Fixed 160 mm
	Туре	Binocular (Siedentopf)
	Tube Inclination	30°
ving Head	Rotation	360°
	Interpupillary	2.2-3 in / 55-75 mm
	Dioptric	-5 to +5
	Magnification	10X
	Туре	Widefield
epieces	Field Diameter	0.7 in / 18 mm
	High-Eyepoint	Yes
	Rubber Eyecups	Yes
	Positions	4
epiece	Pitch	Reverse
	Illumination Technique	Brightfield
	Туре	Achromatic
ectives	Magnification (N.A.)	4X (0.10), 10X (0.25), 40X (0.60), 100X (1.25)
	DIN Standard	Yes
	Туре	Mechanical
	Shape	Rectangular
	Dimensions (L x W)	5.5 x 4.3 in / 140 x 110 mm
ige	X-Y Control	Drop-down, Right Hand, Rack and Pinion
<b>3</b> ·	Cross Travel Range	3 x 1.6 in / 75 x 40 mm
	Slide Holder	Removable spring clip
	Finish	Chemical Resistant
	Туре	Coaxial
	Controls	Coarse / Fine
	Tension Adjustment	Yes
using	Safety Auto-Stop	Yes
	Focusing Range	0.6 in / 15 mm
	Precision	0.002 mm
	Туре	Abbe 1.25 N.A. with Iris Diaphragm
ndenser	Focusing Control	Rack and Pinion
	Type	20W Variable Quartz Halogen
umination	Koehler	Field Diaphragm
	Filter Holder	On Collector
	Base (L x W)	8.3 x 7.1 in / 210 x 180 mm; Anti-Skid Rubber Feet
	Body	Ergonomic; Stain-Resistant Enamel Finish
ecifications	Dimensions (L x W x H)	8.3 x 7.1 x 15 in / 210 x 180 x 380 mm
Johnoudono	Weight	7.0 kg
	Voltage	110 VAC
arranty	5-Year Optical/Mechanical; 1-Year	
pplied With		Spare Lamp, Fuses and Dust Cover

# **VanGuard 1300-Series Compound Microscopes**

Ultra-low focus control, angled viewing head, and interpupillary / dioptric adjustments make this microscope comfortable to use

- The compact frame allows easy arm placement to the controls. Plan objectives provide exceptionally high optical flatness and visual clarity
- Choose from binocular or trinocular viewing heads, and Brightfield or Phase Contrast configurations
- Accessories and parts are available; contact VEE GEE for more information

Specifications		1321BRi	1331BRi	1323PHi	1333PHi					
otical System		VIOS Infinity Correct	ed							
	Туре	Binocular (Siedentopf)	Trinocular (Siedentopf)	Binocular (Siedentopf)	Trinocular (Siedentopf)					
	Tube Inclination	30°								
ewing Head	Rotation	360°								
	Interpupillary	2.2-3 in / 55-75 mm	3 in / 55-75 mm							
	Dioptric	-5 to +5								
amera Port		N/A	20/80 Spilt	N/A	20/80 Spilt					
	Magnification	10X								
	Туре	Widefield								
epieces	Field Diameter	0.8 in / 20 mm								
	High-Eyepoint	Yes								
	Rubber Eyecups	Yes								
epiece	Positions	4								
	Pitch	Reverse								
	Illumination	Brightfield		Brightfield / Darkfiel	d / Phase Contrast					
	Туре	Plan Achromatic		Plan Phase						
jectives	Magnification (N.A.)	4X (0.10), 10X (0.25 (1.25)	), 40X (0.60), 100X	10X (0.25), 20X (0.4 (1.25)	0), 40X (0.60), 100X					
	DIN Standard	Yes								
	Туре	Mechanical								
	Shape	Rectangular								
	Dimensions (L x W)	5.5 x 4.3 in / 140 x								
age	X-Y Control		and, Rack and Pinion							
	Cross Travel Range	3 x 2 in / 75 x 50 m								
	Slide Holder	Removable spring c	ip							
	Finish	Chemical Resistant								
	Туре	Coaxial, Low-Position	n							
	Controls	Coarse / Fine								
ocueina	Tension Adjustment	Yes								
ocusing	Safety Auto-Stop	Yes								
	Focusing Range	0.4 in x 10 mm								
	Precision	0.002 mm								
ndenser	Туре	Abbe 1.25 N.A. with	Iris Diaphragm							
Huensei	Focusing Control	Rack and Pinion			1					
	Туре	20W Variable Quartz	Halogen							
mination	Koehler	Field Diaphragm								
	Filter Holder	On Collector								
	Base (L x W)	8.9 x 6.9 in / 225 x	175 mm; Anti-Skid Rubb	er Feet						
	Body	Ergonomic; Stain-Re	sistant Enamel Finish							
pecifications	Dimensions (L x W x H)	14.6 x 7.5 x 15.7 in 370 x 190 x 400 mr	14.6 x 7.5 x 18.9 in 370 x 190 x 480 mm	14.6 x 7.5 x 15.7 in 370 x 190 x 400 mn	14.6 x 7.5 x 18.9 in 370 x 190 x 480 mm					
	Weight	6.5 kg		6.8 kg						
	Voltage	110 VAC								
arranty	5-Year Optical/Mechanic	al; 1-Year Electrical								
Supplied With	Clear Diffusing Filter, Blu	e Filter, Spare Lamp, F	uses and Dust Cover	Clear Diffusing Filter Green Filter, Telesco Spare Lamp, Fuses a	oing Eyepiece,					

# **VanGuard 1400-Series Compound Microscopes**

Built around a full-size frame, the 1400-Series eliminates the fatigue and strain associated with working long hours over a microscope

- Low-position focusing controls and long drop-down X-Y controls ensure simple, consistent, and comfortable use. Plan objectives provide exceptionally high optical flatness and visual clarity
- Choose from binocular or trinocular viewing heads, and Brightfield or Phase Contrast configurations
- Accessories and parts are available; contact VEE GEE for more information

-	VIOS Infinity Corrected
-	•
Туре	Binocular
	(Siedentopf)
Tube Inclination	30°
	360°
· · · · ·	2.2-3 in / 55-75 mm
Dioptric	-6 to +6
	N/A
Magnification	10X
Туре	Widefield
Field Diameter	0.8 in / 20 mm
High-Eyepoint	Yes
	Yes
	5
	Reverse
	Brightfield / Darkfield / 1423PHi
munimuuvii	Phase Contrast
Туре	Plan Phase
Magnification (N.A.)	10X (0.25), 20X (0.40), 40X (0.60), 100X (1.25)
DIN Standard	Yes
Туре	Mechanical
Shape	Rectangular
Dimensions (L x W)	6.1 x 4.3 in / 155 x 110 mm
X-Y Control	Drop-down, Right Hand, Rack and Pinion
Cross Travel Range	3 x 2 in / 75 x 50 mm
	Removable spring clip
	Chemical Resistant
	Coaxial, Low-Position
	Coarse / Fine
	Yes
•	Yes
	0.8 in / 20 mm
	0.002 mm
	Abbe 1.25 N.A. with Iris Diaphragm
	Rack and Pinion
	20W Variable Quartz Halogen
Koehler	Field Diaphragm
Filter Holder	On Collector
Base	10.2 x 6.7 in / 260 x 170 mm; Anti-Skid Rubber Feet  Fragonomic: Stain-Resistant Enamel Finish
Body	Ergonomic; Stain-Resistant Enamel Finish
Dimensions (L x W x H)	14.9 x 8.6 x 16.7 in 380 x 220 x 425 mm 380 x 220 x 530 mm
Weight	9.3 kg
Voltage	110 VAC
5-Year Optical/Mechan	ical; 1-Year Electrical
	Rotation Interpupillary Dioptric  Magnification Type Field Diameter High-Eyepoint Rubber Eyecups Positions Pitch Illumination Type Magnification (N.A.) DIN Standard Type Shape Dimensions (L x W) X-Y Control Cross Travel Range Slide Holder Finish Type Controls Tension Adjustment Safety Auto-Stop Focusing Range Precision Type Koehler Filter Holder Base Body Dimensions (L x W x H) Weight

# VanGuard 1353SL Stereo Microscope

VanGuard 1353SL Stereo Microscope features precision optics to ensure the highest level of performance for inspection and dissection applications

- Upper and lower LED lighting provides bright and efficient illumination
- Large stage area gives you room to work with bigger objects
- Accessories and parts are available; contact VEE GEE for more information

Specifications		1353SL			
Optical System		Fixed 160 mm			
	Туре	Binocular			
	Tube Inclination	45°			
Viewing Head	Rotation	360°			
	Interpupillary	2.1-3.3 in / 53-83 mm			
	Dioptric	-5 to +5			
	Magnification	10X			
Funciana	Туре	Extra Widefield			
Eyepieces	Field Diameter	0.8 in / 20 mm			
	Rubber Eyecups	Yes			
Stereo Objectives		1X and 3X			
	Plate	Translucent Acrylic			
Stage	Clips	Stainless Steel, Spring-Loaded			
Focusing Movement		Dual Knobs; Rack and Pinion Assembly			
Vertical Travel		125 mm			
	Upper	LED			
Illumination	Lower	LED			
	Base	11.4 x 10.2 in / 290 x 260 mm			
	Body	Ergonomic; Stain-Resistant			
Specifications	Dimensions (L x W x H)	11.4 x 10.2 x 19.7 / 290 x 260 x 500 mm			
.,	Weight	4.5 kg			
	Voltage	110 VAC			
Warranty	5-Year Optical/Mechanical;	5-Year Optical/Mechanical; 1-Year Electrical			
Supplied With	Eyecups (2) and Dust Cover				




# **USB Color Digital Microscope Camera Systems**

Microscope cameras produce high-definition live images and video for documentation of clinical, life sciences, and material science applications

- Cameras feature still and video capture modes, automatic / manual exposure control, manual gamma control, and automatic / manual white balance control
- No separate power cable is required; video, power, and full camera control are all enabled through a single USB cable
- The imaging software allows the user to control all of the advanced settings and has a simple but powerful graphic interface
- Software requires: Windows 7 or later, Intel Core 2 Duo or later, 2 GB of RAM, 2GB hard disk free space, and USB 2.0 / 3.0 port

Item No.	Camera Type	For microscope Models
1400-CDPC-3	USB Digital Camera, 3 Megapixels	1220CM, 1300-Series, 1400-Series
1400-CDPC-5	USB Digital Camera, 5 Megapixels	1220CM, 1300-Series, 1400-Series
1400-CDPC-10	USB Digital Camera, 10 Megapixels	1220CM, 1300-Series, 1400-Series
M1397-K005-1	USB Digital Camera, 3 Megapixels	1353SL
M1397-K007-1	USB Digital Camera, 5 Megapixels	1353SL
M1397-K009-1	USB Digital Camera, 10 Megapixels	1353SL



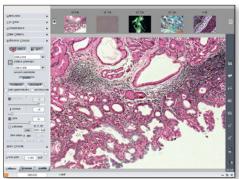

Tcapture: Camera Software Screen



# **100-Place Premium Plus Slide Box**

Modern design for maximum functionality

#### **Optimal Organization**

- Organize by "project" and "date" molded in on lid
- Inventory card index corresponds to slide numbers and is removable for trouble-free writing on flat surface
- Vibrant translucent colors are easy to identify and enable coding options

#### **Functional Design**

- Flexible storage, store vertically or horizontally to efficiently use space
- Molded feet allow secure horizontal stacking of multiple units
- Unique side storage can accommodate a pen, desiccant or other small items needed for working
- Hinged lid smoothly opens 180° sitting flat on surfaces and prevents misplacement

#### Safe & Secure

- EVA foam lines the bottom of the box helping to hold slides securely in place and protect against damage
- Stylish wave surface texture is easy to hold for safe sample transport

# 100-Place Premium Plus Slide Box

Pack of 5

Item No.	Color	Item No.	Color	
120576	Blue	120579	Red	
120577	Orange =	120592	Purple	
120578	Green	120593	Assorted*	

\*Assorted colors (Purple, Blue, Green, Red, Orange)

Lining	Lock	L x W x H in	L x W x H cm
EVA Foam	Sliding Latch	8.2 x 7.7 x 1.3	20.7 x 19.5 x 3.2

# Replacement Inventory Card

Pack of 10

Item No.	L x H in	L x H cm
120618	8.5 x 11	21.6 x 28

# **Premium Microscope Slide Boxes**

Slide boxes are perfect for storage or transportation of standard size slides (76 x 26 mm and 3 x 1 in)

#### Organize & Identify

- Lining and slots are numbered to correspond to the ABS plastic boxes available with either slide inventory sheet
- Organize by "Project" and "Date" with pre-molded labeling area
- Available in multiple colors for project coding

#### Durable & Safe

- cork or foam lining
- Designed to securely store slides in grooved slots
- Molded feet secure stacking for easy transport and storage




# 25-Place Slide Box

Easily access slides with friction-fit cover



Item No.	Color	Lining	Lock	L x W x H in	L x W x H cm
HS15989A	Blue				
HS15989B	Green	Corle	Friction Fit	5.6 x 3.5 x 1.4	14.1 x 8.8 x 3.5
HS15989C	Red	Cork	FIICUOII FIL	5.6 X 3.5 X 1.4	14.1 X 0.0 X 3.3
HS15989E	White				
HS15989M	Blue				
HS15989P	Red	Foam	Friction Fit	5.6 x 3.5 x 1.4	14.1 x 8.8 x 3.5
HS15989R	White				

# **50-Place Slide Box**

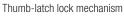



Item No.	Color	Lining	Lock	L x W x H in	L x W x H cm
HS15996A	Blue				
HS15996B	Red	Cork	Thumb-Latch	8.2 x 3.4 x 1.4	20.8 x 8.6 x 3.5
HS15996C	White				

## **100-Place Slide Box**

Closes securely with a rust-resistant nickelplated clasp and pin hinges



Item No.	Color	Lining	Lock	L x W x H in	L x W x H cm
HS15994A	Blue				
HS15994B	Green				
HS15994C	Red		Niekal platad		
HS15994D	Yellow	Cork	Nickel-plated	8.2 x 6.9 x 1.3	20.8 x 17.5 x 3.4
HS15994E	White		Clasp		
HS15994F	Grey				
HS15994G	Black				
HS15994M	Blue		Niekal plated		
HS15994P	Red	Foam	Nickel-plated	8.2 x 6.9 x 1.3	20.8 x 17.5 x 3.4
HS15994R	White		Clasp		

# 100-Place Tall **Slide Box**

Additional height for 3 x 2 in slides Stackable with 100-Place Premium Slide Boxes



Item No.	Color		Lining	Lock	LxWxH in	L x W x H cm
HS15979A	Blue	•	Cork	Nickel-plated Clasp	8.2 x 6.9 x 2.4	20.8 x 17.5 x 6



# Storage Rack for Microscope Slide Boxes

Create a file system for 50-Place and 100-Place slide boxes

#### Storage rack makes organization of samples easy

- Clear acrylic construction allows the labels to be seen on the slide boxes
- Place boxes in rows or columns by color to create an easy coding system

Item No.	Color	Capacity	L x W x H	L x W x H
HS159890B	Clear	15 Bins	28 x 6.7 x 9	71.1 x 17.1 x 22.8

# **Standard Microscope Slide Boxes**

Constructed of durable polypropylene



#### Slide boxes hold either 25 or 100 standard slides (76 x 26 mm and 3 x 1 in)

#### Improved Identification

- Molded slot numbers correspond directly to the slide inventory card
- Organize with "Project" and "Date" pre-molded labeling area

#### Secure & Stackable

- Hinged lids are held securely closed with a latch lock
- Molded feet secure stacking for easy transport and storage

#### 25-Place Slide Box

Item No.	Color	Lining	Lock	LxWxH in (cm)
HS15990A	Blue	None		5.6 x 3.6 x 1.2
HS15990B	Red	None	Latch Lock	(14.1 x 9.2 x 3)
HS15990C	White			(14.1 x 9.2 x 3)
120740	25 Dlac	Poplacoment Invente	5.3 x 3.1	
120/40	20-Pla0	e Replacement Invento	(13.4 x 8)	

#### 100-Place Slide Box

Item No.	Color	Lining	Lock	L x W x H in (cm)
HS15991A HS15991B	Blue Red	- None	Latch Lock	8.3 x 6.7 x 1.5
HS15991C	White	_		(21 x 16.9 x 3.7)
120741	100-Place Replace	3.9 x 6.1 (10 x 15.5)		

# **Low Temperature Microscope Slide Boxes**

Store in ultra-low temperatures with confidence



#### Holds standard slides (76 x 26 mm and 3 x 1 in)

#### **Durable Construction**

- Made of polycarbonate to store in ultra-low temperatures
- Uniform grooved slide slots separate slides preventing slide contact
- Molded feet enable secure stacking for easy transport and storage

#### **Durable Construction**

- Low temperature resistant numbered foam lining corresponds with identification sheet
- Organize with "Project" and "Date" pre-molded labeling area



# Low Temperature Microscope Slide Boxes

Item No.	Color		Capacity	Lock	L x W x H in (cm)
120316	Blue	•	25-Place	Friction Fit	5.6 x 3.5 x 1.4 (14.1 x 8.8 x 3.5)
HS15988H	Purple	•	100-Place	Stainless Steel Clasp	8.2 x 6.9 x 1.3 (20.8 x 17.5 x 3.4)



# **Cardboard Slide Tray**

Cardboard slide trays hold standard slides (76 x 26 mm and 3 x 1 in)

Cardboard slide folders offer a convenient solution for safely working, transporting and short-term storing of slides

#### Organize & Protect

- Index on the front cover corresponds with the number grid inside for identification
- Recessed slide compartments and inside covers for specimen protection
- Slides cannot be accessed when cover is closed reducing contamination

#### Sturdy Construction

- Compartment holder is made from a double layer of densely pressed cardboard
- Cover is single layer cardboard
- Front protection flaps are attached with durable tape

Item No.	Color	Capacity	L x W x H in	L x W x H cm
HS9902	Black			
HS9913	Blue	_		
HS9914	Green	20-Place	8.2 x 13.4 x 0.4	20.8 x 34.1 x 0.9
120196	Red			
120197	Yellow			
HS9930	Black	30-Place	8.2 x 12.9 x 0.4	20.8 x 32.8 x 1

# **Slide Tray Clasps**

Effortlessly locate, retrieve and keep slide folders organized

 Assorted color clasps are ideal for color coding identification and to keep slide folder flaps securely closed  Clasps can be labeled with a permanent marker or a label can be attached

#### Slide Tray Clasps

 Item No.
 Color
 L x W x H in
 L x W x H cm

 120653
 Assorted\* □ ■ ■ ■ 2.2 x 1.7 x 0.6
 5.6 x 4.4 x 1.6

\*Assorted colors (Clear, Grey, Blue, Green, Red)

Pack of 100



# Have a Question?

Our experienced and responsive CUSTOMER SERVICE and TECHNICAL SUPPORT teams are ready and available to help.





# **Slide Mailers**

Safe and secure transport of standard microscope slides  $(76 \times 26 \text{ mm})$  and  $3 \times 1 \text{ in})$ 

#### **Cardboard Slide Mailers**

Item No.	Color	Capacity	UOM	L x W x H in	L x W x H mm
HS9903	Natural	1-Place	Pack of 25	3.7 x 1.7 x 0.2	95 x 42 x 5
HS9904	Natural	2-Place	Pack of 36	4.1 x 3.1 x 0.2	103 x 79 x 5
120557	White	4-Place	Pack of 20	4.7 x 5.7 x 0.3	120 x 140 x 6

- Heavy cardboard mailers for safe transport
- Cover that opens for viewing and thumb grooves for easy slide removal
- Securely holds slides with or without cover glass
- Flat surface area is ideal for labeling or writing necessary information



#### **Plastic Slide Mailers**

Item No.	Color	Capacity	UOM	LxWxH in	LxWxH mm
120556	Natural	1-Place	Pack of 50	3.3 x 1.6 x 0.2	83 x 41 x 5
HS15983B	Natural				
HS15983A	Assorted*	2-Place	Pack of 25	3.7 x 2.7 x 0.2	93 x 68 x 5

\*Assorted colors (Blue, Purple, Red, Orange, Yellow)

- Polypropylene protectors ideal for secure shipping and storing of slides
- Large, easy-to-secure snaps, provide safe storage for valuable slides
- Variety of colors brings some fun to the lab and enables color coding



1-Place

## Plastic 5-Place Slide Mailers

Item No.	Color	Capacity	UOM	L x W x H in	L x W x H mm
HS15982	Natural	5-Place	Pack of 25	3.2 x 0.7 x 1.2	81 x 19 x 31
HS15983G	Assorted*	5-Place	Pack of 25	1.7 x 0.9 x 3.5	43 x 24 x 88
HS15986	Natural				

\*Assorted colors (Blue, Purple, Red, Orange, Yellow)

- Polypropylene containers ideal for secure shipping, storing or staining slides
- Polypropylene containers hold five standard slides (76 x 26 mm and 3 x 1 in)
- Inside wall features slotted channels to keep slides separated and to make removal easier
- Lid closes securely to keep slides safe

Side Opening HS15982



HS15983G

2-Place

### **Slide Holders**

Item No.

HS15981A

HS15981B

HS15981C

Color


Natural

Natural

Natural

# Hold slides during evaluation, cleaning and drying

- Made from chemical-resistant polyethylene making the unit resist stains and easy to clean
- Solid, strong construction with rubber feet for benchtop stability
- Angled slide holder slots allow easy access to individual slides and viewing of slide labels

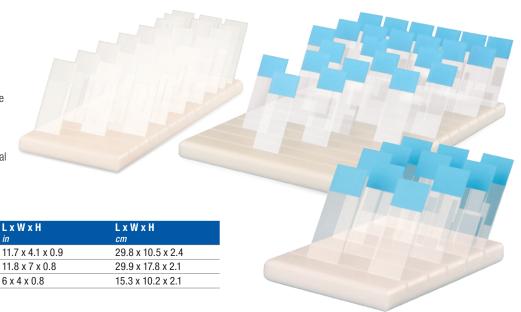



Capacity

40-Place

78-Place

15-Place

### 8-Slide Staining Tray Set

Reduce possible cross-contamination and messy work areas

This rugged staining tray is an economical alternative to metal versions

- Made from a polypropylene and polyethylene blend for stability
- Stain, rinse, and dry slides on a single working tray
- A clean and fresh working surface area every time
- Up to eight slides fit comfortably onto the base and the deep well holds up to 38 mL
- Compact size and recessed handles for easy transport from work area to sink
- Pack includes 4 trays and 1 lid

### 8-Place Slide Staining Tray

for guick disposal of unwanted liquid waste

Item No.	Color	Capacity	L x W x H in	L x W x H cm
HS15951A	Black ■	8-Place	12 x 5 x 1.2	30.5 x 12.7 x 3

# **18-Place Slide Staining Tray Set**

**Ideal for Clinical and Research laboratory staining** 

Tray holds standard microscope slides (76 x 26 mm) securely in place for staining, rinsing and drying applications

### Robust design

- Durably constructed from solvent resistant polypropylene blend
- Non-skid Silicone feet keep tray in place

#### Convenient to use

- Two ways to empty stain liquid
  - Pour spouts on each corner
  - Drain hole with plug
- Equipped with handle cut-outs for easy handling

#### Versatile for multiple testing formats

- Black lid included for light sensitive applications
- Humidity retention wells for specialized staining needs



### 18-Place Slide Staining Tray

Item No.	Color	Capacity	L x W x H in	LxWxH cm
120879	Black	18-Place	16.1 x 8.9 x 1.6 in	41 x 22.5 x 4 cm

# Safety Products & Dispensers

Protect and comply with our range of safety products and PPE dispensing solutions



Wash Bottles
Biohazard Containment
PPE Dispensers & Holders

# **Wash Bottles Color-Coded**

**Color-coded for easy identification** 

The closure and stem are molded in one piece to make them leak-proof

The angled neck and tip provide a uniform stream of liquid and the tip can be cut back if required to increase the flow

- Color-coding for user, application or location identification
- Wide bottle mouth makes filling easier Bottle is made from LDPE, with a PP closure

			Pack of 6
Item No.	Cap/Lid Color	Capacity	
120245	White		
120246	Red <		
120247	Yellow	500 mL	
120248	Green		
120249	Blue		



# Create a safer working surface with our reusable Silicone Lab Mats

- Anti-skid to keep glassware in place
- Protect counter tops from stain, spills and wear
- Available in multiple colors and sizes



# Wash Bottles Safety Labeled and Self Venting

### Clear and informative product safety information

The information provided on the front of each bottle is specific to the product contained and includes the following — the Chemical Abstract Service (CAS) No., National Fire Protection Association (NFPA) codes and symbols, Protective Equipment Symbol (PES), Target Organ Information and Hazard Codes

There is also the option of using other liquids within the bottles and creating the appropriate safety information on the outside of the bottle by using the self-labeling versions

The closure and stem are molded in one piece to make them leak-proof, the angled neck and tip provide a uniform stream of liquid and the tip can be cut back to increase the flow

- Wide bottle mouth makes filling easier
- Sodium Hypochlorite bottle made from natural LDPE for light protection and has a HDPE closure
- All other bottles are made from LDPE, with a PP closure
- Non-metallic vent that prevents the risk of solvent leakage
- Six named chemicals, a self-labeling\*\* and assorted pack of bottles to choose from

Pack of 6

Item No.	Cap/Lid C	olor	Description	Capacity
120250	Red		Acetone	
120251	Clear		Distilled Water	
120252	Yellow		Isopropanol	
120253	White		Ethanol	500 mL
120254	Green		Methanol	JUU IIIL
120255	Blue		Sodium Hypochlorite	
120256	White		Self-labeling**	
120257	Assorted*		Assorted*	

<sup>\*</sup> Assorted Pack includes one of each: acetone, distilled water, isopropanol, ethanol, sodium hypochlorite & a self-labeling bottle

<sup>\*\*</sup>Self-labeling bottles must only be used with chemicals that are safe and compatible with LDPE.

Labeling stickers included



# smartSlicer™ Laboratory Safety Razor Blade

Say no to dangerous sharps

- At just 0.27 mm the US Patent Pending product is as thin as a single edge industrial razor blade, but much safer
- Cleanly slice agarose gels right on the light box without scratching the surface
- Ideal for cutting soft tissue and plant matter
- Adjust the size required by simply cutting with scissors
- Made from food grade PET

Pack of 100

Item No.	L x W x H in	L x W x H
120259	0.8 x 1.5 x 0.01	19 x 38 x 0.3
120260	1 x 2 x 0.01	25 x 50 x 0.3

smartSlicer™ is a trademark of LevGo Inc.

# Safely seal PP and PE bags, sheets or flat tubing with our Heat Sealer

- Heating element measures 31 cm long and is only hot when sealing
- Remove excess material safe with integrated sliding blade



Heat Sealer with Cutter, Page 156



# **Lead Rings**

Stabilize vessels and labware on the bench or in water baths

- Vinyl-coated lead rings will not mark or scratch benchtops and labware
- C-shape format allows some manual size variation for a customized fit

Item No.	Shape	Color	Flask Sizes mL	<b>ID</b> in	ID cm	<b>OD</b> in	OD cm	Weight g
HS8881A	С	- Blue <b>=</b>	125-500	1.9	4.2	3	7.5	215
HS8882A	0	Diue =	125-500	1.9	4.8	3	7.7	210
HS8881B	С	Dluo =	e ■ 250–1000	2	5.1	3.5	9	465
HS8882B	0	Blue		2	5.1	3.5	9	490
HS8881C	С	Dluo	Blue ■ 500–2000	2.1	5.4	3.9	9.8	595
HS8882C	0	blue		2.2	5.7	3.9	9.8	610
HS8881D	С	Dluo -	1000 4000	2.6	6.6	4.4	11.3	890
HS8882D	0	Blue	1000-4000	2.8	7	4.5	11.5	935

\*Weights may vary by +/- 3%



Safety Bottle Carriers, Page 85

# Safety Bottle Carrier

- Handle swings down to either side for easier loading and unloading
- Carriers come in four colors with designated labeling areas for easy identification

# **Biohazard Disposal Bags**

For personal or general lab use, sizes to suit all needs

- Polypropylene
- A range of biohazard collection bags printed with biohazard warning symbol and precautions in English, Spanish and French
- Includes an indicator patch that darkens when autoclaved



### **Biohazard Disposal Bags**

Item No.	Color	Thickness mil	Capacity	Packaging	L x W in	L x W cm
HS10320			2 L	500	8 x 12	20.3 x 30.5
HS10321	- - Red ■	1.0	9 L	200	14 x 19	35.6 x 48.3
HS10322	- Reu -	1.6	19 L	200	19 x 23	48.3 x 58.4
HS10323	_		40 L	200	25 x 33	63.5 x 83.8

### **HDPE Coated Steel Wire Bag Holder**

Item No.	LxWxH in	L x W x H cm
HS10044H	3.6 x 5 x 8.4	9.1 x 12.8 x 21.3

\*Holds Biohazard Disposal Bag HS10320





# **Autoclavable Disposal Bags**

Clear bags imprinted with blue biohazard labeling designed for safe disposal

- Constructed of 2 mil thick polypropylene, Front write-on panel allowing for repeated autoclaving
- Strong side seams for securely holding items
- 4 inner packs of 50 bags



Case of 200

Item No.	Color	Capacity	L x W in	L x W cm
HS1002A		14 L	12.2 x 26	31 x 66
HS1002B	Clear	26 L	16.3 x 23.6	41.5 x 60
HS1002C		 69 L	24 x 31.9	62 x 81

# **Sharps Containers**

Easy to use and transport

- Clearly labeled with the Biohazard warning on one side
- Provides safe disposal of syringes, needles and other small sharps
- Stackable to reduce storage space before use
- 1 L Capacity (quart) with a clear lid to view fill level
- Slide access lid provides temporary and final closures for safety and security
- Built-in needle removal port, prevents direct contact with sharps
- Compatible with the Droplet blood collection tray, Totemax® and blood collection tray, Pages 63 and 64

Item No.	Color	L x W x H in	L x W x H cm
120177	Red <	4.5 x 4.5 x 7.5	11.5 x 11.5 x 19.2
120178	Yellow	4.3 X 4.3 X 7.3	



# **PPE Dispensers & Holders**

### **Soft Cover Dispenser**

Store and dispense with ease

Translucent acrylic dispensers provide easy access to safety accessories

- Hinged lid protects contents from dust and debris
- Set on counter or wall mount (screws included)
- Great for holding other personal items



Item No.	Color	L x W x D in	L x W x D cm
HS1040C	Clear	8 x 8 x 16	20.3 x 20.3 x 40.6



# **Earplug Dispenser**

Translucent acrylic dispensers provide easy access to safety accessories

- Holds approximately 200 pairs of foam earplugs
- Set on counter or wall mount (screws included)



Item No.	Color	L x W x D	L x W x D cm
HS1040B	Neon Green	8 x 8 x 16	20.3 x 20.3 x 40.6

### 3-in-1 Dispenser

Three separate lift-up lids make loading and dispensing supplies easy

Three-compartment dispenser can be used to organize safety glasses, earplugs, soft covers (including bouffant caps, beard covers, shoe covers) or any small items that need to be kept within reach

The unit includes screws for wall mounting or can be set on a counter

- Middle compartment has reach-in hole for dispensing
- Durable acrylic is see-through so product levels are always visible

Item No.	Color	L x W x D	L x W x D cm
HS1042	Grev	16 x 16 x 8	40.6 x 40.6 x 20.3



# **Safety Glasses Dispensers**

**Ergonomically angled for easy dispensing** 

Translucent acrylic dispensers provide easy access to safety accessories

- Hinged lid protects glasses from dust and debris
- Set on counter or wall mount (screws included)
- Large dispenser holds approx. 20 glasses, small approx. 10 glasses
- Great for holding other personal items

Item No.	Color	L x W x D in	L x W x D cm
HS1041	Neon Blue	8 x 6 x 6	20.3 x 15.2 x 15.2
HS1040A	Neon Blue	8 x 8 x 16	20.3 x 20.3 x 40.6



# **PPE Dispensers & Holders**

# **Workstation Storage Bin**

The essential combination accessory holder

- Organize your most commonly used laboratory disposables in one central location
- Acrylic unit has room to store gloves, Kimwipes® and Parafilm® M
- Use bins to store tubes, tips, etc.

Item No.	Color	L x W x H in	L x W x H cm
HS234523	Clear	10.4 x 5.4 x 17.3	26.5 x 13.7 x 44

Parafilm® is a registered trademark of Bemis Inc. Kimwipes® is a registered trademark of Kimberly-Clark





# **Lab Wipe Holders**

Holds Kimwipes® and other major brands

Clear acrylic holders fit boxes of small or large lab wipes

- The smaller version attaches with double-sided adhesive tape and the larger version with screws
- Mount on a shelf or wall
- Mounting materials included

Item No.	Description	Color	L x W x H in (cm)
HS234521	Small Holder	Clear $\square$	4.6 x 4.8 x 3.3 (11.8 x 12.2 x 8.4)
HS234522	Large Holder	Clear $\square$	12.6 x 5 x 3.8 (32 x 12.7 x 9.7)

Kimwipes® is a registered trademark of Kimberly-Clark

# **Benchtop Dispensing Bins**

Work great with gloved hands, bright colors

Ergonomically designed with large, angled openings and oversized lip

- Lid remains open when tilted back
- Nonskid rubber feet and pre-drilled holes for wall mounting
- Green bin includes a removable divider
- Mounting screws included

Item	No.	Color	Compartments	L x W x H in	L x W x H cm
HS2	3411	Neon Green	2	7.1 x 6.1 x 6.5	18 x 15.5 x 16.5
HS2	3412	Neon Blue	1	6.1 x 6.1 x 6.7	15.5 x 15.5 x 17



# **Clear Acrylic Glove Box Holders**

Keep gloves organized & ready for immediate use

- Free up bench-top space by mounting holders on wall for easy dispensing of gloves
- Universal fit for most major brands of glove boxes
- Mounting screws included
- Includes a polyester foam insert to securely hold a range of glove box sizes—snug fit ensures easy removal of gloves from the box

Item No.	Capacity	L x W x H in	L x W x H cm
HS23471A	1	10.3 x 5.9 x 4.5	26.3 x 14.6 x 11.4
HS23471B	2	11.4 x 10.3 x 4.5	29.2 x 26 x 11.4
HS23471C	3	11.6 x 15.2 x 4.3	29.4 x 38.8 x 11.6





### **HDPE Coated Wire Glove Box Holders**

Keep gloves organized & ready for immediate use

- HDPE-coated steel wire racks are an economical solution
- Free up bench-top space by mounting holders on wall for easy dispensing of gloves
- Universal fit for most major brands of glove boxes
- Mounting screws included

Item No.	Capacity	L x W x H in	L x W x H cm
HS23457	1	5.5 x 8.2 x 4.3	14 x 20.9 x 10.8
HS23458	3	9.8 x 17.9 x 4.3	25 x 45.4 x 10.8

# **Modular Glove Box Holders**

Available infused with Microban® for anti-microbial protection

The anti-microbial version (blue) is ideal for use in clean rooms, hospitals, and microbiological work areas. The additive prevents growth of molds and bacteria

- High-impact, durable resin is resistant to harsh cleaning solutions
- Modular design allows single boxes or horizontal stacking for multiple box holders
- Universal fit for most major brands of glove boxes
- Mounting screws included

US Patent No. 6669157 Pack of 3

Item No.	Description	Color	L x W x H in	L x W x H cm
120258	Standard	White	— 10 x 6.1 x 3.8	25.4 x 15.6 x 9.7
HS23491A	Anti-microbial	Blue	— 10 X 0.1 X 3.0	20.4 x 10.0 X 9.7





### **Flexible Glove Box Holder**

Glove box holder ideal for gloves, face masks, shoe covers and more

The adjustable holding points provide true flexibility to hold up to four different boxed protective products and other lab essentials (size dependent)

- Ideal for holding different boxed sizes of gloves, mop caps, ear plugs, face masks, shoe covers
- Adjustable universal fit
- Easy to clean, with a low profile
- Wall- or door-mountable
- Mounting hardware included

#### Adaptable

Store different size boxes all in one holder

Item No.	Color	L x W x D in	L x W x D cm
120107	White	21.9 x 5.2 x 2.7	55.8 x 13.2 x 6.8

# **Delivery**

We stock our warehouse to ensure product availability and on-time delivery.

By offering low minimum order quantities, we help dealers efficiently maintain lower stocking levels.



# **Sealing & Labeling**

Sample protection and identification for a wide range of tubes and vials



Heat Sealer
Parafilm & Dispensers
Films, Tapes & Foils
Labels

### **Heat Sealer with Cutter**

Quickly and efficiently seal polyethylene or polypropylene bags, sheets or flat tubing from 0.01 to 0.8 mm thick

### Sliding blade removes excess material 19 mm from the seal

- Can seal a bag up to 12 in (30.5 cm) long
- Heating element measures 31 cm long and is only hot when sealing
- Includes one replacement part kit
- Additional replacement part kits can be ordered separately
- For use with PE and PP material items only. Not compatible with combination plastic/paper pouches



#### Replacement Part Kit:

- 2 extra blades,
- 2 extra heating elements,
- 2 extra PTFE cloths





Item No.	Power Supply		L x W X H in	L x W X H
HS28270	220-230 VAC, UK			
HS28271	220-230 VAC, Cont. Europe/Israel	<u>••</u>	17.6 x 3.4 x 6.5	44.6 x 8.6 x 16.6
HS28272	110-120 VAC, US/Taiwan/Japan	11		
HS28280	Replacement Part Kit			

# Parafilm® M Sealing Film

Protect the integrity of your laboratory and research applications

# Parafilm® M is a semi-transparent, flexible, thermoplastic film with two-way stretch

- Protect the integrity of research applications with a tight moisture-barrier that holds a strong seal and doesn't release from glassware, despite off-gassing
- Reduce the risk of tampered data with a covering that is odorless, colorless and semi-transparent
- Prevents moisture and volume loss due to material composition and tight seal
- Flexible and self-sealing material allows it to cling around irregular shapes and surfaces
- Stretches more than 200% of original length
- Resistant up to 48 hours against many polar substances such as saline solutions, inorganic acids and alkaline solutions
- Two color options available in 2" rolls (natural and purple) for color coding purposes



Item No.	Color	Film Width		Roll Length	
		in	ст	ft	m
HS234526A	Natural	2	5	250	76.2
120762	Purple	2	5	250	76.2
HS234526B	Natural	4	10	125	38.1
HS234526C	Natural	4	10	250	76.2
120265	Natural	20	50.8	50	15.2

Parafilm® M is a registered trademark of Amcor Limited

# Box Top Dispenser for Parafilm® M Sealing Film

### **Easy application**

Simply place the dispenser on top of the Parafilm® M box and pull the sealing film through, there is no need to remove rolls from the packaging

- Accommodates two 50 mm or one 100 mm Parafilm® M rolls
- Built in safety blade slides and cuts with no danger to fingertips and eliminates the need for scissors
- Supplied with two safety blades, additional replacement blades are available for purchase
- Made with ABS plastic making the unit easy to clean

Item No.		L x W x H in	L x W x H cm
120667	Blue		
120668	Purple	4.9 x 6.8 x 3.9	12.3 x 17.2 x 9.9
120758	Natural $\square$		
Item No.			
120722	Replacemen	t Blades for Box Top Dis	penser, <i>Pack of 3</i>

Parafilm® M is a registered trademark of Amcor Limited





# **ABS Dispenser for Parafilm® M Sealing Film**

Store, dispense, and cut Parafilm® M, tape, labels, and Tough-Tags®

- ABS plastic dispenser designed with integrated safety razor that slides smoothly creating a clean cut
- Supplied with two safety blades, additional replacement blades are available for purchase
- Dispenser holds one or two 50 mm (2-inch) rolls or one 100 mm (4-inch) roll

Item No.		L x W x H in	L x W x H cm
HS234525B	Blue	47×61×67	12 x 15.6 x 17.1
HS234525C	Green	— 4.7 x 6.1 x 6.7	12 X 13.0 X 17.1

<b>ПЭZ343Z3</b> U	areen =	
Item No.		
HS234525R	Replacement Blades for ABS Dispenser, Pack of 5	

Parafilm® M is a registered trademark of Amcor Limited Tough-Tags® is a registered trademark of Diversified Biotech

# **Acrylic Dispenser for Parafilm® M Sealing Film**

Store, dispense, and cut Parafilm® M, tape, labels, and Tough-Tags®

- Strong, crystal-clear acrylic dispenser makes your film supply visible and features a serrated cutting blade
- Rubber feet secures in place while cutting
- Dispenser holds one or two 50 mm (2-inch) rolls or one 100 mm (4-inch) roll

Item No.	L x W x H in	L x W x H cm	
HS234524	4.7 x 6.6 x 6.9	11.9 x 16.8 x 17.5	

Parafilm® M is a registered trademark of Amcor Limited Tough-Tags® is a registered trademark of Diversified Biotech



# **PetriSEAL™** and ContainerSEAL™

### Sealing tape for petri dishes, microplates, lids and caps

- Tape stretches, conforms and seals to irregular perimeters
- Waterproof and chemically resistant to most solvent and caustic agents
- Seal remains flexible even with intermittent exposure to 100°C

#### PetriSEAL™


Item No.	Color	Film Width in	ст	Roll Length ft	т
<b>HEAPSEAL-108B</b>	Blue				
HEAPSEAL-108C	Clear				
HEAPSEAL-108R	Red	0.5	1.3	108	33
HEAPSEAL-108Y	Yellow				
HEAPSEAL-108W	White				





### **ContainerSEAL™**

Item No.	Color/Size	Film Width in	ст	Roll Length ft	т
HEACSEAL-75R	Red				
<b>HEACSEAL-75W</b>	White	8.0	1.9	108	33
HEACSEAL-75Y	Yellow				
HEACSEAL-XLR	Red XL				
HEACSEAL-XLW	White XL	1	2.5	108	33
HEACSEAL-XLY	Yellow XL				

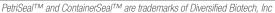





### **Foil Rolls**

### Designed for easy use and minimal waste

- Two sizes available for covering small surface such as beakers, flasks, tubes and Petri dishes
- Withstands autoclaving, incubating and freezing temperatures

Item No.	Thickness mm	Roll Width in	ст	Roll Length ft	т
HS23534A	0.015	4	10.2	500	152.4
HS23536A	0.015	6	15.3	500	152.4

# **Cryo-Babies® and Cryo-Tags® Low Temperature Labels**

For secure uncompromised labeling

- Chemically inert
- Adheres to most plastics (100% virgin PP tubes recommended), glass and metals without cracking, peeling or degrading
- Flexible polyolefin material and acrylic adhesive expand and contract with temperature fluctuations







### Cryo-Babies® on a Roll

oryo-bables on a no			
Ideal Use	0.5 mL Tubes	1.5 mL Tubes	General Purpose
	LABEL	LABEL	LABEL
Label Dimensions	0.94 x 0.50 in (24 x 13 mm)	1.28 x 0.50 in (33 x 13 mm)	1.50 x 0.75 in (38 x 19 mm)
UOM	1,000/roll	1,000/roll	1,000/roll
Item No.			
White	HEACRY-0500	HEACRY-1000	HEACRY-2000
Blue	-	HEACRY-1000-B	HEACRY-2000-B
Green	-	HEACRY-1000-G	HEACRY-2000-G
Orange	-	HEACRY-1000-0	HEACRY-2000-0
Red	-	HEACRY-1000-R	HEACRY-2000-R
Yellow	-	HEACRY-1000-Y	HEACRY-2000-Y

Cryo-Babies® and Cryo-Tags® are registered trademarks of Diversified Biotech, Inc

# For write-on labels, see our selection of Sharpie® Pens




**Cryo-Babies® and Cryo-Tags® Low Temperature Labels** 

For secure uncompromised labeling

- Chemically inert
- Adheres to most plastics (100% virgin PP tubes recommended), glass and metals without cracking, peeling or degrading
- Flexible polyolefin material and acrylic adhesive expand and contract with temperature fluctuations





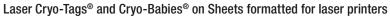



Ideal Use	0.5 mL Tubes	1.5 mL Tubes	1.5 mL Tubes	Cryogenic Vials	Bottles & Beakers	Microplates
Label Dimensions	0.94 x 0.50 in (24 x 13 mm)	1.28 x 0.50 in (33 x 13 mm)	1.28 x 0.50 in (33 x 13 mm)	1.69 x 0.75 in (43 x 19 mm)	3.0 x 2.0 in (76 x 51 mm)	1.50 x 0.25 in (38 x 6 mm)
Labels/Sheet	119	51	85	52	12	156
Labels/pk	2,380	1,020	1,700	1,040	240	3,120
UOM	20/pk	20/pk	20/pk	20/pk	20/pk	20/pk
Item No.						
White	HEALCRY-2380	HEAMINI-1700	HEALCRY-1700	HEALCRY-1100	HEALCRY-3020	HEASIDE-1000
Blue	HEALCRY-2380-B	-	HEALCRY-1700-B	HEALCRY-1100-B	-	HEASIDE-1100
Green	HEALCRY-2380-G	-	HEALCRY-1700-G	HEALCRY-1100-G	-	HEASIDE-1200
Orange	HEALCRY-2380-0	-	HEALCRY-1700-0	HEALCRY-1100-0	-	HEASIDE-1300
Red	HEALCRY-1380-R	-	HEALCRY-1700-R	HEALCRY-1100-R	-	HEASIDE-1400
Yellow	HEALCRY-1380-Y	-	HEALCRY-1700-Y	HEALCRY-1100-Y	HEALCRY-3020-Y	HEASIDE-1500
Rainbow*	HEARNBW-2300	-	HEARNBW-2200	HEARNBW-2000	=	HEARNBW-2400

Laser Cryo-Tags® and Cryo-Babies® on Sheets formatted for laser printers - continued

Ideal Use	Microscope Slides	Racks & Boxes	Racks & Boxes	Racks & Boxes	General Purpose	General Purpose
Label Dimensions	1.0 x 1.0 in (25 x 25 mm)	2.50 x 0.50 in (64 x 13 mm)	2.625 x 1.0 in (67 x 25 mm)	1.125 x 1.75 in (29 x 45 mm)	8.5 x 11.0 in (216 x 279 mm)	1.50 x 0.75 in (38 x 19 mm)
Labels/Sheet	80	60	30	42	1	60
Labels/pk	1,600	1,200	600	840	20	1,200
UOM	20/pk	20/pk	20/pk	20/pk	20/pk	20/pk
Item No.						
White	HEALCRY-1400	HEALCRY-2500	HEALCRY-1258	HEALCRY-1300	HEALCRY-1000	HEALCRY-1200
Blue	-	-	HEALCRY-1258-B	-	-	HEALCRY-1200-B
Green	-	=	HEALCRY-1258-G	-	-	HEALCRY-1200-G
Orange	-	-	HEALCRY-1258-0	-	-	HEALCRY-1200-0
Red	-	-	HEALCRY-1258-R	-	-	HEALCRY-1200-R
Yellow	-	-	HEALCRY-1258-Y	HEALCRY-1300-Y	-	HEALCRY-1200-Y
Rainbow*	-	-	HEARNBW-2500	-	-	HEARNBW-2100

\*Rainbow 4 sheets of each Blue, Green, Orange, Red, Yellow Cryo-Babies® and Cryo-Tags® are registered trademarks of Diversified Biotech. Inc

# Tough-Spots® Tube Labels Pre-cut peel off labels for tube tops

- Adheres to most plastics, glass, and metals without cracking, peeling or degrading
- Flexible polyolefin material and acrylic adhesive



### Tough-Spots® on a Roll

Ideal Use	0.2 mL Tubes	0.5 mL Tubes	1.5/2.0 mL Tubes	1.5/2.0 mL Tubes	15-50 mL Tubes	50 mL Tubes
Label Diameter	0.1875 in/4.8 mm	0.375 in/9.5 mm	0.4375 in/11 mm	0.50 in/13 mm	0.75 in/19 mm	1.0 in/25 mm
UOM	1,500/roll	1,000/roll	1,000/roll	1,000/roll	500/roll	500/roll
Item No.						
White	HEATS-TNY	HEAT-SPOTS	HEATSPOTS-44	HEAT-SPOTS-50	HEATSPOTS-75	HEAT-SPOTS-100
Blue	-	HEAT-SPOTS-B	-	HEAT-SPOTS-50-B	-	-
Coffee	-	HEAT-SPOTS-C	-	HEAT-SPOTS-50-C	-	-
Gold	-	HEAT-SPOTS-D	-	HEAT-SPOTS-50-D	-	-
Green	-	HEAT-SPOTS-G	-	HEAT-SPOTS-50-G	-	-
Lavender	-	HEAT-SPOTS-L	-	HEAT-SPOTS-50-L	-	-
Maroon	-	HEAT-SPOTS-M	-	HEAT-SPOTS-50-M	-	-
Orange	-	HEAT-SPOTS-0	-	HEAT-SPOTS-50-0	-	-
Pink	-	HEAT-SPOTS-P	-	HEAT-SPOTS-50-P	-	-
Red	-	HEAT-SPOTS-R	-	HEAT-SPOTS-50-R	-	-
Silver	-	HEAT-SPOTS-S	-	HEAT-SPOTS-50-S	-	-
Yellow	-	HEAT-SPOTS-Y	-	HEAT-SPOTS-50-Y	-	-





- Label sheets are top coated and heat-resistant to prevent jams in laser printers
- Printer formatting instruction included



### Laser Tough-Spots® on Sheets formatted for laser printers

Ideal Use	0.2 mL Tubes	0.5 mL Tubes	1.5/2.0 mL Tubes	1.5/2.0 mL Tubes	15-50 mL Tubes
	00000000000000000000000000000000000000				
Label Diameter	0.375 in/9.5 mm	0.4375 in/11 mm	0.50 in/13 mm	0.75 in/19 mm	1.0 in/25 mm
Labels/sheet	192	234	192	108	63
Labels/pk	3,840	4,680	3,840	2,160	1,260
UOM	20/pk	20/pk	20/pk	20/pk	20/pk
Item No.					
White	HEASPOT-1000	HEASPOT-5000	HEASPOT-2000	HEASPOT-3000	HEASPOT-4000
Blue	HEASPOT-1100	-	HEASPOT-2100	-	-
Green	HEASPOT-1200	-	HEASPOT-2200	-	-
Orange	HEASPOT-1300	-	HEASPOT-2300	-	-
Red	HEASPOT-1400	-	HEASPOT-2400	-	-
Yellow	HEASPOT-1500	-	HEASPOT-2500	-	-
Rainbow*	HEARNBW-1000	-	HEARNBW-1100	-	-

# **Tough-Tags® Tube Labels** Pre-cut labels for identification

- Adheres to most plastics, glass, and metals without cracking, peeling or degrading
- Flexible polyolefin material and acrylic adhesive





### Tough-Tags® on a Roll

Ideal Use	0.2 mL Tubes	0.5 mL Tubes	1.5/2.0 mL Tubes	General Purpose	Microplates
Label Diameter	0.81 x 0.28 in/21 x 7 mm	0.94 x 0.50 in/24 x 13 mm	1.28 x 0.50 in/33 x 13 mm	1.50 x 0.75 in/38 x 19 mm	1.50 x 0.28 in/38 x 6 mm
UOM	1,500/roll	1,000/roll	1,000/roll	1,000/roll	1,000/roll
Item No.					
White	HEATT-TNY	HEATTSW-1000	HEATTLW-1000	HEATTGP-1000	HEATT-SWALL
Blue	-	HEATTSB-1000	HEATTLB-1000	-	-
Green	-	HEATTSG-1000	HEATTLG-1000	-	-
Orange	-	HEATTSO-1000	HEATTLO-1000	-	-
Red	-	HEATTSR-1000	HEATTLR-1000	-	=
Yellow	-	HEATTSY-1000	HEATTLY-1000	-	-
Clear	-	HEATTSC-1000	HEATTLC-1000	-	-

- Label sheets are top coated and heat-resistant to prevent jams in laser printers
- Printer formatting instruction included

Tough-Tags® is a registered trademarks of Diversified Biotech, Inc.

### Laser Tough-Tags® on Sheets formatted for laser printers

Ideal Use	0.5 mL Tubes	1.5/2.0 mL Tubes	General Purpose	Microscope Slides
Label Dimensions	0.94 x 0.50 in/24 x 13 mm	1.28 x 0.50 in/33 x 13 mm	1.50 x 0.75 in/38 x 19 mm	0.94 x 0.94 in/22 x 22 mm
Labels/Sheet	119	85	60	96
Labels/pk	2,975	2,125	1,500	2,400
UOM	25/pk	25/pk	25/pk	25/pk
Item No.				
White	HEATTSW-2240	HEATTLW-2016	HEATTLW-1050	HEAMISL-1000
Blue	HEATTSW-2240-B	HEATTLW-2016-B	HEATTLW-1050-B	HEAMISL-1000-B
Green	HEATTSW-2240-G	HEATTLW-2016-G	HEATTLW-1050-G	HEAMISL-1000-G
Orange	HEATTSW-2240-0	HEATTLW-2016-0	HEATTLW-1050-0	HEAMISL-1000-0
Red	HEATTSW-2240-R	HEATTLW-2016-R	HEATTLW-1050-R	HEAMISL-1000-R
Yellow	HEATTSW-2240-Y	HEATTLW-2016-Y	HEATTLW-1050-Y	HEAMISL-1000-Y
Rainbow*	HEARNBW-3100	HEARNBW-3000	HEARNBW-3200	HEARNBW-3300

\*Rainbow 4 sheets of each Blue, Green, Orange, Red, Yellow

### **Label Station**

Organize your labels for easy accessibility

Portable label station holds up to six standard size boxes of Tough-Tags®, Cryo-Babies® and Tough-Spots®

Just drop the box in the label station and feed the paper through the slot

- Acrylic station includes a writing surface, a cutting edge for easy dispensing and non-skid rubber feet
- Counter-sunk holes provided for wall mounting

Item No.	Color	L x W x H in	L x W x H cm	
HS23460	Blue	11 x 4.7 x 3.5	28 x 12 x 8.9	

Tough-Tags®, Cryo-Babies® and Tough-Spots® are registered trademarks of Diversified Biotech, Inc.



# **Charts & References**



**Plastic Material References** 

# **Plastic Material References**

# **Chemical Resistance of Plastics**

- **G** Good Resistance: continuous exposure to the substance does not cause damage within 30 days
- F Fair Resistance: continuous exposure to the substance causes minor damage within 7-30 days
  P Poor Resistance: not suitable for continuous exposure to the substance. Immediate damage may occur

	LD	PE	HC	PE	F	P		MP PX)	P	VC	P	C	ABS	/PS	SA	AN	PM	MA	PT	FE	P	DM
TEMPERATURE °C	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50
Acetaldehyde	G	Р	G	F	F	Р	F	Р	Р	Р	F	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Acetic Acid (Glacial)	F	Р	G	G	G	F	F	F	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	Р	Р
Acetic Anhydride	Р	Р	F	F	F	F	G	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	Р	Р
Acetone	F	Р	Р	Р	F	Р	G	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Ammonium Chloride (10%)	G	G	G	G	G	G	G	G	G	F	F	F	G	G	G	G	F	F	G	G		
Ammonium Hydroxide (30%)	G	G	G	G	G	G	G	G	G	Р	Р	F	F	F	F	F	F	F	G	G	G	G
Amyl Acetate	F	Р	G	F	F	Р	G	F	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G	G	G
Aniline (Phenylamine)	G	F	G	G	G	G	G	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	F	F
Aqua Regia	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	F	Р	Р	Р	G	G	G	G
Benzaldehyde	G	F	F	Р	Р	F	G	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Benzene	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	F
Benzoic Acid	G	G	G	G	G	F	G	G	G	F	G	F	F	F	G	G	G	G	G	G		
Boric Acid (10%)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G
Butyl Acetate	F	F	F	F	F	F	F	F	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	F
Butyric Acid (Butanoic acid)	Р	Р	F	Р	Р	Р	Р	Р	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G		
Calcium Hydroxide (Saturated)	G	G	G	G	G	G	G	G	G	G	Р	Р	F	F	G	G	F	Р	G	G	G	G
Carbon Disulphide	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Carbon Tetrachloride	F	Р	F	Р	F	Р	Р	P	G	Р	Р	Р	P	Р	F	Р	P	P	G	G	F	F
Chloroform	F	Р	F	Р	Р	Р	Р	Р	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	Р	Р
Citric Acid (1M)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G		G	G		
Cresol	Р	Р	F	Р	F	F	Р	Р	Р	Р	Р	Р	F	Р	Р	Р	Р	Р	G	G		
Cyclohexane	F	Р	F	Р	F	Р	Р	Р	Р	Р	G	F	Р	Р	F	Р	F		G	G	G	G
Dibutyl Phthalate	F	Р	F	Р	F	Р	F	F	Р	Р	F	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Dichlorobenzene	F	P	F	P	F	P	Р	P	P	P	Р	Р	P	Р	Р	Р	P	Р	G	G	Р	Р
Diethyl Ether	Р	Р	F	Р	F	Р	Р	Р	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Diethylene Glycol	G	G	G	G	G	G	G	G	F	Р	F	F	G	G	G	G	P	Р	G	G	G	G
Dimethyl Formamide (DMF)	G	G	G	G	G	G	G	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Dimethyl Sulfoxide (DMSO)	G	G	G	G	G	G	G	G	P	Р	Р	Р	-	F	Р	Р	P	Р	G	G	-	-
Dioxane	F	F	F	F	P	Р	F	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	F	F
Ethyl Acetate (Ethyl Ester)	G	G	G	G	F	Р	F	Р	Р	Р	Р	Р	F	Р	Р	Р	Р	Р	G	G		
Ethyl Alcohol (Absolute Ethanol)	G	F	G	G	G	G	G	F	G	G	G	F	G	F	Р	Р	Р	Р	G	G	G	G
Ethyl Chloride (Chloroethane)	F	Р	Р	Р	F	Р	F	Р	Р	Р	Р	Р	F	Р	Р	Р	Р	Р	G	G		
Ethylene Chloride	F	Р	Р	Р	F	Р	F	Р	Р	Р	Р	Р	F	Р	Р	Р	Р	Р	G	G	Р	Р
Ethylene Oxide (Pure)	F	F	F	F	F	Р	F	Р	Р	Р	F	Р	Р	Р	Р	Р	G	G	G	G	G	G
Ethylene Oxide (Gas)	F	F	G	G	G	G	F	F	F	F	G	G	G	F	F	F	G	G	G	G		
Formaldehyde (Formalin) 40%	G	F	G	F	G	G	G	G	G	G	G	G	F	F	F	Р	G	F	G	G	G	G
Formic Acid (50%)	F	F	G	G	G	F	G	G	G	Р	F	F	F	F	G	F			G	G		
Formic Acid (100%)	F	F	G	G	G	F	G	G	G	Р	F	Р	F	F	Р	Р			G	G	Р	Р
Glycerine (Glycerol)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	F	F
Hexane	Р	Р	F	F	F	F	F	Р	G	Р	F	Р	Р	Р	G	G	G		G	G	G	G
Hydrobromic Acid (69%)	G	G	G	F	G	F	G	G	G	G	F	F	F	F	F	F	Р	Р	G	G	Р	Р
Hydrochloric Acid (5%)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	F		G	G	Р	Р
Hydrochloric Acid (35%)	G	G	G	G	G	F	G	F	G	G	F	Р	G	G	G	F	Р	Р	G	G	Р	Р
Hydrofuoric Acid (48%)	G	G	G	G	G	F	G	F	F	F	F	Р	Р	Р	G	F	Р	Р	G	G	Р	Р
Hydrogen Peroxide (30%)	G	G	G	G	G	F	G	F	G	G	G	G	G	F	G	G	F	Р	G	G	G	Р
Lactic Acid (85%)	G	G	G	G	G	F	G	F	G	F	G	F	G	G	G	G	G	G	G	G	G	Р
Methyl Acetate	G	Р	G	G	G	F	F	Р	Р	Р	Р	Р	Р	Р	Р	Р			G	G		
Methyl Alcohol (Methanol)	G	F	G	G	G	G	G	F	G	G	F	F	F	F	Р	Р	Р	Р	G	G	G	G
Methyl Ethyl Ketone (Butanone)	Р	Р	Р	Р	G	G	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G		
Methylene Chloride (Dichloro Methane)	Р	Р	G	Р	G	Р	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	Р	Р
Mineral Oil	F	G	G	G	G	F	G	F	G	G	G	G	G	G	G	G	G	G	G	G	G	G
Nitric Acid (10%)	G	G	F	F	F	F	F	F	F	G	G	F	F	G	F	F	G		G	G	Р	Р
Nitric Acid (70%)	G	P	G	P	P	P	G	P	P	Р	P	P	P	Р	P	P	P	Р	G	G	P	P
Nitrobenzene	P	P	Р	P	P	P	G	P	P	P	P	P	P	P	P	P	P	P	G	G	F	P
Oxalic Acid (10%)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G
Perchloric Acid (70%)	F	P	F	P	F	P	F	P	F	P	P	P	F	F	F	F	G	G	F	F	P	P
Phenol (100%)	P	Р	P	Р	Р	Р	G	P	P	Р	Р	Р	Р	P	Р	P	Р	P	G	G	Р	Р
Phosphoric Acid (85%)	G	P	G	G	G	F	G	F	G	F	G	F	G	F	G	G	F	P	G	G	G	P
า แบอคแบบเด พดเน (00%)	u	1	u	u	u	"	u	[	u	[	u	Г	u	Г	u	u	Г	۲	u	u	u	"

	LD	PE	HD	PE	Р	Р		MP PX)	P	/C	Р	C	ABS	S/PS	SA	AN	PM	MA	PT	FE	PC	OM
TEMPERATURE °C	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50
Picric Acid	Р	Р	Р	Р	Р	Р	G	G	Р	Р	Р	Р	F	F	G	G			G	G		
Potassium Hydroxide (30%)	G	G	G	G	G	G	G	G	G	G	Р	Р	F	F	G	G	G	G	G	G	G	G
Potassium Permanganate	G	G	G	G	G	F	G	G	G	F	G	G	G	F	F	F	F	Р	G	G	F	F
Propylene Glycol	G	G	G	G	G	G	G	G	F	F	F	F	G	G	G	G			G	G	G	G
Pyridine	Р	Р	Р	Р	G	G	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	F
Salicylic Acid	G	G	G	G	G	G	G	G	F	F	G	F	G	F	G	F			G	G	Р	Р
Silver Nitrate	G	F	G	G	G	G	G	G	G	F	G	G	G	G	G	F			G	G	F	F
Sodium Hydroxide (50%)	F	F	G	G	G	G	G	G	G	G	Р	Р	G	F	G	F	G	Р	G	G	G	G
Sodium Hypochlorite (15%)	G	F	G	F	F	Р	G	G	G	G	F	F	G	F	G	G	G	G	G	G		
Sulphuric Acid (20%)	G	G	G	G	G	G	G	G	G	F	G	F	G	G	G	F	G	G	G	G		
Sulphuric Acid (60%)	G	F	G	F	F	F	G	F	G	F	F	F	G	F	Р	Р	Р	Р	G	G	Р	Р
Sulphuric Acid (98%)	F	F	F	F	F	G	F	F	F	Р	Р	Р	F	F	Р	Р	Р	Р	G	G	Р	Р
Tetrahydrofuran (THF)	F	Р	F	G	F	F	F	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	F	F
Toluene	F	Р	Р	Р	Р	Р	F	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	Р	Р
Trichloroacetic Acid	F	Р	F	Р	F	F	G	G	F	Р	F	Р	F	Р	Р	Р			G	G		
Trichloroethylene	Р	Р	Р	Р	Р	Р	G	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	Р	Р
Turpentine	F	Р	F	Р	F	Р	F	Р	F	F	F	Р	Р	Р	F	Р	F	Р	G	G	G	G
Xylene	Р	Р	F	Р	Р	Р	F	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	Р	G	G	G	G
Zinc Chloride (10%)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	F
Zinc Sulphate (10%)	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	G	F	Р

# **Physical Properties of Plastics**

	LDPE	HDPE	PP	PMP (TPX)	PC	PS	PMMA	PTFE
Max Usage Temp. °C	80	120†	135	145	135	90	50	260
Min Usage Temp. °C	-50	-100	-20‡	-180	-135	-40	-60	-200
Transparency	Translucent	Translucent	Translucent	Clear	Clear	Clear	Clear	Opaque
Flexability	Flexible	Rigid	Rigid	Rigid	Rigid	Rigid	Rigid	Rigid
Specific Gravity	0.92	0.95	0.9	0.83	1.2	1.05	1.2	2.2
Microwavable	Yes	No	Yes	Yes	Marginal	No	No	Yes
Sterilization								
Autoclavable	No	No	Yes	Yes	Yes	No	No	Yes
Gas	Yes	Yes	Yes	Yes	Yes	Yes	No	Yes
Dry Heat	No	No	No	Yes	No	No	No	Yes
Radiation (Gamma)	Yes	Yes	No	No	Yes	Yes	Yes	Yes
Disinfection	Yes	Yes	Yes	Yes	Yes	Some	Some	Yes

<sup>†</sup> Polymer may be malleable at temperatures above 80°C if the product is under structural stress

#### Note:

The information provided in the Chemical Resistance of Plastics chart is intended to provide guidelines for reference only. Observe the tolerated temperature ranges of plastic when autoclaving laboratory plasticware. Remove any stoppers, fittings or caps from the plasticware prior to autoclaving. Plastic vessels should be autoclaved separately from their closures and other fittings. Autoclaving vessels with their closures in place can lead to deformation and destruction of the vessels. All statements are advisory only and imply no liability on the part of Heathrow Scientific.

#### **Sterilization:**

- Autocalving (121°C, 16 psi for 15 minutes)—Clean and rinse items with distilled water before autoclaving. (Always completely disengage thread before autoclaving.) Certain chemicals which have no appreciable effect on resins at room temperature may - Radiation—gamma irradiation at 25 kGy (2.5 MRad) with unstabilized plastic. cause deterioration at autoclaving temperatures unless removed with distilled water beforehand.
- Gas—Ethylene Oxide, formaldehyde, hydrogen peroxide.

- Dry heat (160°C, 120 minutes)
- Disinfectants—Benzalkonium chloride, formalin/formaldehyde, ethanol, etc.

<sup>‡</sup> Material may become brittle at low temperatures



